


Unison UK-C1
 Unison UK-C18
 Cadenza CD-C18

50 x 3 mm

Application

Tizanidine (muscle-relaxant) impurity

チザニジン(筋緊張緩和剤) 不純物

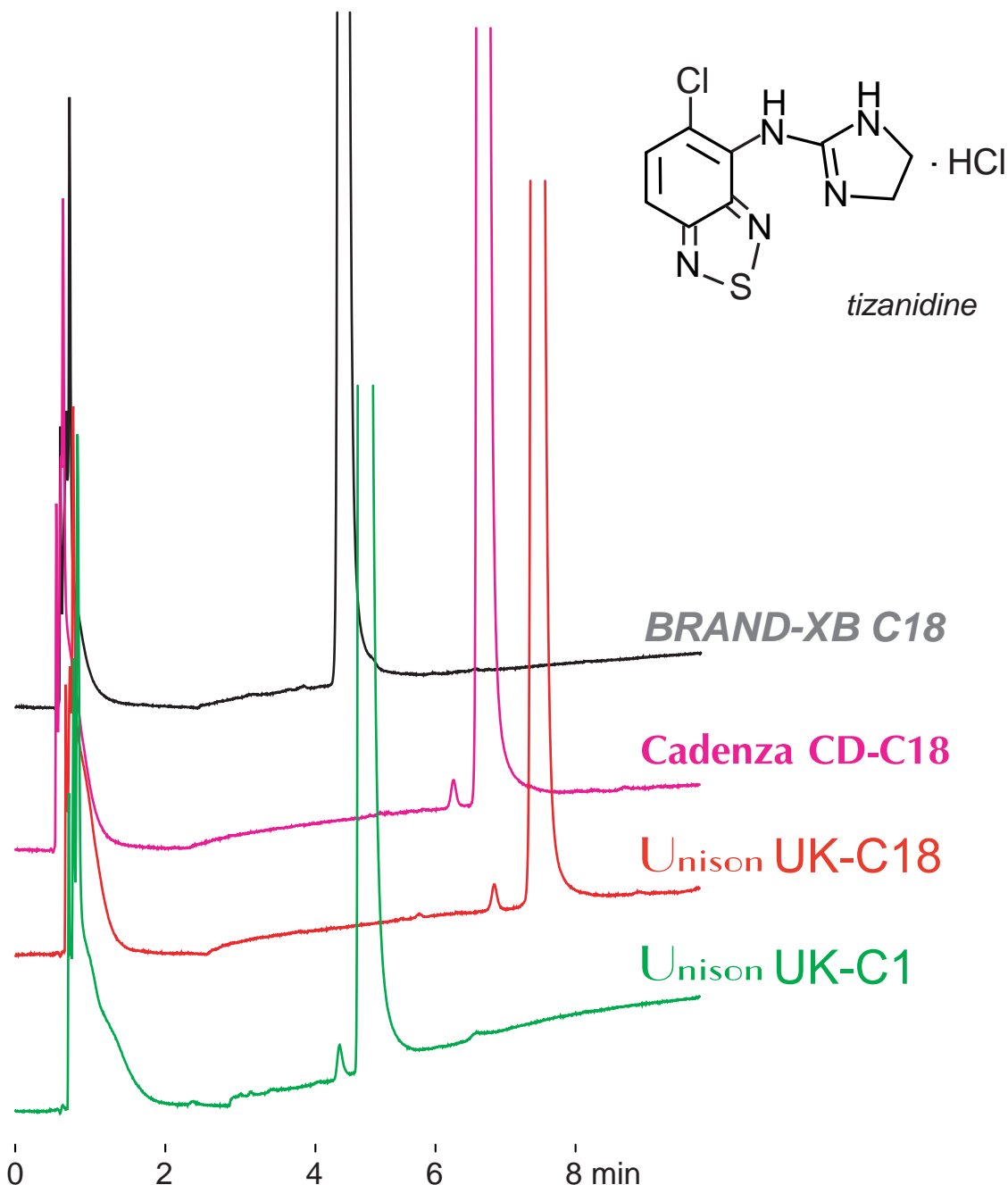

50 x 3 mm
 A: 10mM HCOONH₄
 B: 50mM HCOONH₄ / acetonitrile = 80 / 20
 0 - 10%B (0 - 10min)
 0.4mL/min (5MPa), 37°C
 2 μL (0.1mg/mL, DMSO), 320nm