

Cadenza CD-C18 Cadenza CW-C18



Unison UK-C18 Unison UK-C8

Cadenza CX-C18 Unison UK-Phenyl 50 x 2 mm

Application

## LC-MS: Omega-3 fatty acid ethyl esters

オメガ-3脂肪酸エチルのLC-MS分析

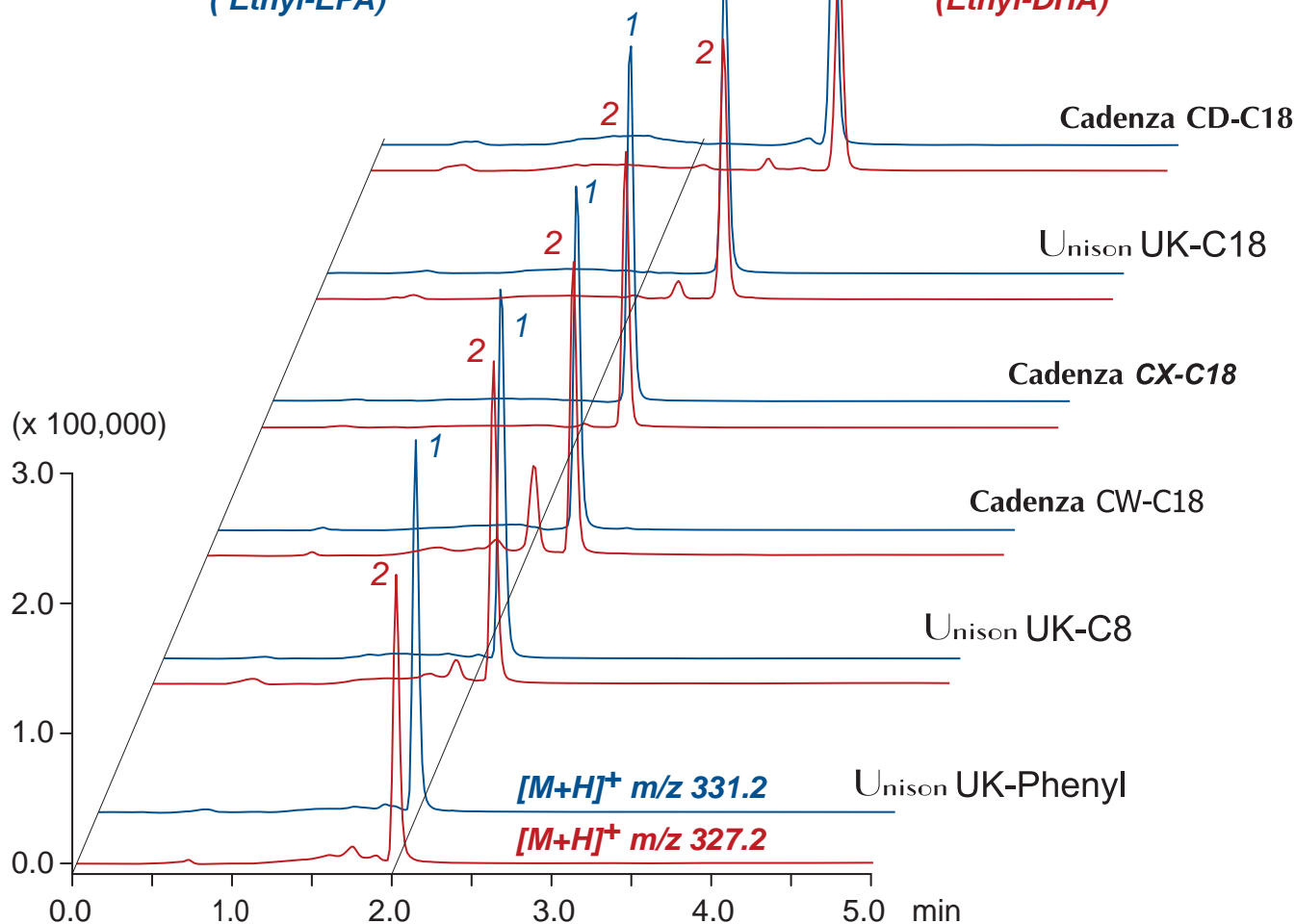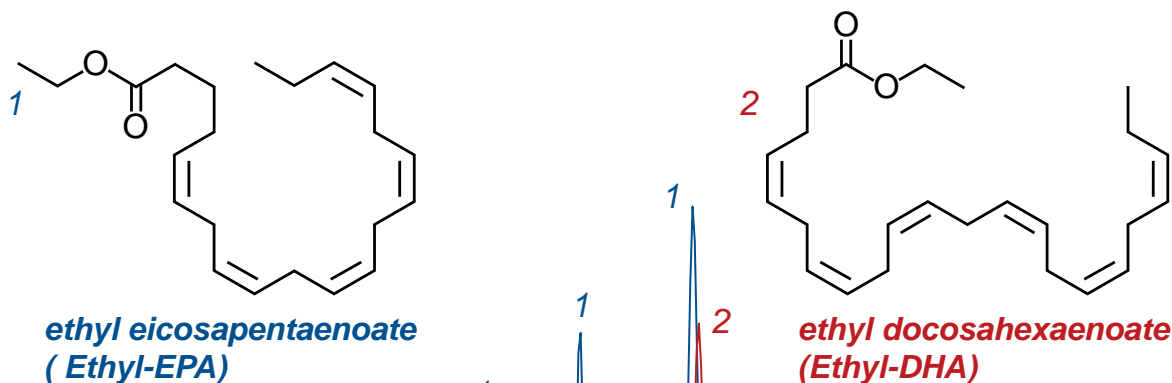

50 x 2 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

30-100%B (0-1min), 100%B (1-5min)

0.2mL/min (4MPa), 37°C, MS (APCI, Positive)

5μL (20 μg/mL, EtOH)