

**Scherzo SS-C18**  
Unison UK-Amino

50 x 3 mm

Application

### Amantadine

(Antiviral, Neuropsychiatric agent, Transporter substrate)

アマンタジン

(抗ウイルス薬、精神神経用薬、トランスポーター基質薬)

#### RP+ IEX

Scherzo SS-C18, 50 x 3 mm


A: 0.1% HCOOH

B: 100mM HCOONH<sub>4</sub> / ACN = 50 / 50

0-100%B (0-10min)

0.4mL/min (5MPa), 37°C, ELSD

2 μL (1mg/mL, DMSO)

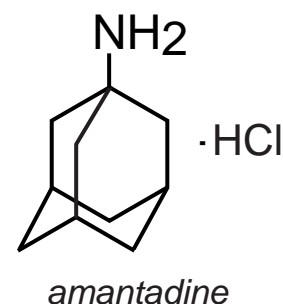

#### NP

Unison UK-Amino, 50 x 3 mm

A: acetonitrile

B: acetonitrile / 100mM HCOONH<sub>4</sub> = 50 / 50

0-100%B (0 -10min)

0.4mL/min (3MPa), 37°C, ELSD

2 μL (1mg/mL, DMSO)

