


Cadenza CD-C18
Cadenza CX-C18
Unison UK-C18

50 x 3 mm


Application

LC-MS: Omarigliptin (Oral antidiabetic drug)

オマリグリプチン(経口糖尿病用剤)のLC-MS分析

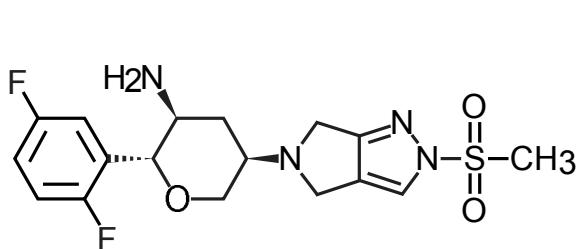

omarigliptin

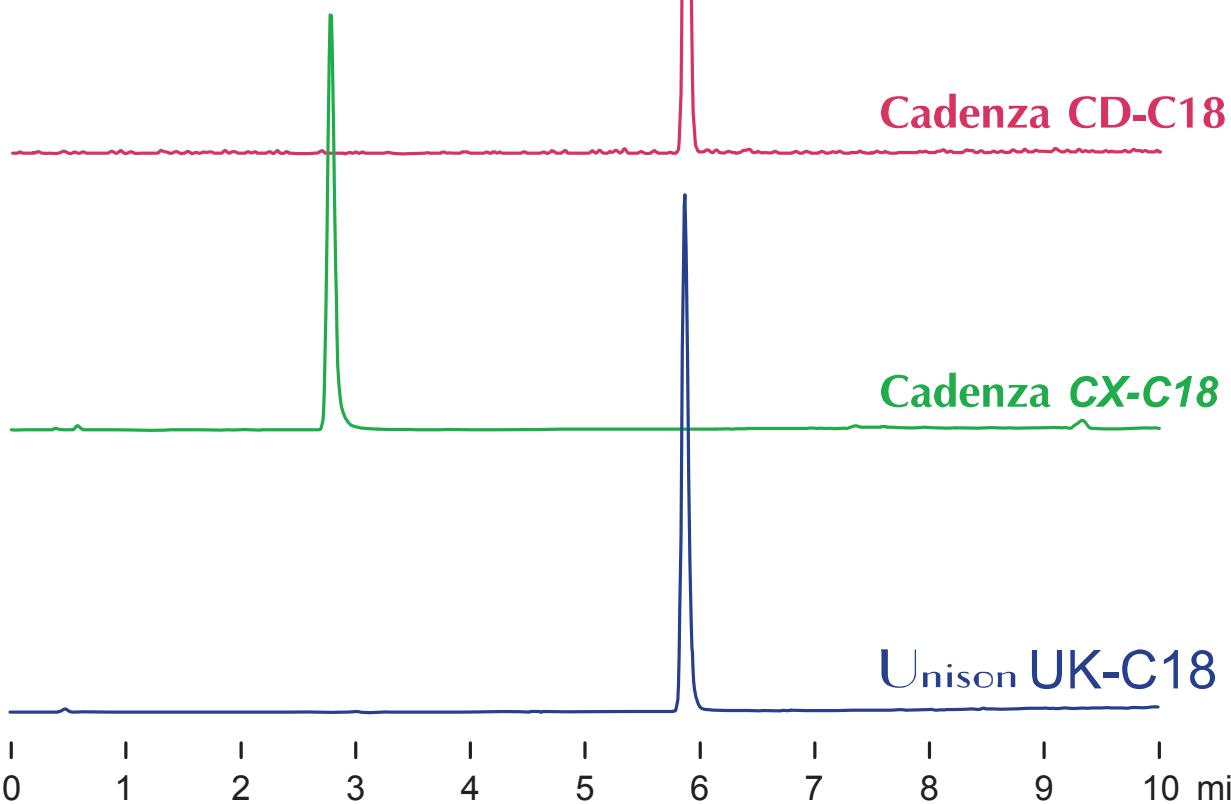 $[M+H]^+$ m/z 399.1

50 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-30% B (0-10min)

0.4mL/min (4MPa), 37°C, MS (ESI, SIM, positive)

5µL (1µg/mL, 0.1% HCOOH)