


**Intrada Organic Acid**

150 x 2 mm

Application

**LC-MS: Organic acids in Soy sauce**

醤油中の有機酸のLC-MS分析

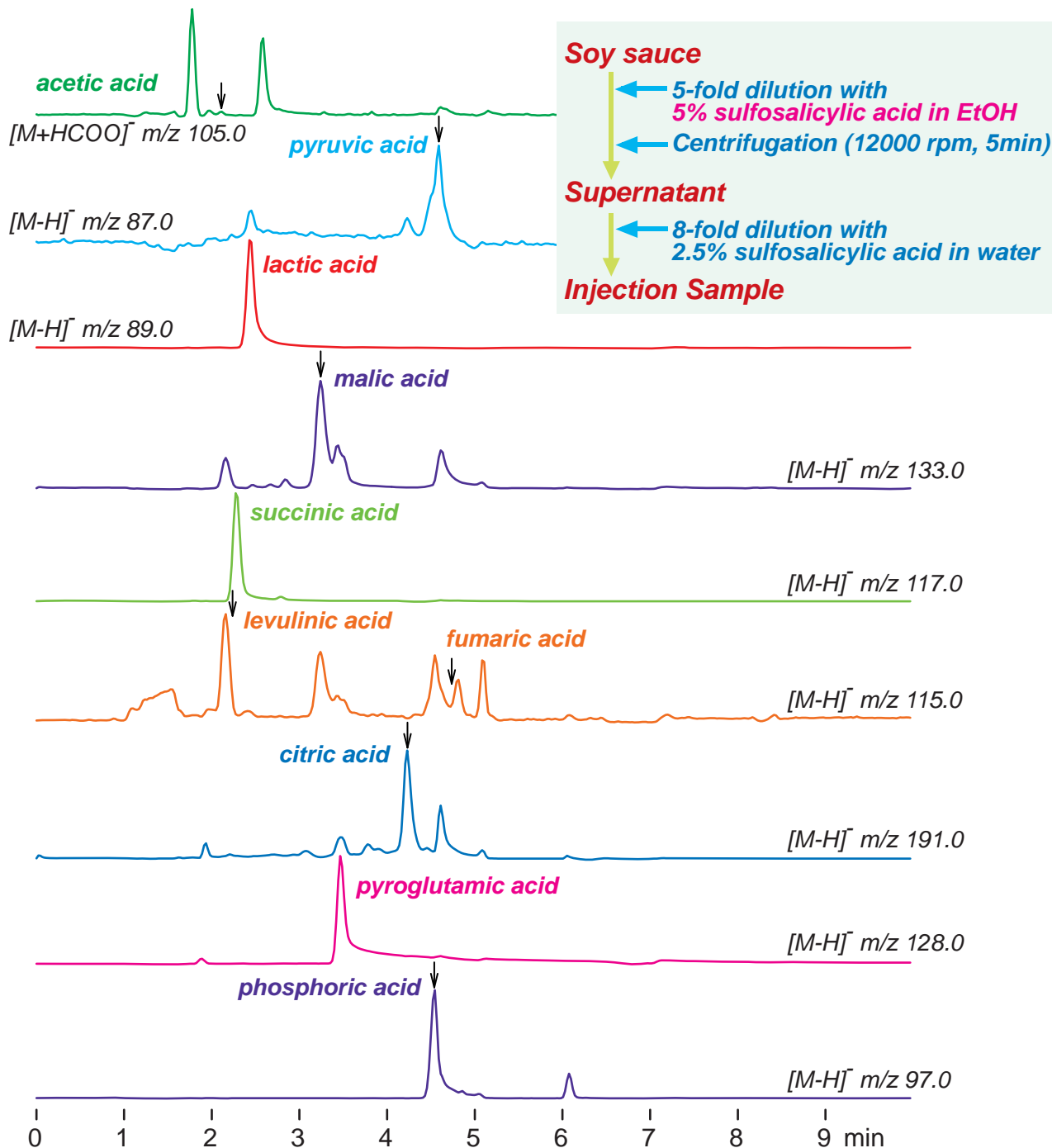

Intrada Organic Acid, 150 x 2 mm  
 A: acetonitrile / water / HCOOH = 10 / 90 / 0.1  
 B: acetonitrile / 100mM HCOONH<sub>4</sub> = 40 / 60  
 0%B (0-1min), 0-100%B (1-7min), 100%B (7-10min)  
 0.2mL/min (12MPa), 37°C, 5μL  
 MS (ESI, SIM, Negative)