




Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

50 x 3 mm

Application

LC-MS: Nelfinavir mesylate (antiretroviral drug)

ネルフィナビル メシル酸塩 (抗HIV薬) の LCMS分析

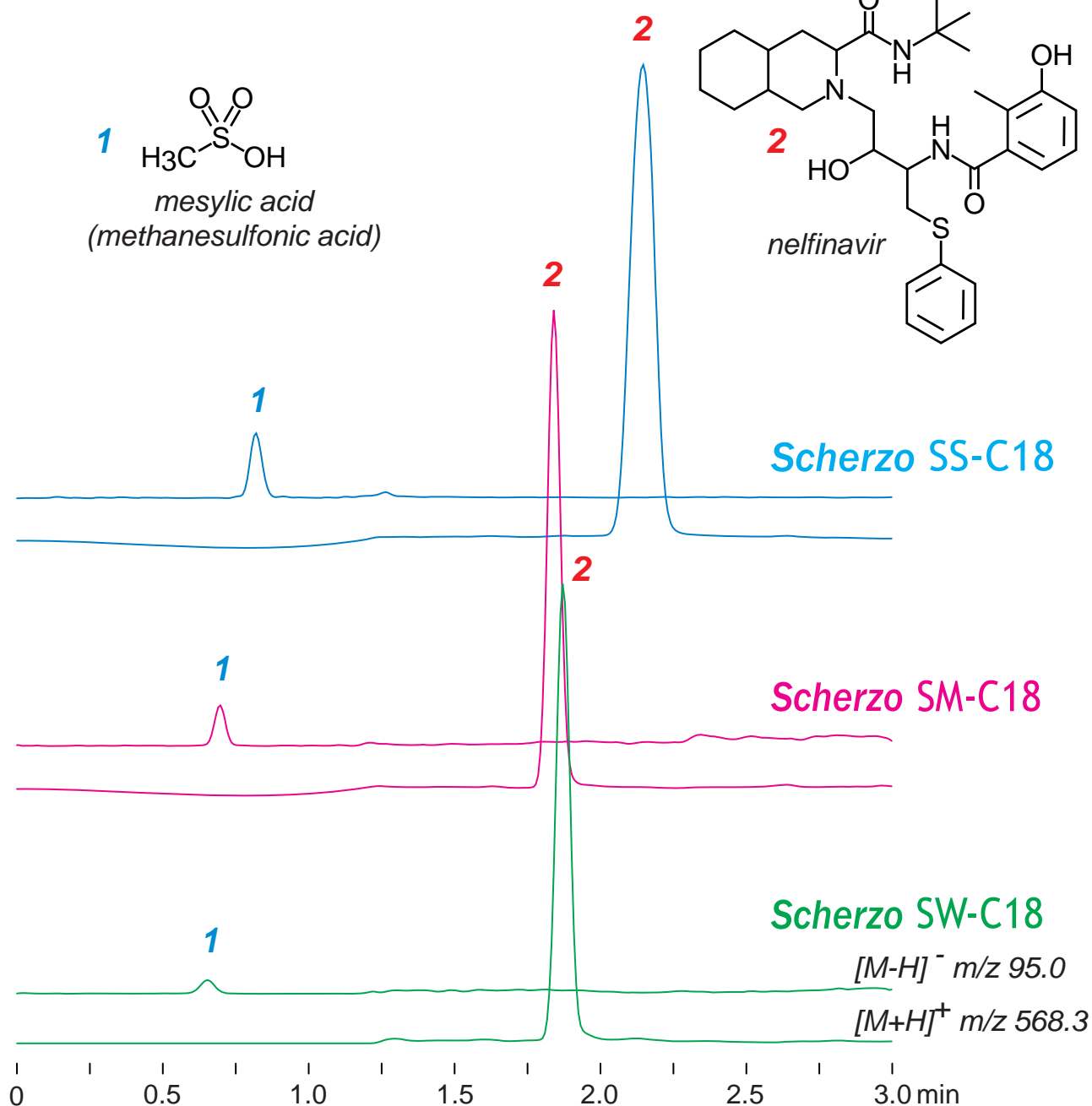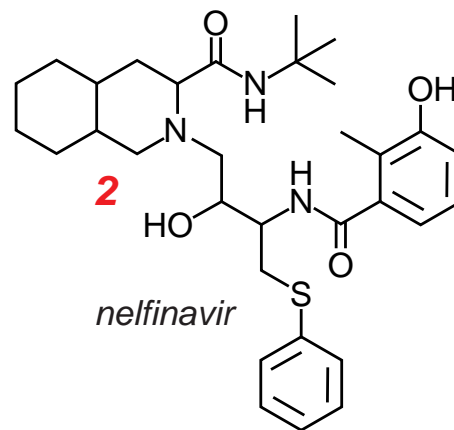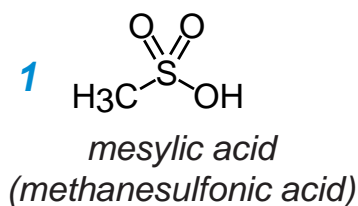

50 x 3 mm

A: 50mM ammonium formate

B: acetonitrile

0-100%B (0-1 min), 100%B (1-3 min)

0.4mL/min (4MPa), 37°C, MS (ESI, SIM)

5µL (2µg/mL, 50mM HCOONH₄ / ACN =90 / 10)