



Unison UK-C1  
 Unison UK-Phenyl  
 Unison UK-C8  
 Unison UK-C18

50 x 3 mm

Application

**LC-MS: Ivermectin (antiparasitic drug)**

イベルメクチン(抗寄生虫薬)のLC-MS分析

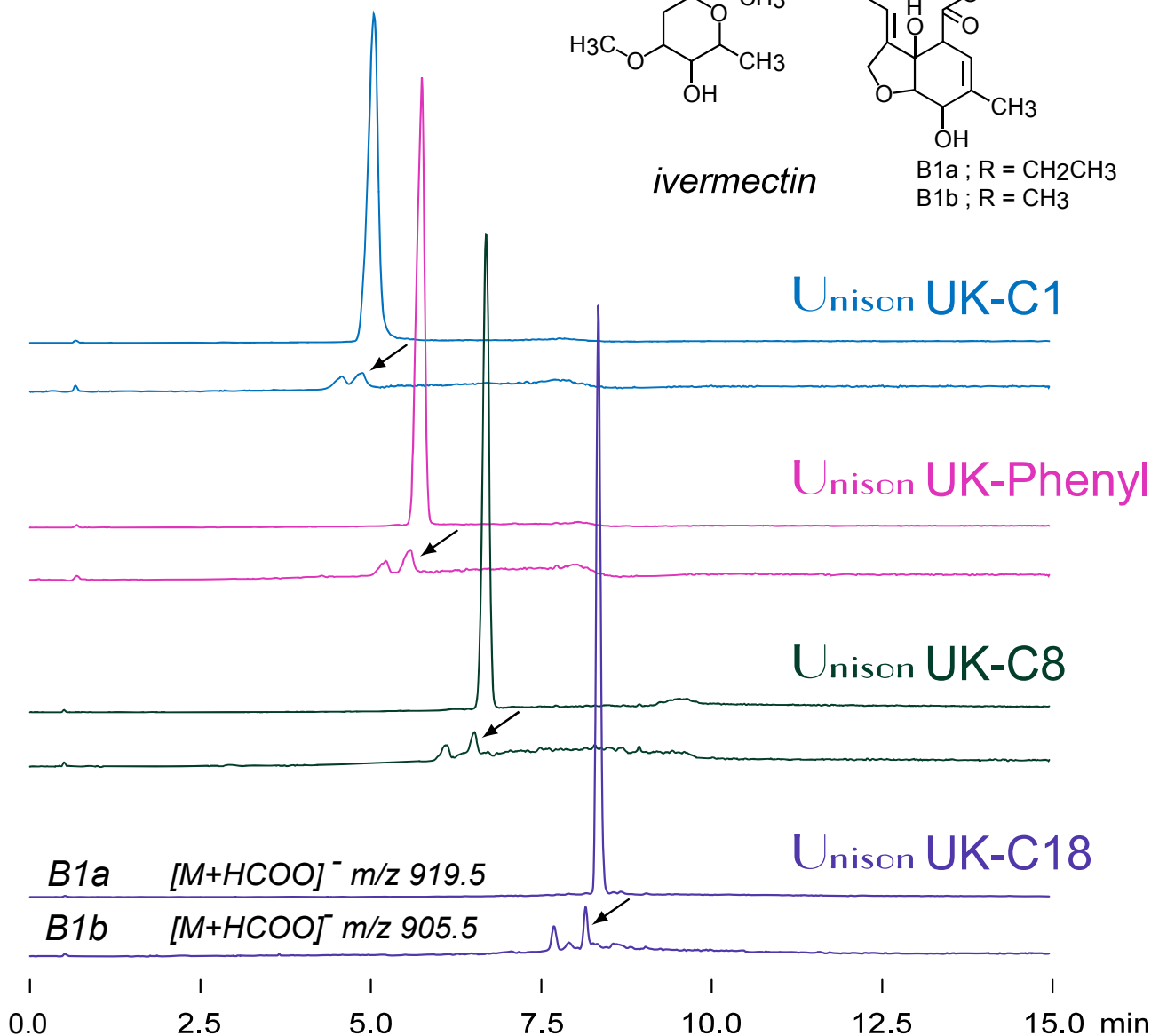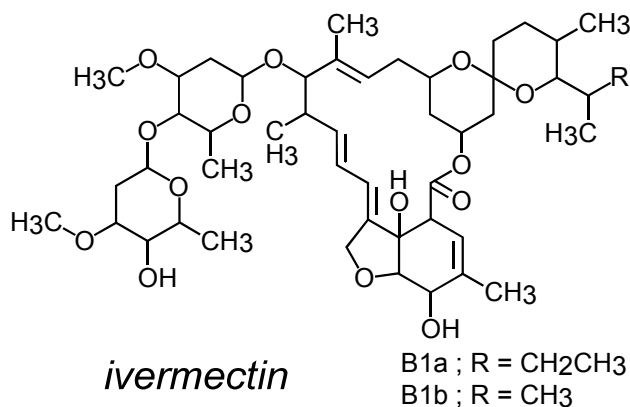

50 x 3 mm

A: water / formic acid = 100 / 0.1

B: methanol / formic acid = 100 / 0.1

50-100%B (0-8 min), 100%B (8-15 min)

0.4mL/min (6MPa) , 37°C, MS (ESI, SIM)

5µL (10µg/mL, EtOH / water / HCOOH = 95 / 5 / 0.1)