



Scherzo SS-C18  
 Scherzo SM-C18  
 Scherzo SW-C18

100 x 4.6 mm

Application

*Oral hypoglycemic drug (metformin)*

ビグアナイド系経口糖尿病治療薬 (メトホルミン)

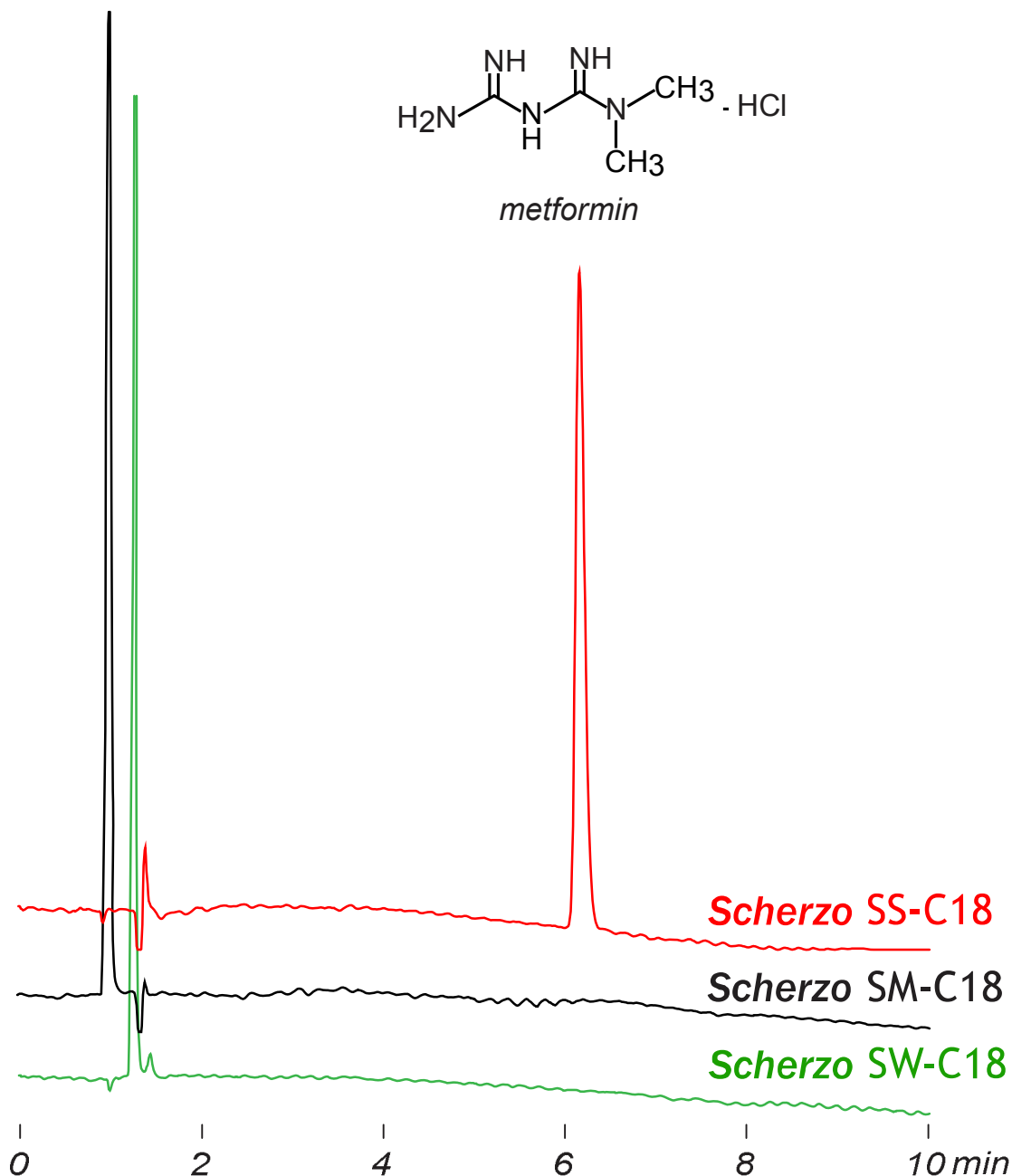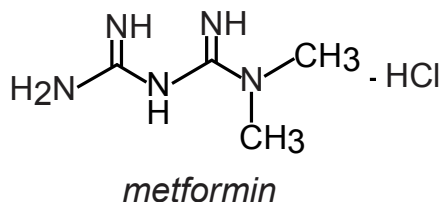

100 x 4.6 mm  
 A: acetonitrile / water / formic acid = 5 / 95 / 0.1  
 B: acetonitrile / 100mM HCOONH4 = 50 / 50  
 40%B (0-1min), 40-60%B (1-10min)  
 1 mL/min (8MPa), 37°C  
 5 µL (0.1mg/mL, water), 230nm