

Unison UK-Amino

Unison UK-Silica



250 x 4.6 mm

Application

Oral hypoglycemic drug (metformin)

ビグアノイド系経口糖尿病治療薬 (メトホルミン)

Unison UK-Amino

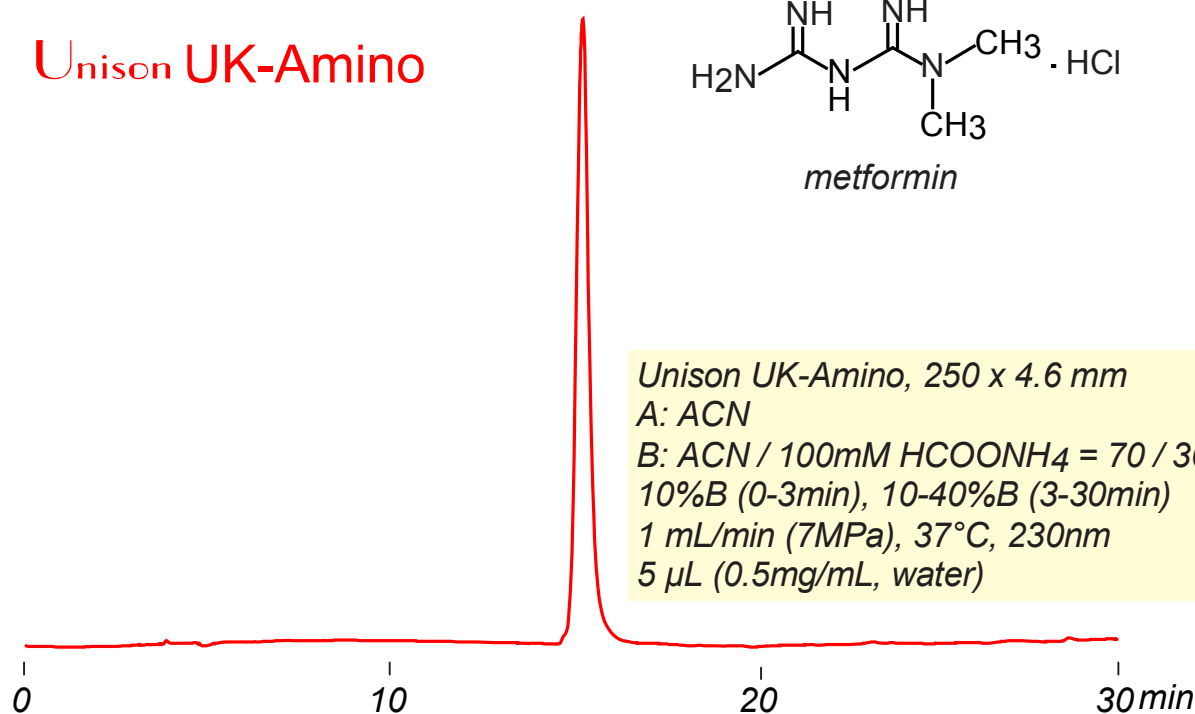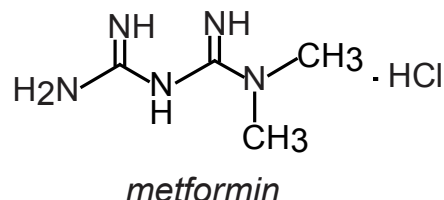

Unison UK-Silica

Unison UK-Silica, 250 x 4.6 mm
 A: ACN / water / HCOOH = 90 / 10 / 0.1
 B: ACN / water / HCOOH = 30 / 70 / 0.15
 15%B (0-3min), 15-80%B (3-30min)
 1mL/min (7MPa), 37°C, 230nm
 5 μL (0.5mg/mL, water)

