



Cadenza **CX-C18**  
Unison **UK-C8**  
Cadenza **CD-C18**  
Cadenza **CW-C18**  
Unison **UK-C18**

75 x 4.6 mm

Application

**Enrofloxacin (fluoroquinolone)**

エンロフロキサシン (キノロン系抗菌薬)

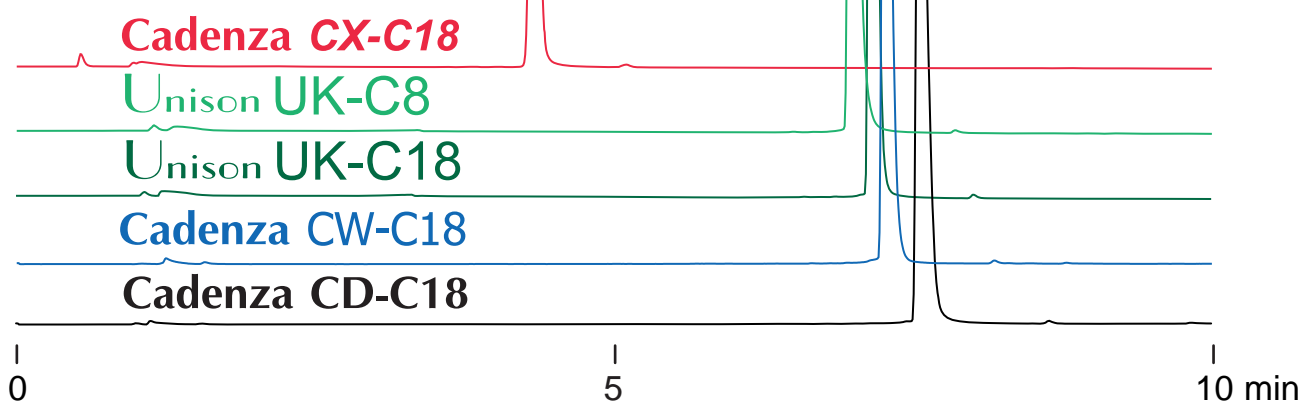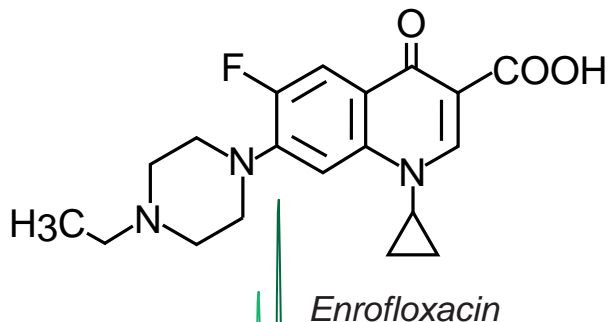

75 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-30%B (0-10min)

1.0 mL/min (7MPa), 37°C, 290nm

10  $\mu$ L (0.1mg/mL, 1%HCCOH)