

**Cadenza CX-C18**

**Cadenza CD-C18**

**Cadenza CW-C18**



**Unison UK-C18**

75 x 4.6 mm

Application

**Lomefloxacin (fluoroquinolone)**

ロメフロキサシン (キノロン系抗菌薬)

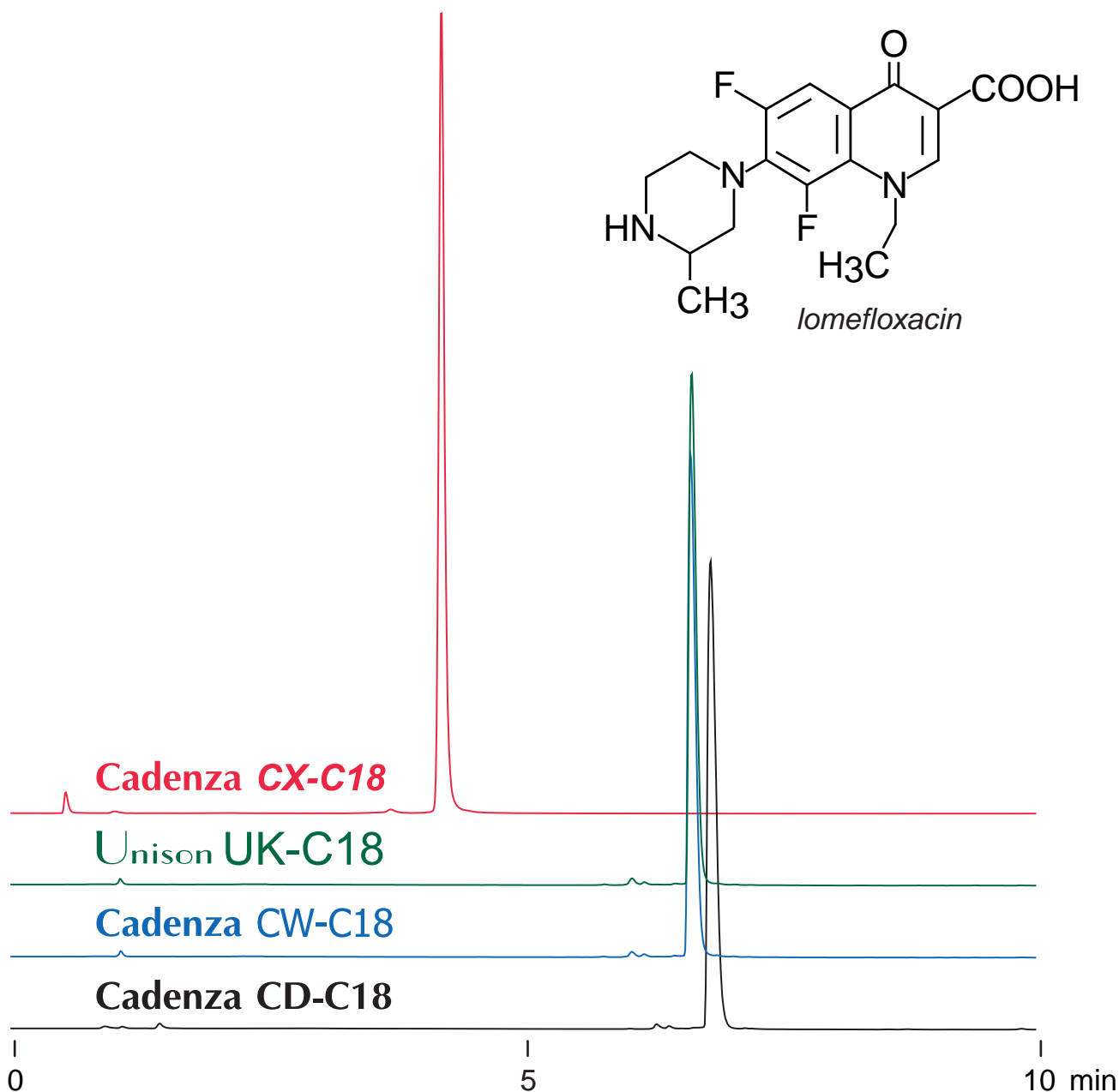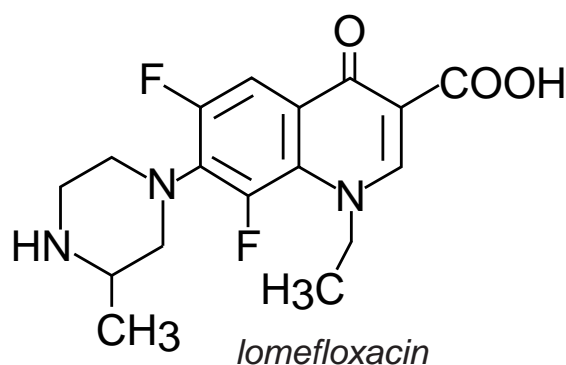

75 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-30%B (0-10min)

1.0 mL/min (7MPa), 37°C, 290nm

10 µL (0.1mg/mL, 1%HCCOH)