


Cadenza CW-C18

250 x 4.6 mm

Application

Galactosylceramide from bovine brain
 スフィンゴ糖脂質 (ガラクトセレブロシド, 牛脳由来)

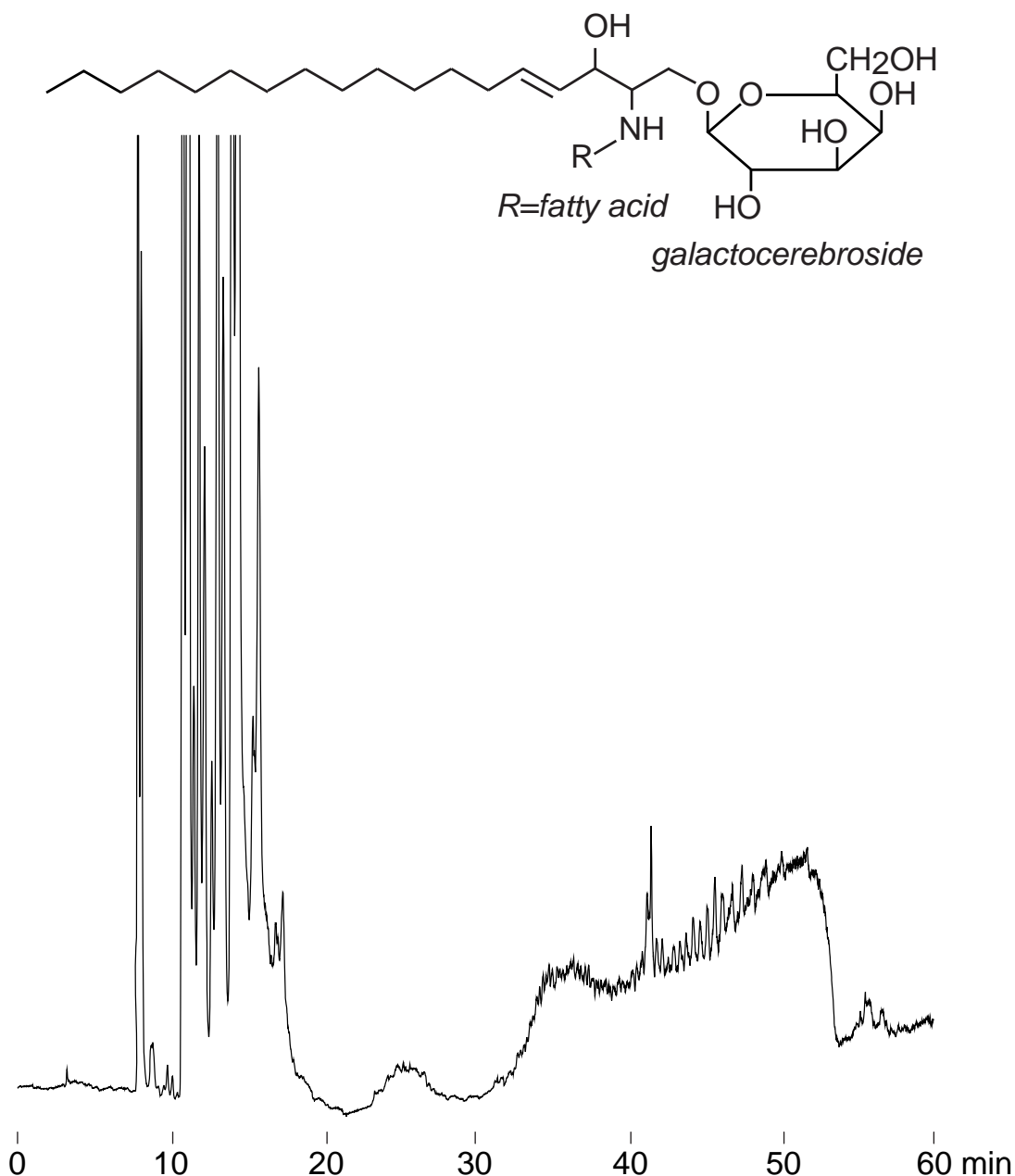

Cadenza CW-C18, 250 x 4.6 mm
 A: acetonitrile / formic acid = 100 / 0.3
 B: tetrahydrofuran
 10-15%B (0-12min), 15-100%B (12-50min), 100%B (50-60min)
 1mL/min (7MPa), 37°C, ELSD (spray chamber 50°C, drift chamber 100°C)
 5µL (4mg/mL, chloroform / methanol = 2 / 1)