

Unison UK-Phenyl

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

Application

## Lomerizine (Antimigraine, Calcium channel blocker)

ロメリジン ( 偏頭痛治療薬 , カルシウム拮抗薬 )

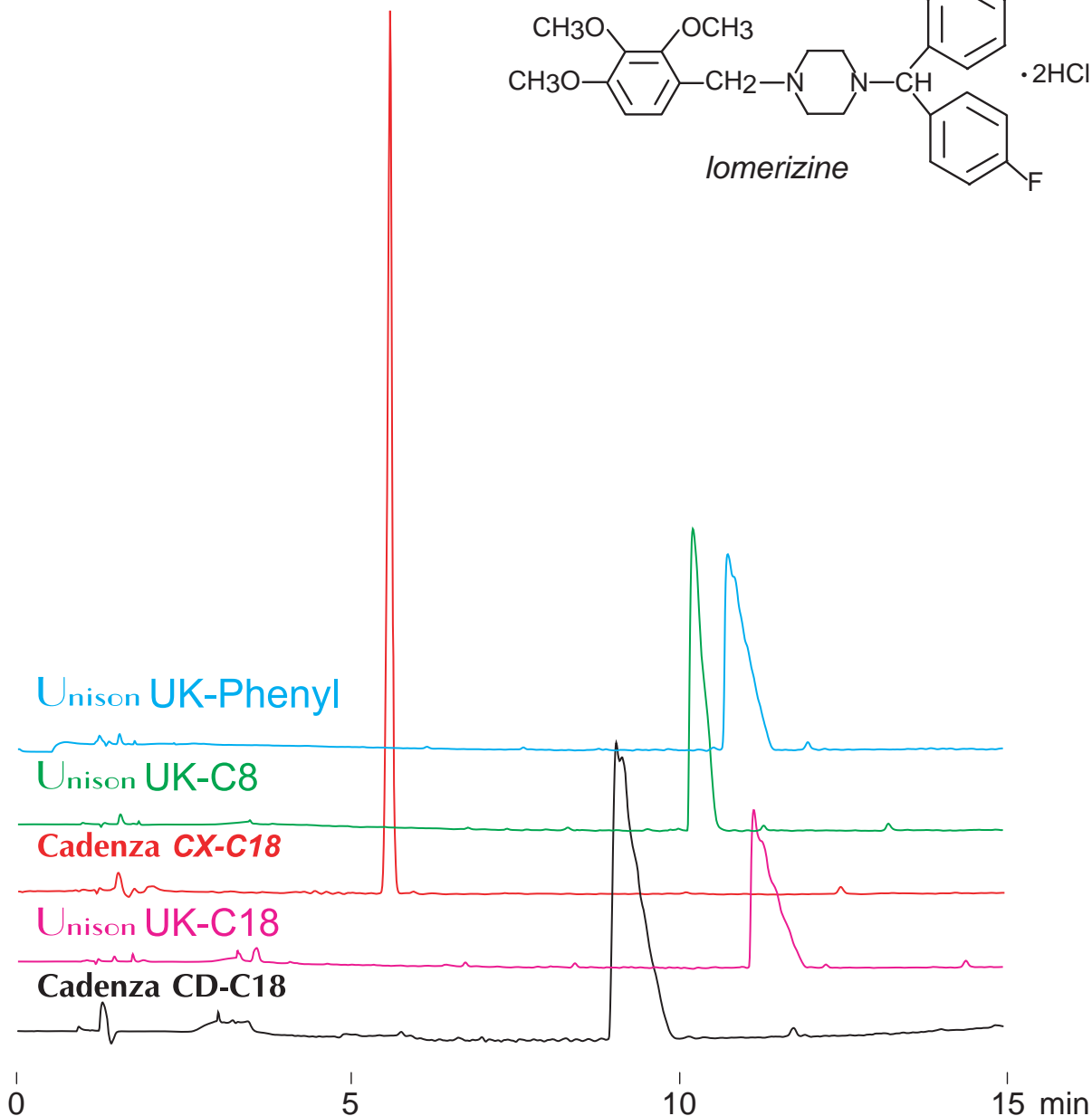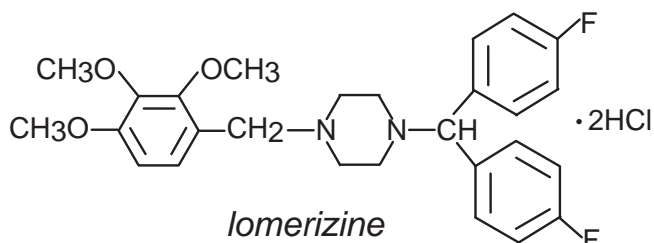

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-75%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (1mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1 )