

Unison UK-Phenyl



Unison UK-C8

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

Application

Lomerizine (Antimigraine, Calcium channel blocker)

ロメリジン (偏頭痛治療薬 , カルシウム拮抗薬)

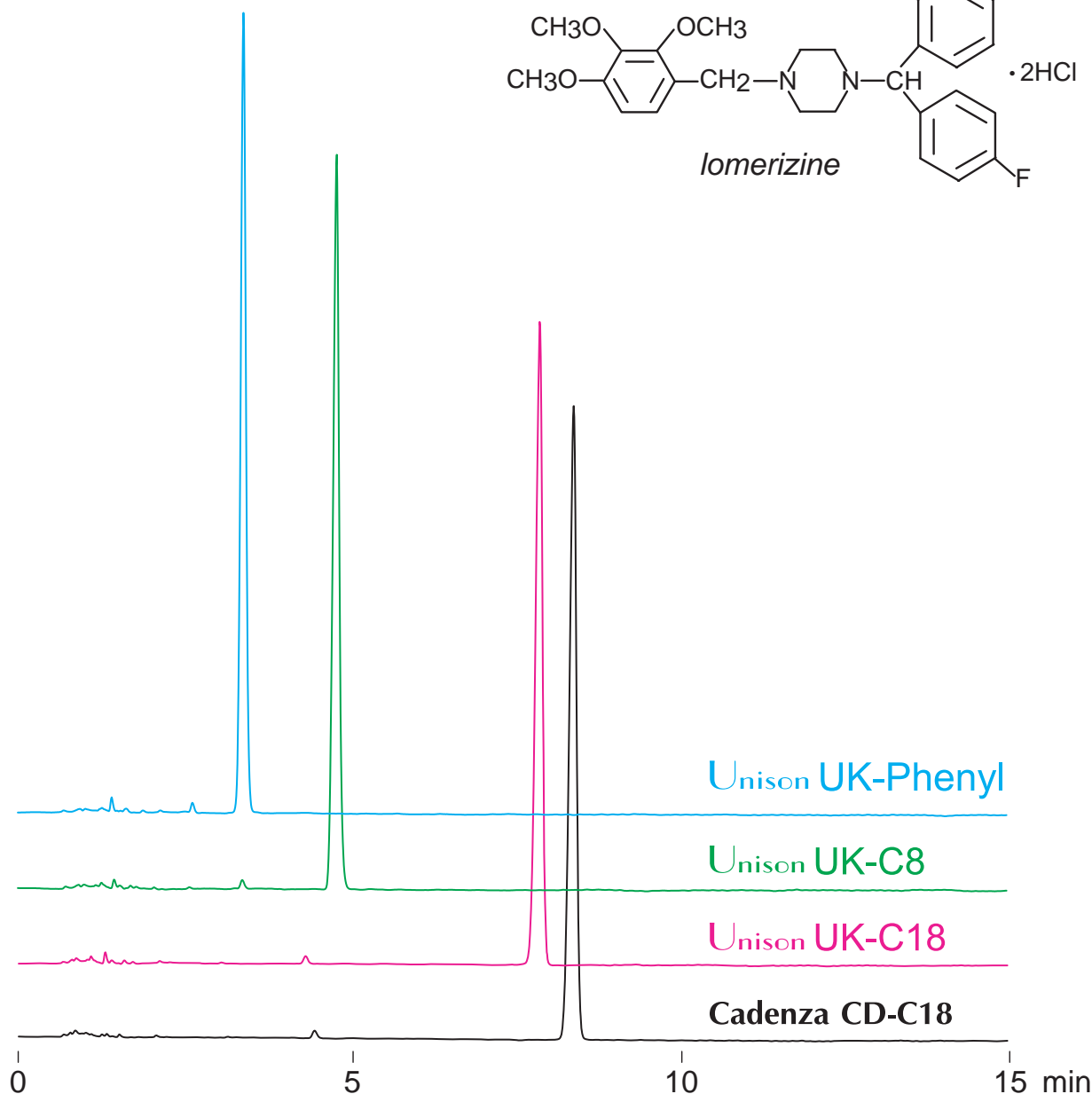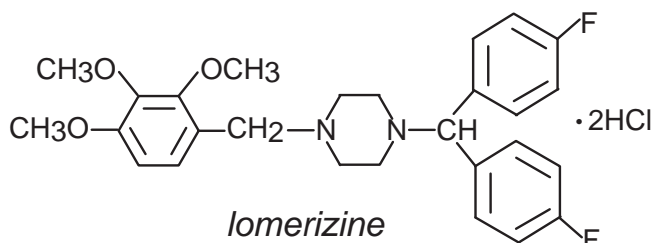

100 x 4.6 mm

A: acetonitrile / 100mM ammonium formate= 50 / 50

B: acetonitrile / 100mM ammonium formate= 90 / 10

20-100%B (0-15min)

1.0 mL/min (6MPa), 37°C, 260nm

10 µL (1mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)