

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


**HPLC COLUMN APPLICATION**

**Clotiazepam ( Antianxiety, GABA-A receptor agonist )**

クロチアゼパム ( 抗不安薬, GABA-A 受容体作動薬 )

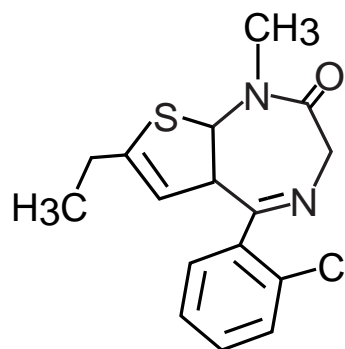

clotiazepam

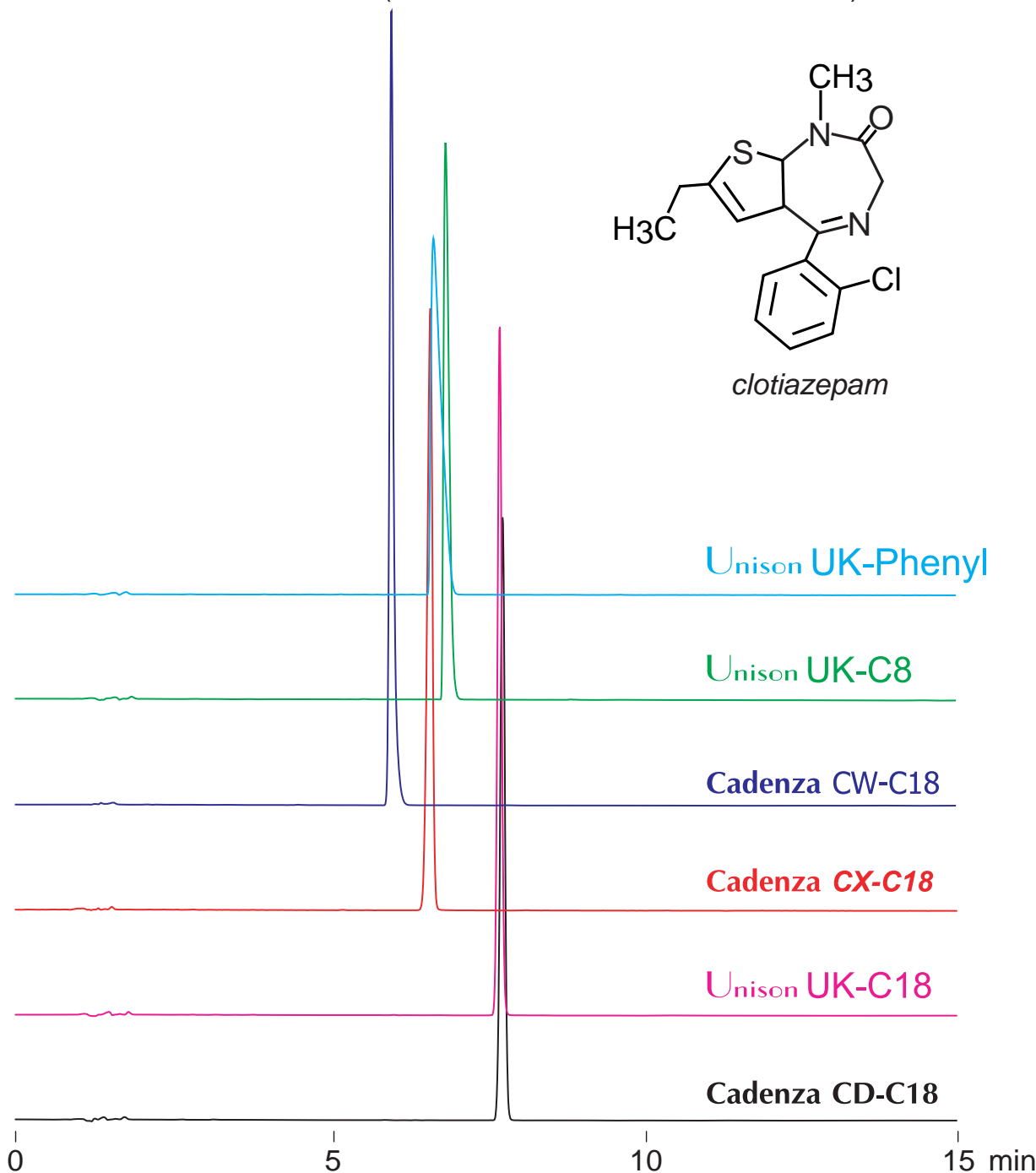

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

25-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (0.25mg/mL, methanol / formic acid = 100 / 1)