

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Raloxifene (Osteoporosis agent)

ラロキシフェン (骨粗しょう症治療薬)

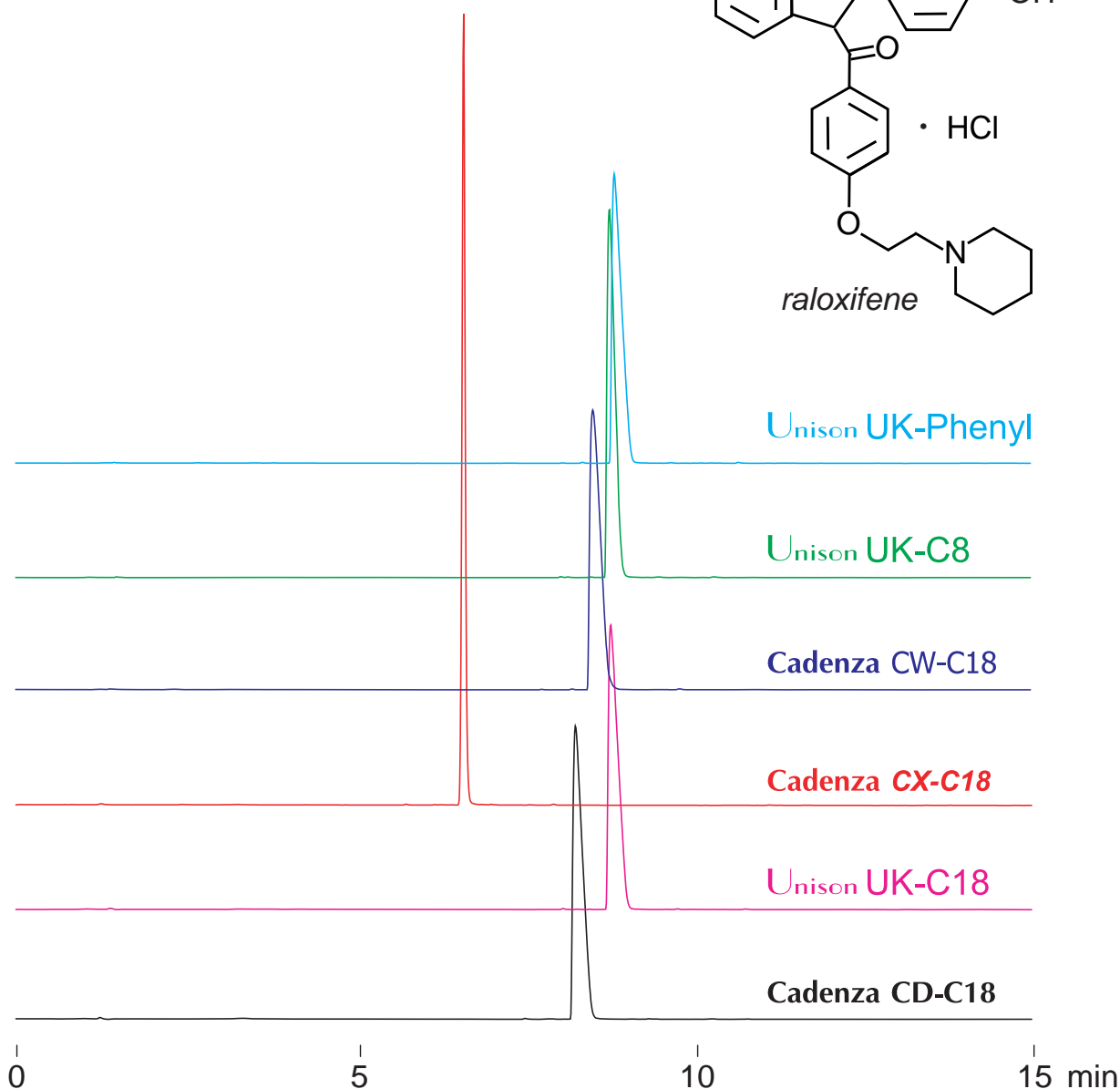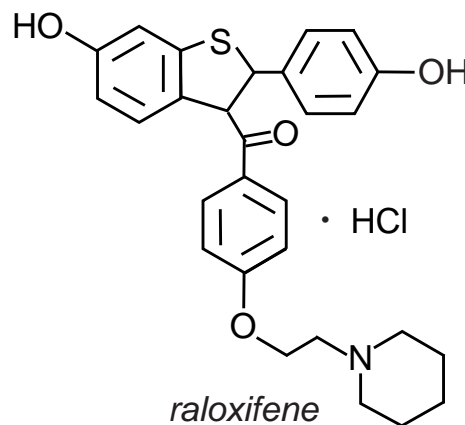

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-80%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.25mg/mL, methanol / water / formic acid = 50 / 50 / 1)