


DACAPO DX-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION


Bromhexine (Expectorant drug)

ブロムヘキシン (去痰薬)

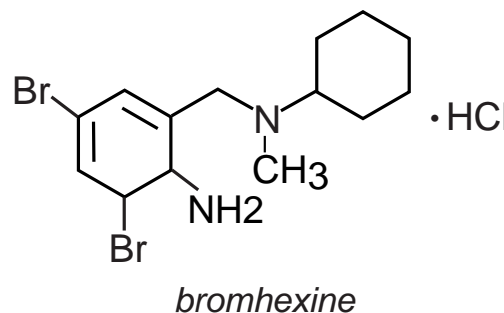

A: 100mM ammonium formate
 B: acetonitrile / 100mM ammonium formate
 = 90 / 10
 80-100%B (0-20min)
 16MPa

A: water / formic acid = 100 / 0.1
 B: acetonitrile / formic acid = 100 / 0.1
 0-100%B (0-20min)
 20MPa

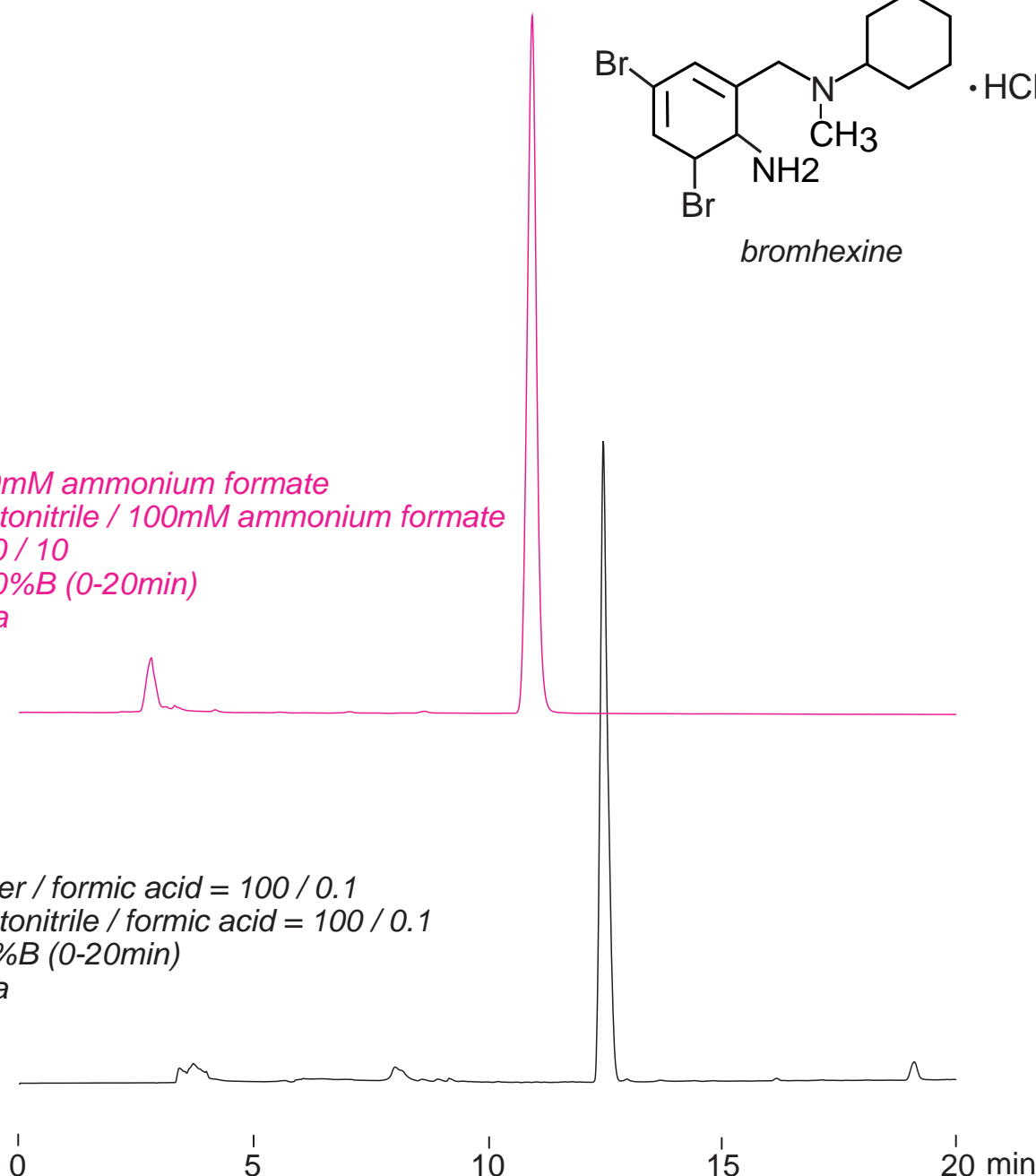

Dacapo DX-C18, 250 x 4.6 mm
 1.0mL/min, 37°C, 260nm
 10µL (1mg/mL, water / formic acid = 100 / 1)