

Unison UK-C8

Unison UK-Phenyl

Cadenza GL-C18

Cadenza CW-C18

Cadenza CD-C18



Unison UK-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Glycyrrhizic acid ( Anti-inflammatory )

グリチルリチン酸 ( 抗炎症薬 )

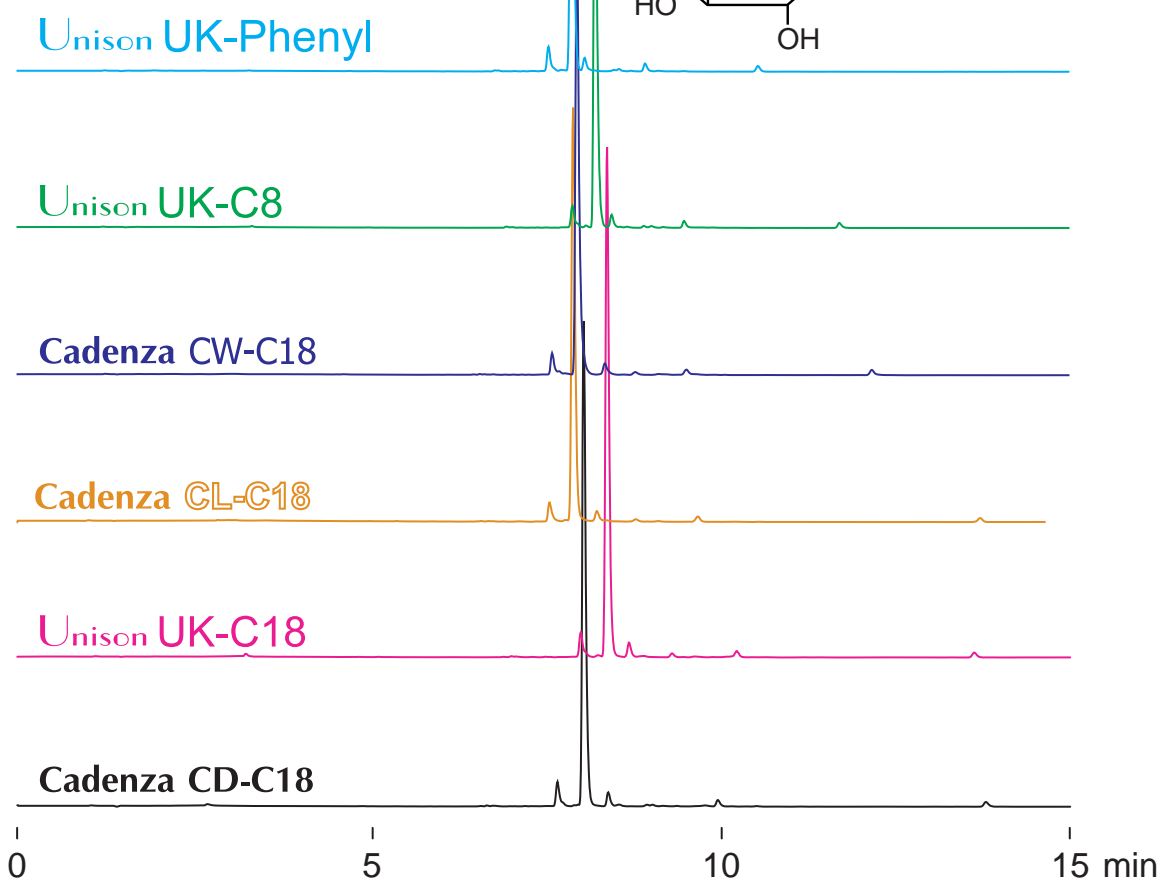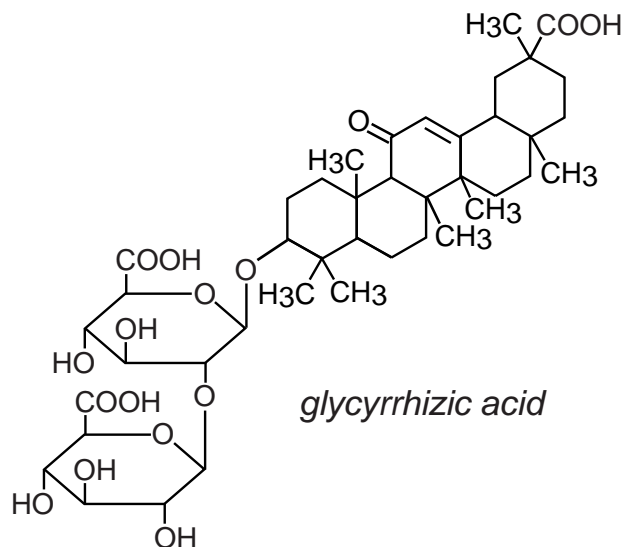

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-15min),

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, 100mM ammonium formate)