

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

### Picosulfate ( Laxative )

ピコスルファート ( 緩下薬 )

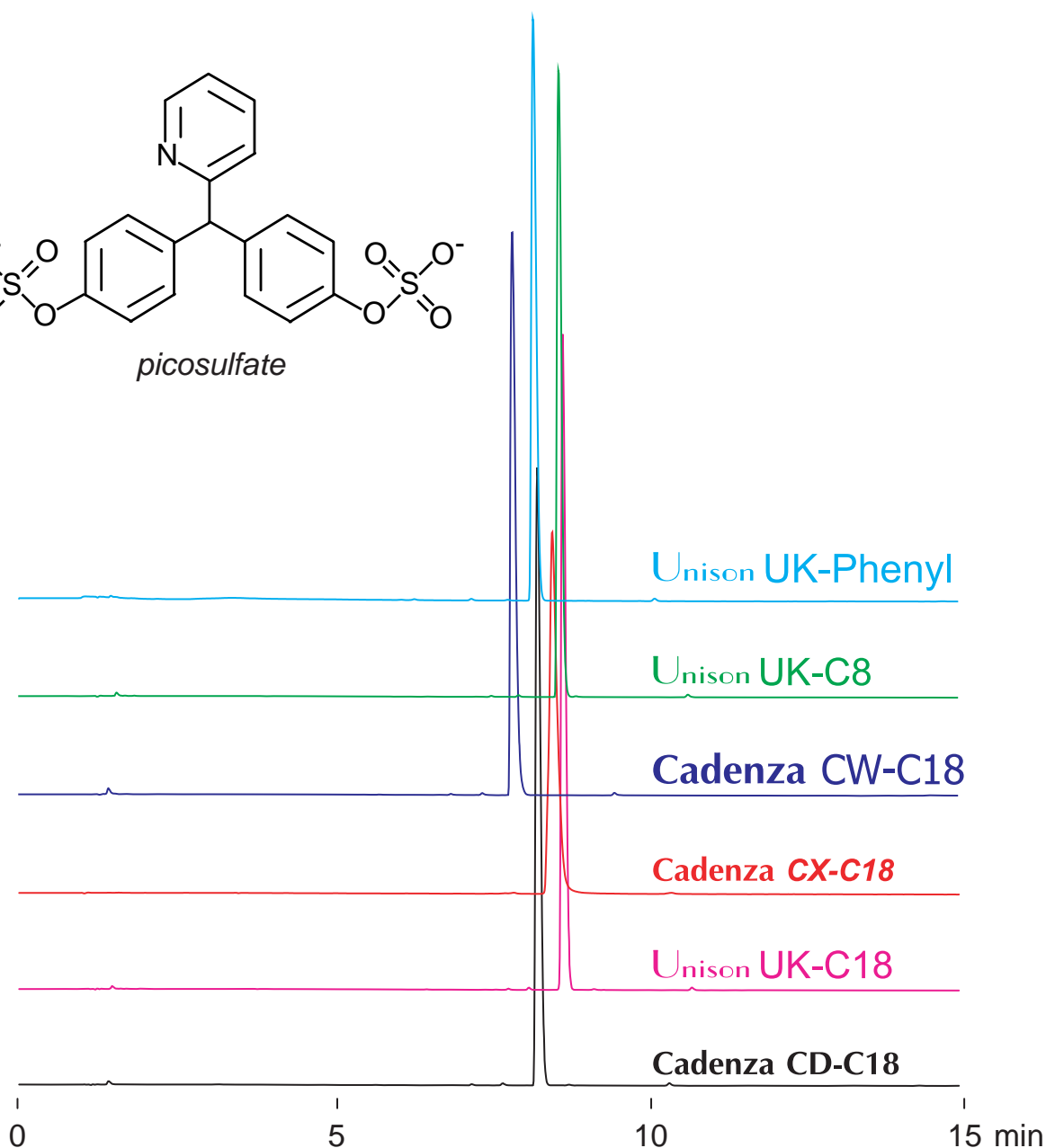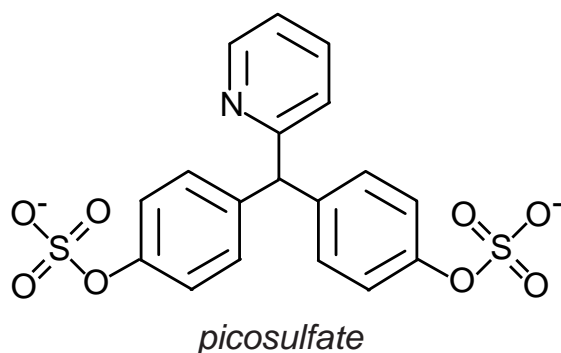

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

0-50%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (1mg/mL, water / formic acid = 100 / 0.1)