

Unison UK-C8

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Nematic, Smectic Liquid-crystal materials

ネマチック , スメクチック液晶材料

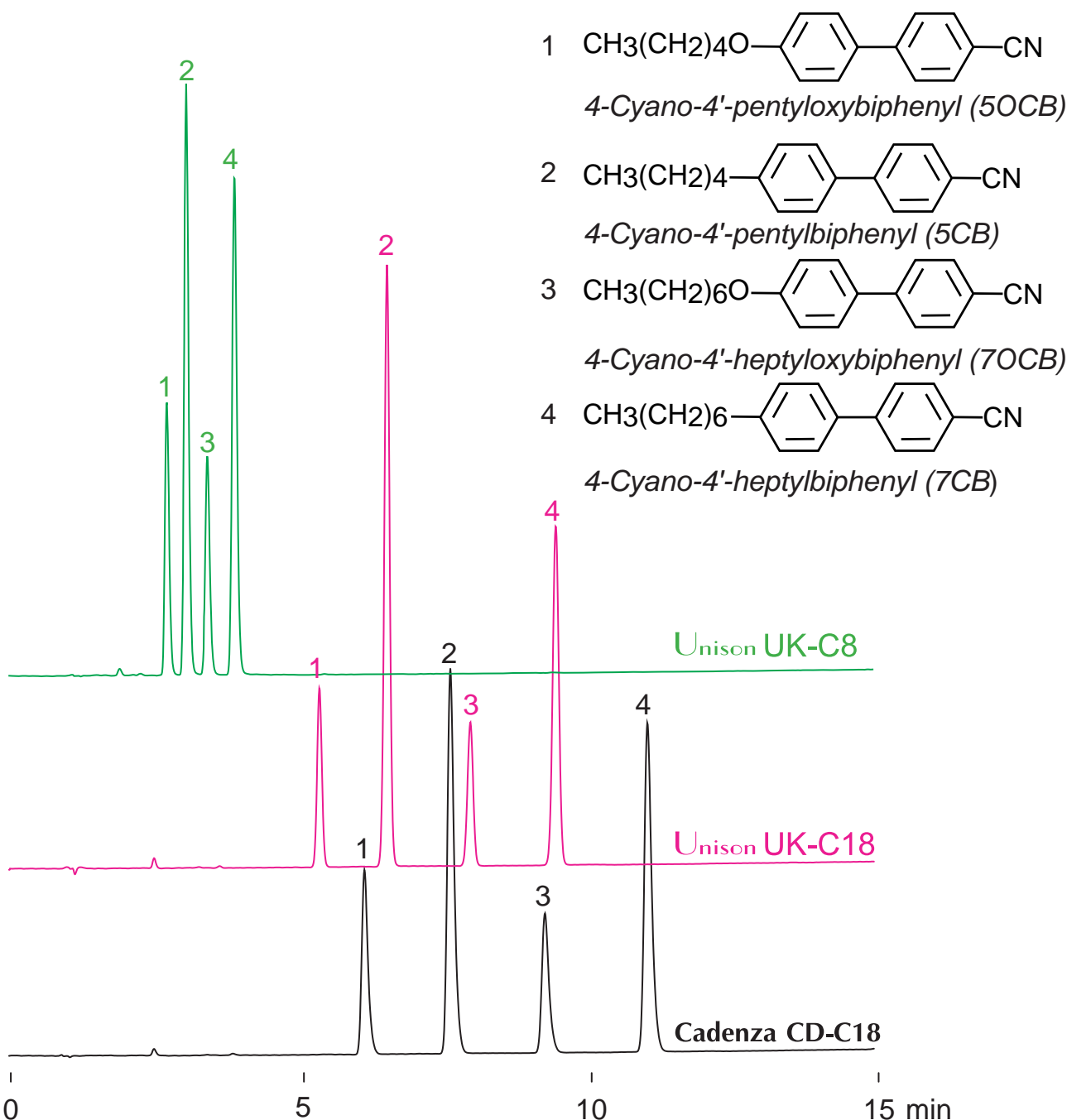

100 x 4.6 mm

A: acetonitrile / water / acetic acid = 50 / 50 / 0.5

B: acetonitrile / acetic acid = 100 / 0.5

50-100%B (0-15min)

1.0 mL/min (6MPa), 37°C, 260nm

10µL (0.1mg/mL, acetonitrile)