

Cadenza CD-C18



100 x 4.6 mm

250 x 4.6 mm

HPLC COLUMN APPLICATION

Teprenone ( Antiulcer drugs )

テプレノン ( 抗消化性潰瘍薬 )

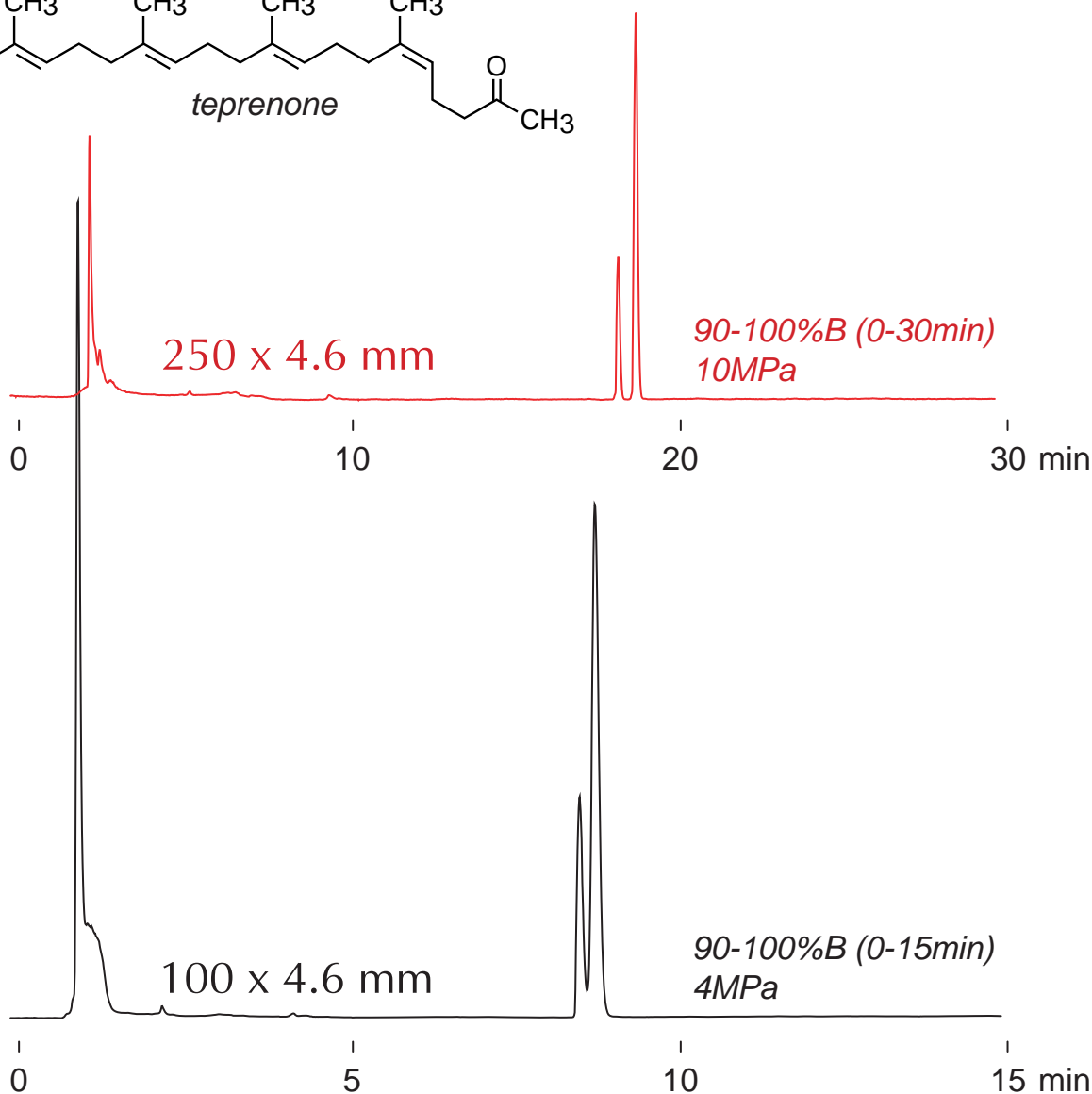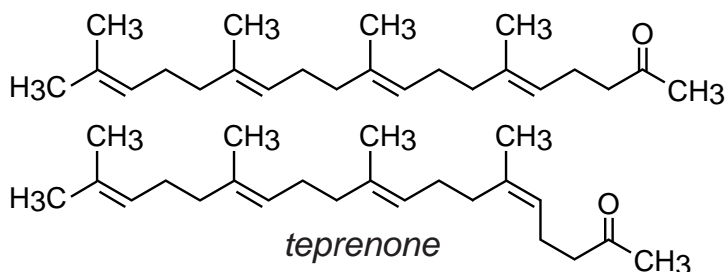

Cadenza CD-C18

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

1.0mL/min, 37°C

ELSD (spray chamber 30°C, drift chamber 60°C)

10µL (1mg/mL, dimethyl sulfoxide)