

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Unison UK-C18

Cadenza CX-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## Fexofenadine ( Allergic rhinitis, Opioid receptor agonist )

フェキソフェナジン ( アレルギー鼻炎薬 , ヒスタミン受容体拮抗薬 )

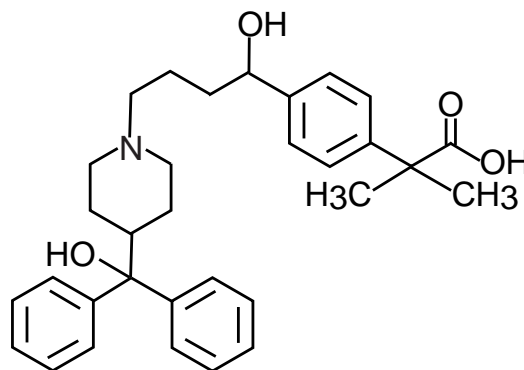

fexofenadine

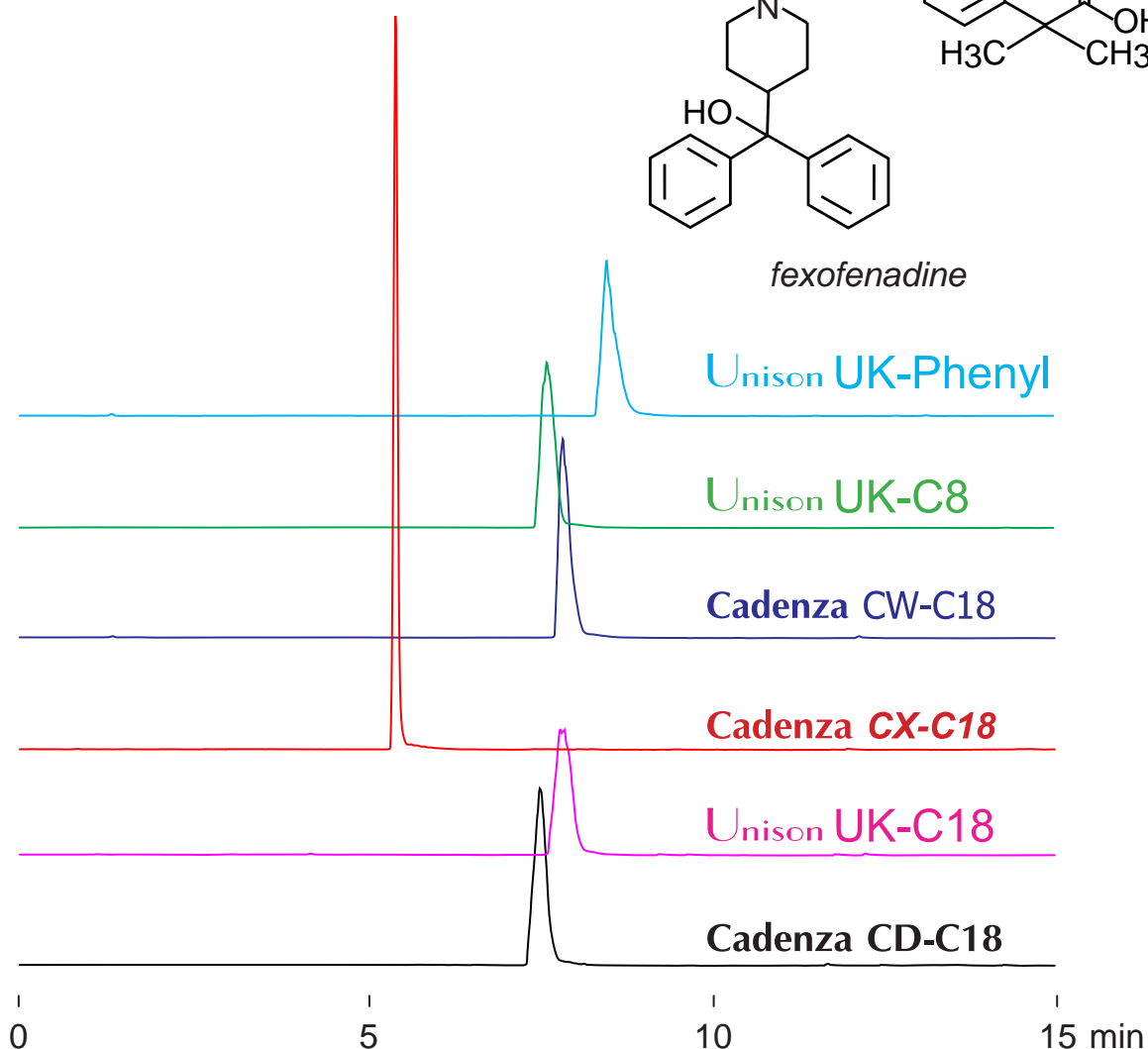

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C

ELSD (spray chamber 30°C, drift chamber 60°C)

10µL (0.25mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)