

Cadenza CW-C18

Cadenza CX-C18

Unison UK-C18

Cadenza CD-C18

Unison UK-Phenyl

Unison UK-C8


DACAPO DX-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Pararosaniline ( Dye )

パラローズアニリン ( 染料 )

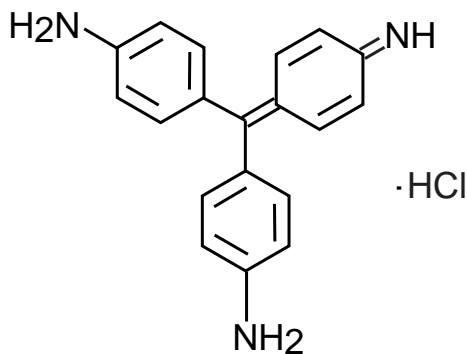

Pararosaniline

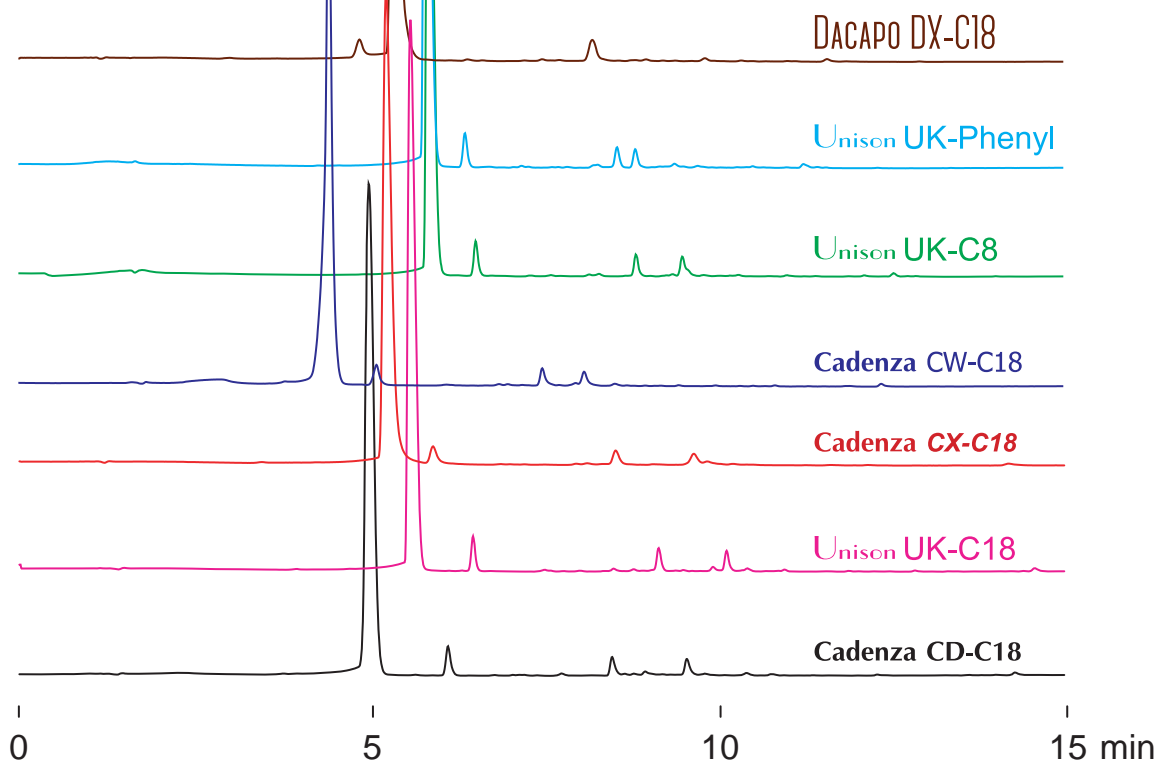

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate= 90 / 10

20-100%B (0-15min)

1.0 mL/min, 37°C, 260nm

10 µL (1mg/mL, ACN / 100mM ammonium formate = 50 / 50)