

Cadenza CX-C18 Unison UK-Phenyl


Unison UK-C18 Unison UK-C8

Cadenza CD-C18 Cadenza CW-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Fursultiamine (Vitamin B1 derivative)

フルスルチアミン (ビタミン B1 誘導体)

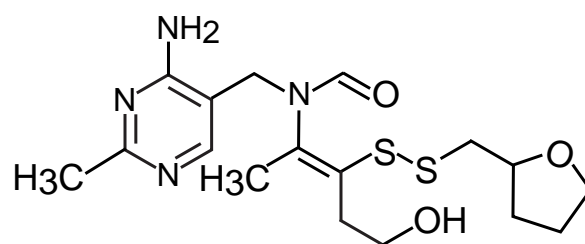

fursultiamine

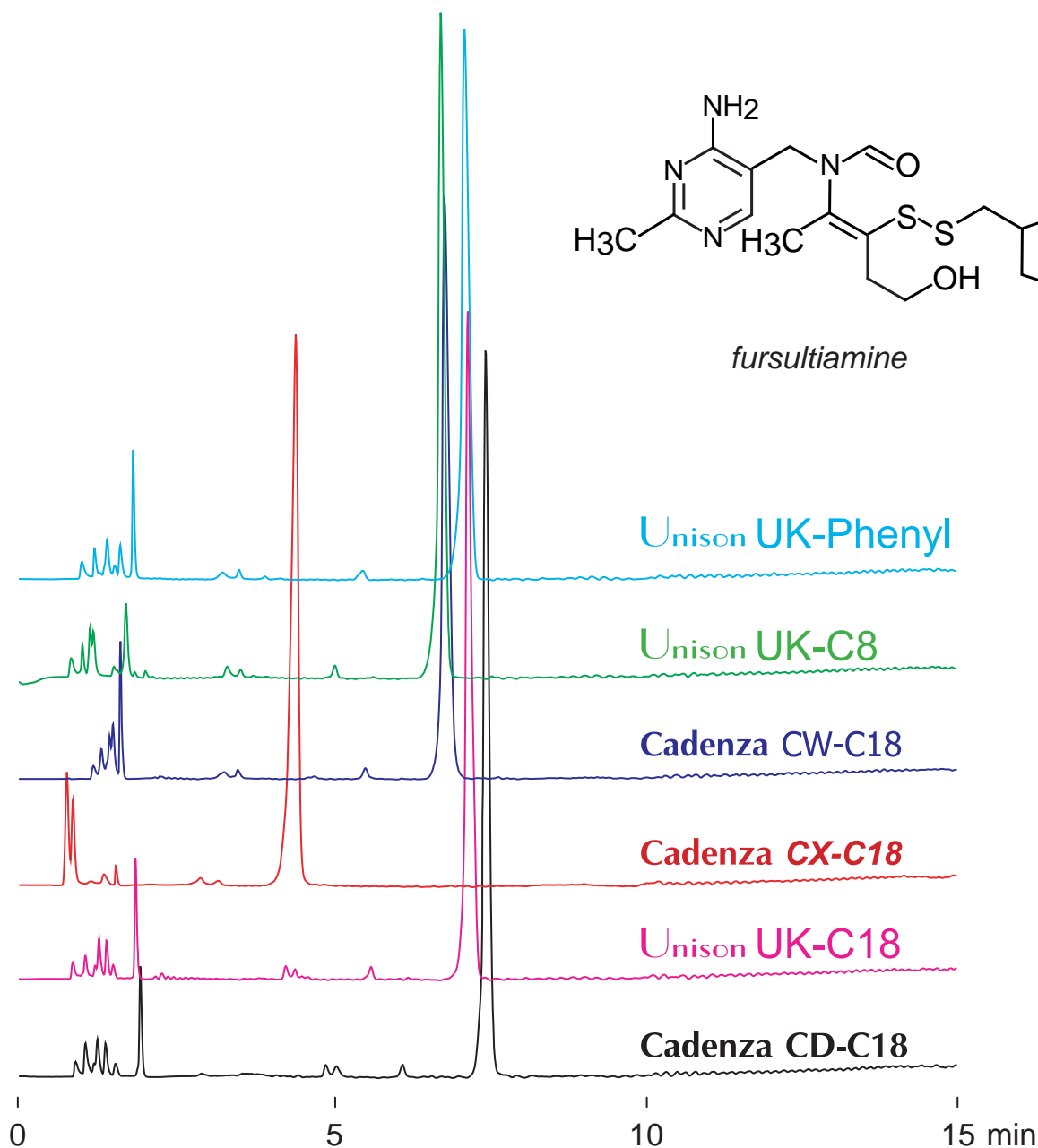

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: methanol / formic acid = 100 / 0.1

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10μL (0.5mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)