

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Ponceau 3R ( Artificial coloring )

ポンソー 3R ( 人工着色料 )

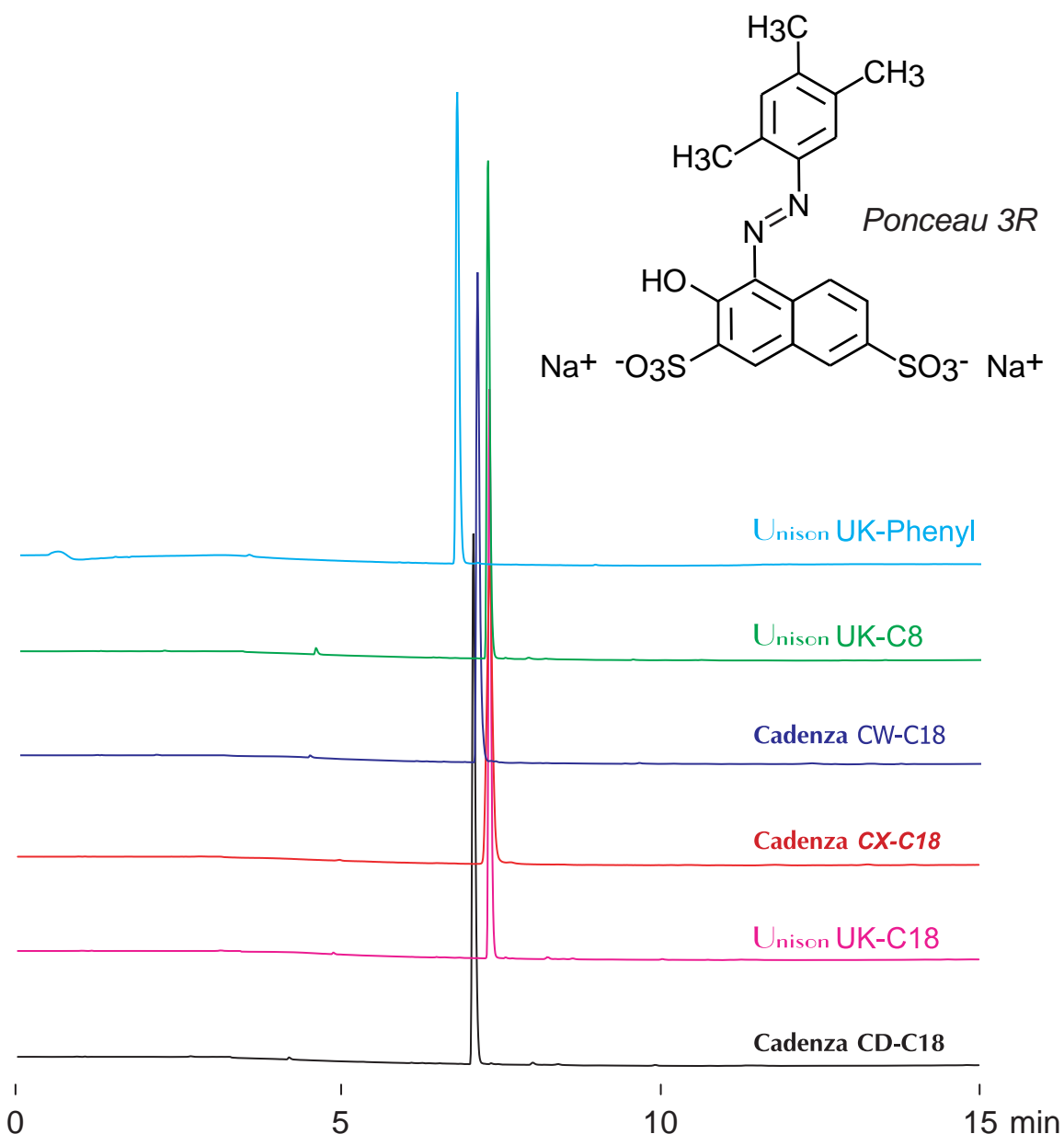

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-100%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, 100mM ammonium formate)