

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Allura Red AC (Food additive)

赤色 40号 (食品添加物)

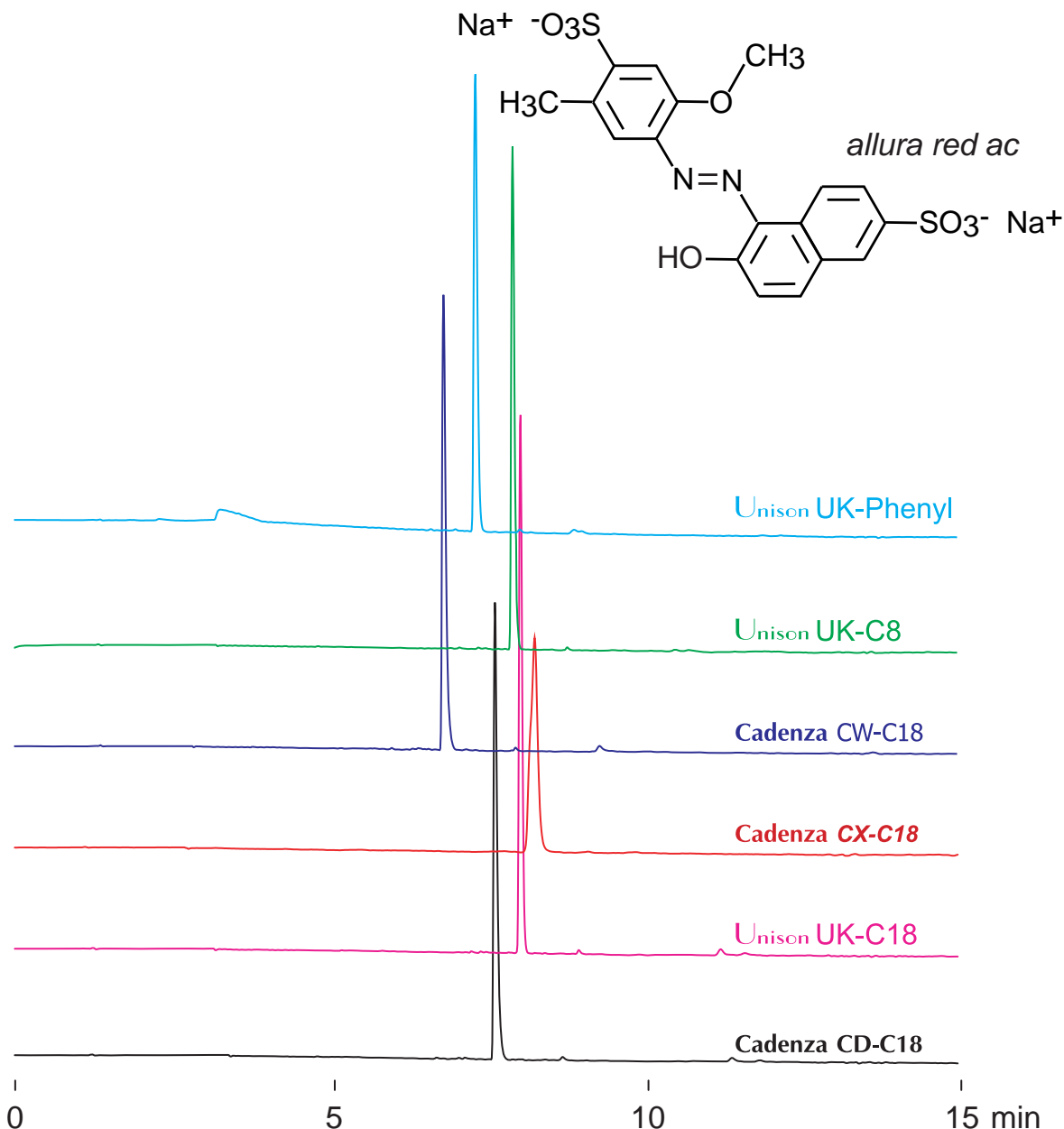

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-60%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate)