

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

## New Coccine ( Food additive )

赤色 102号 (食品添加物 )

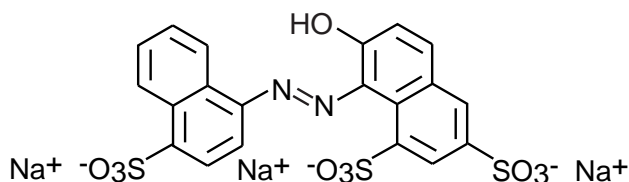

New Coccine

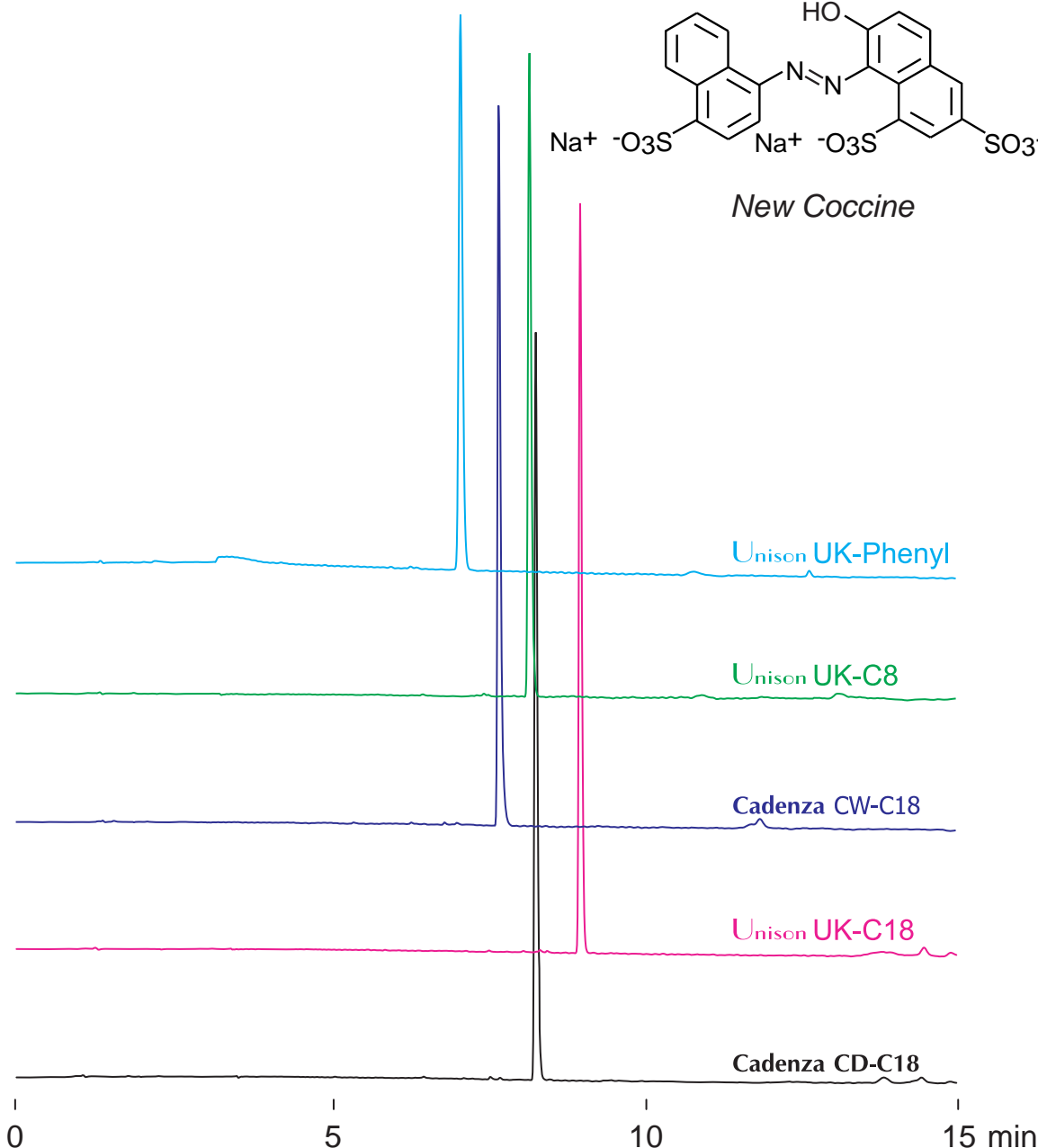

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-40%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate)