

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18



Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Rose Bengal (Food additive)

赤色 105号 (食品添加物)

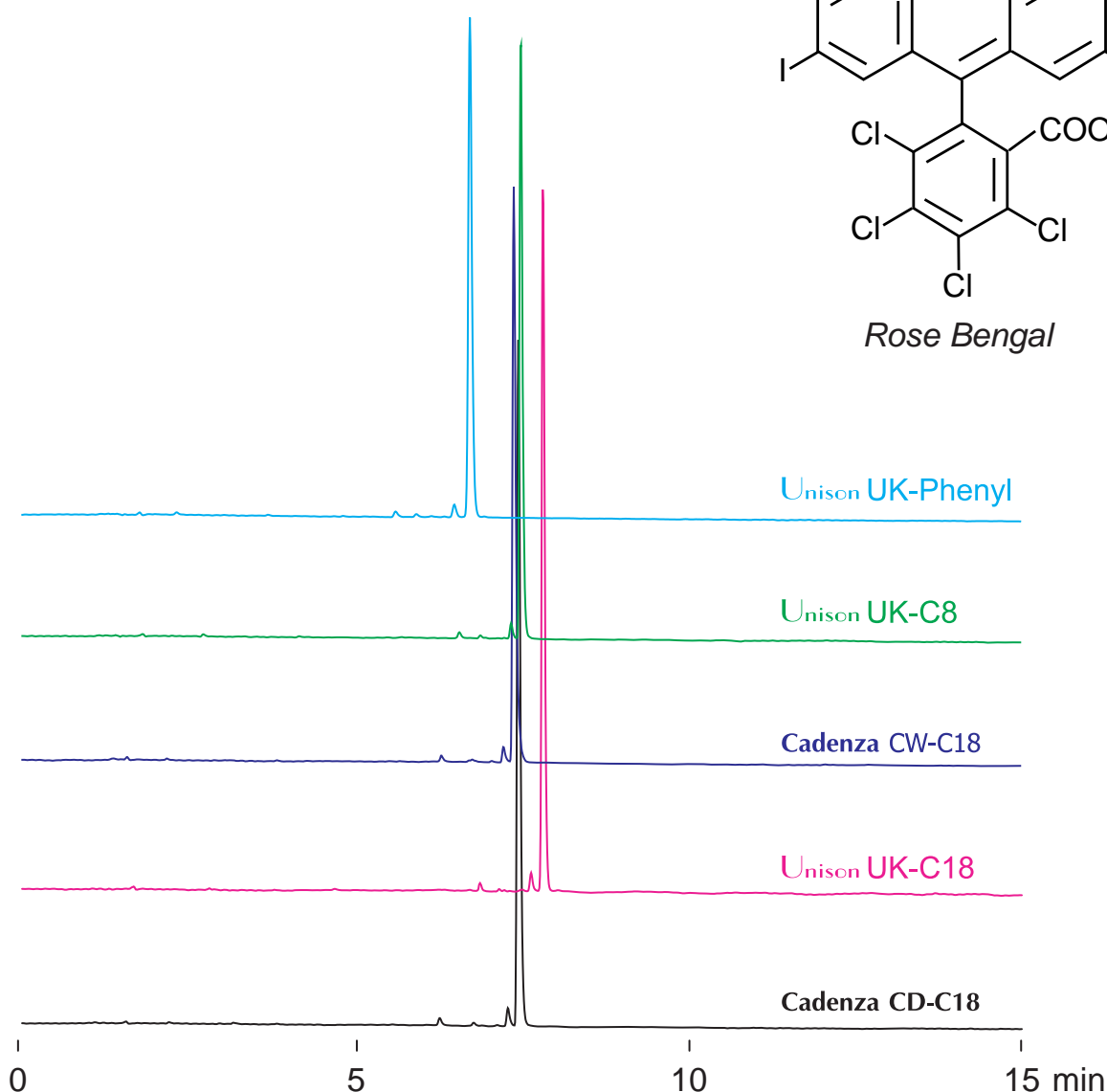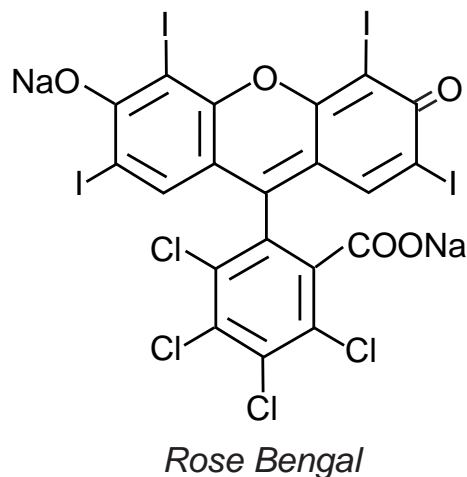

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

20-100%B (0-15min)

1mL/min (8MPa), 37°C, 260nm

10µL (0.2mg/mL, 100mM ammonium formate)