


Unison UK-Phenyl
Unison UK-C8
Unison UK-C18

Cadenza CW-C18
Cadenza CX-C18
Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Light Green SF Yellowish (Artificial coloring)

ライトグリーンSF黄口 (人工着色料)

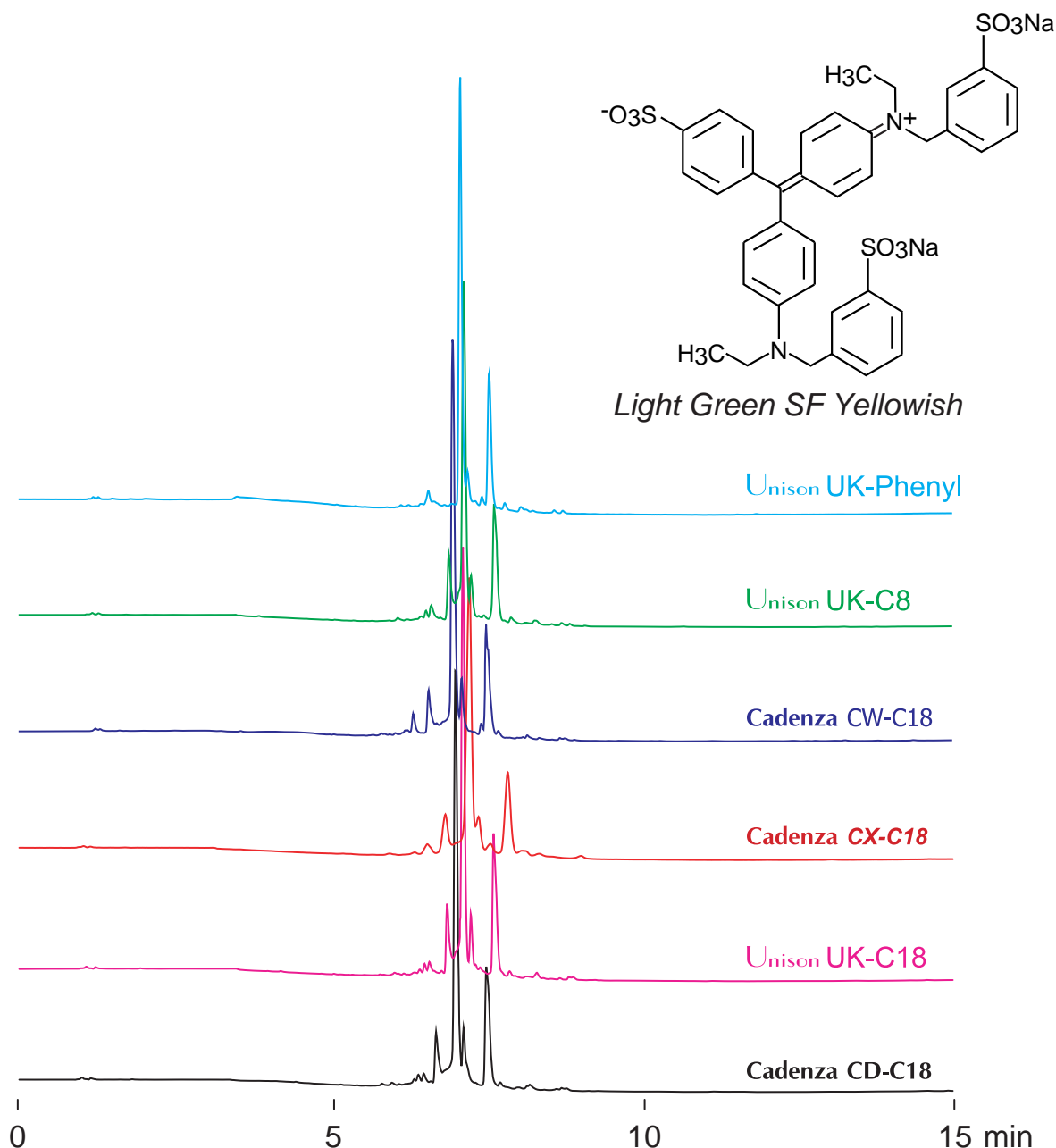

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-100%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10µL (0.5mg/mL, 100mM ammonium formate)