
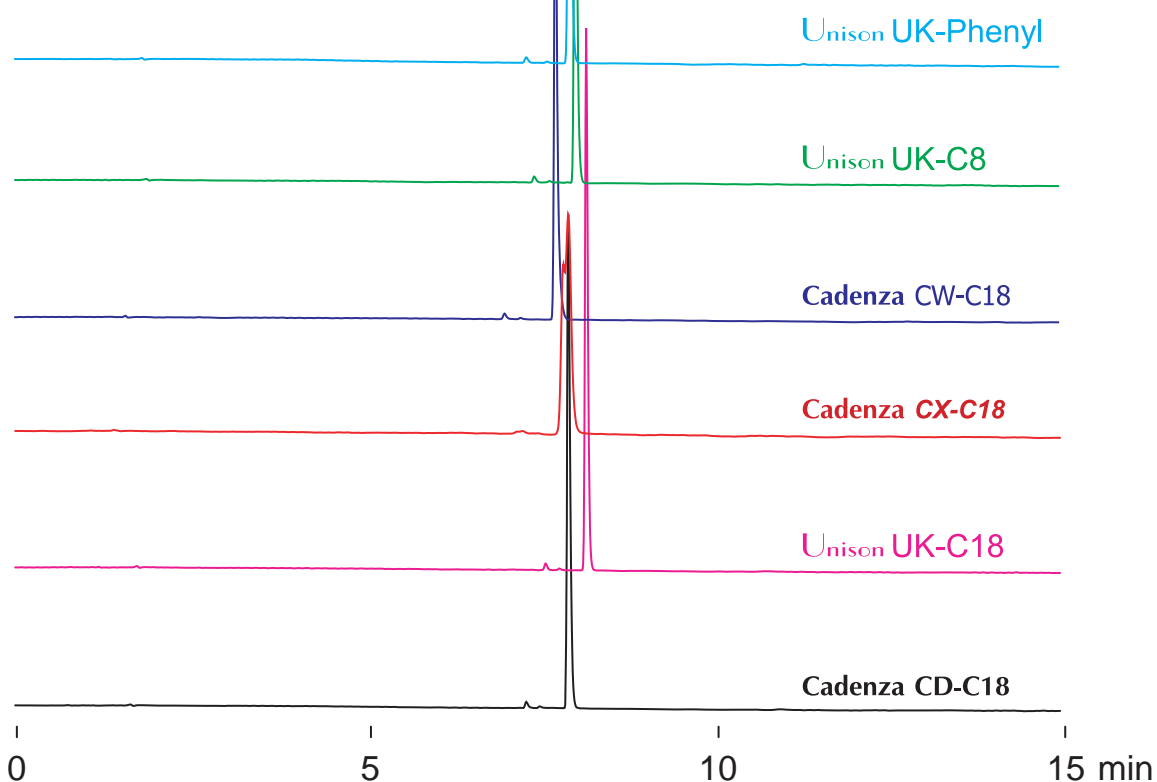


Unison UK-Phenyl Cadenza CW-C18
 Unison UK-C8 Cadenza CX-C18
 Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Acid Red (Food additive)

赤色 106号 (食品添加物)

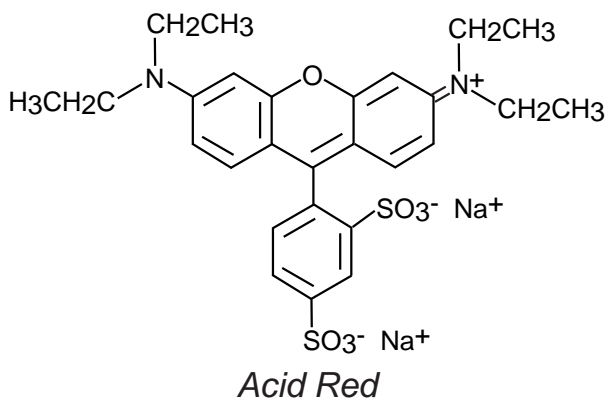

100 x 4.6 mm
 A: 100mM ammonium acetate
 B: acetonitrile / 100mM ammonium acetate = 90 / 10
 10-100%B (0-15min)
 1mL/min (8MPa), 37°C, 260nm
 10µL (0.08mg/mL, 100mM ammonium formate)