



Unison UK-Phenyl  
 Unison UK-C8  
 Unison UK-C18

Cadenza CW-C18  
 Cadenza CX-C18  
 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Acid Green 9 ( Artificial coloring )

ブリリアントグリーン (人工着色料)

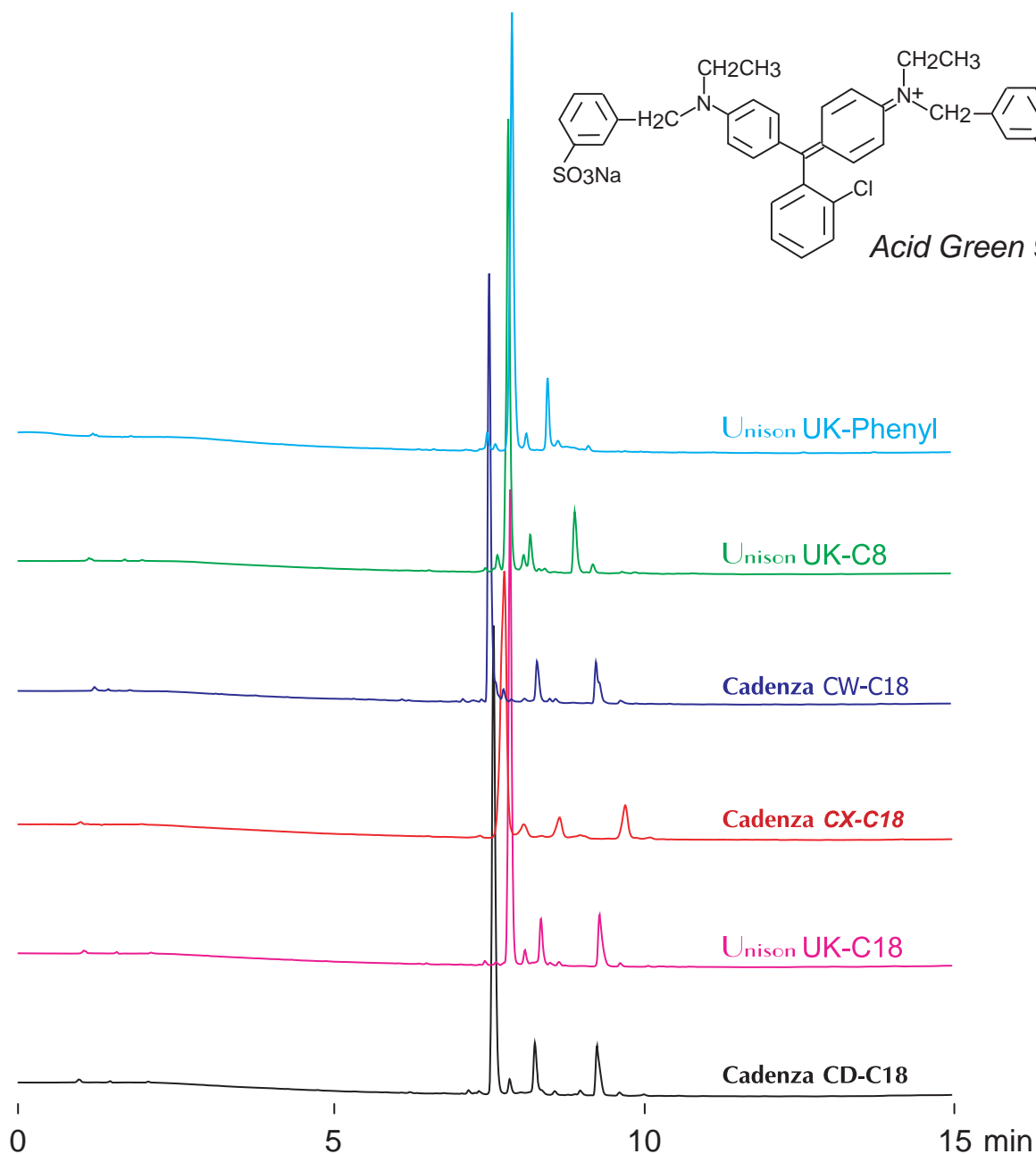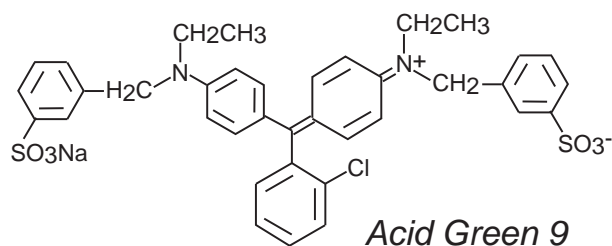

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

10-100%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, 100mM ammonium formate)