
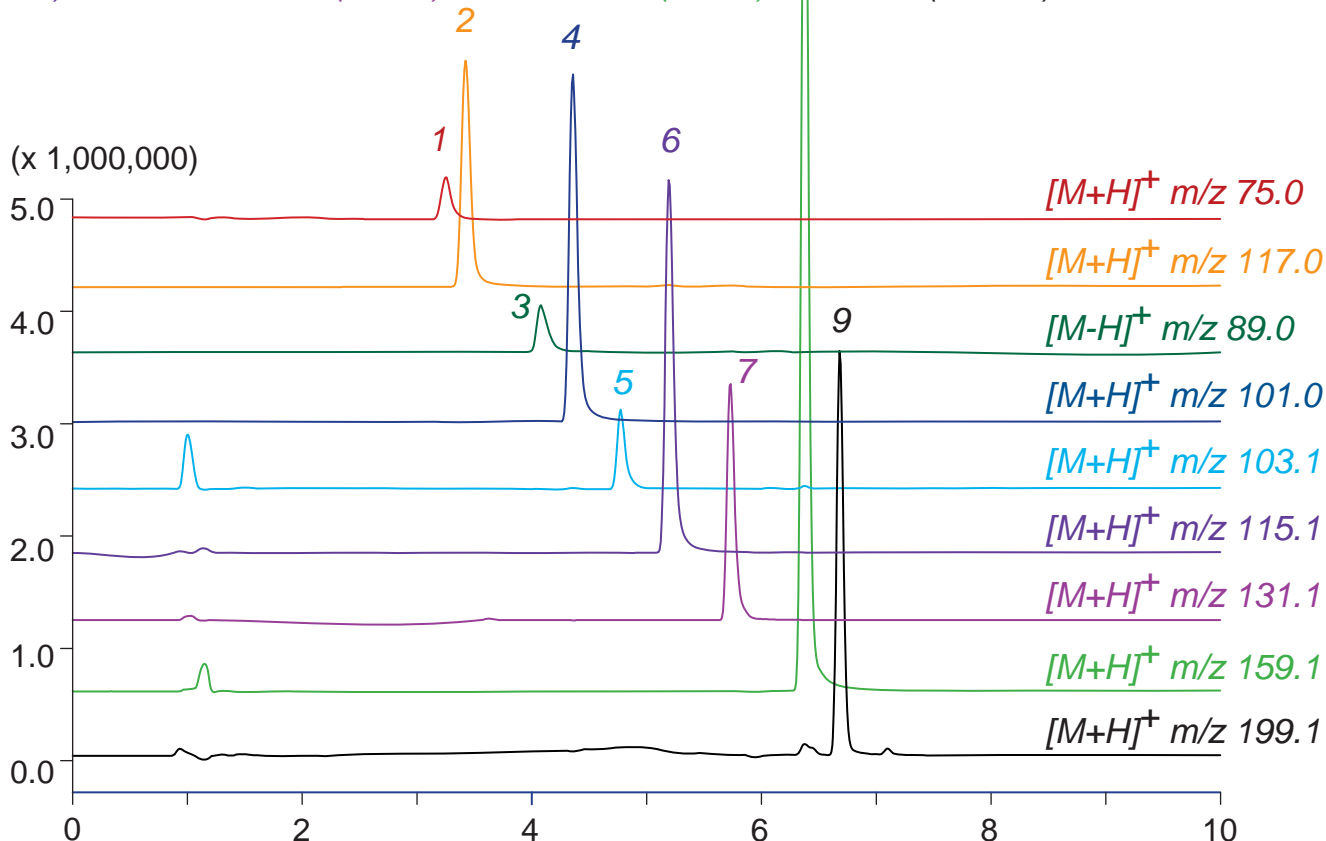


Unison UK-Phenyl MF 100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: N-nitrosoamines

N - ニトロソアミンの LC-MS 分析

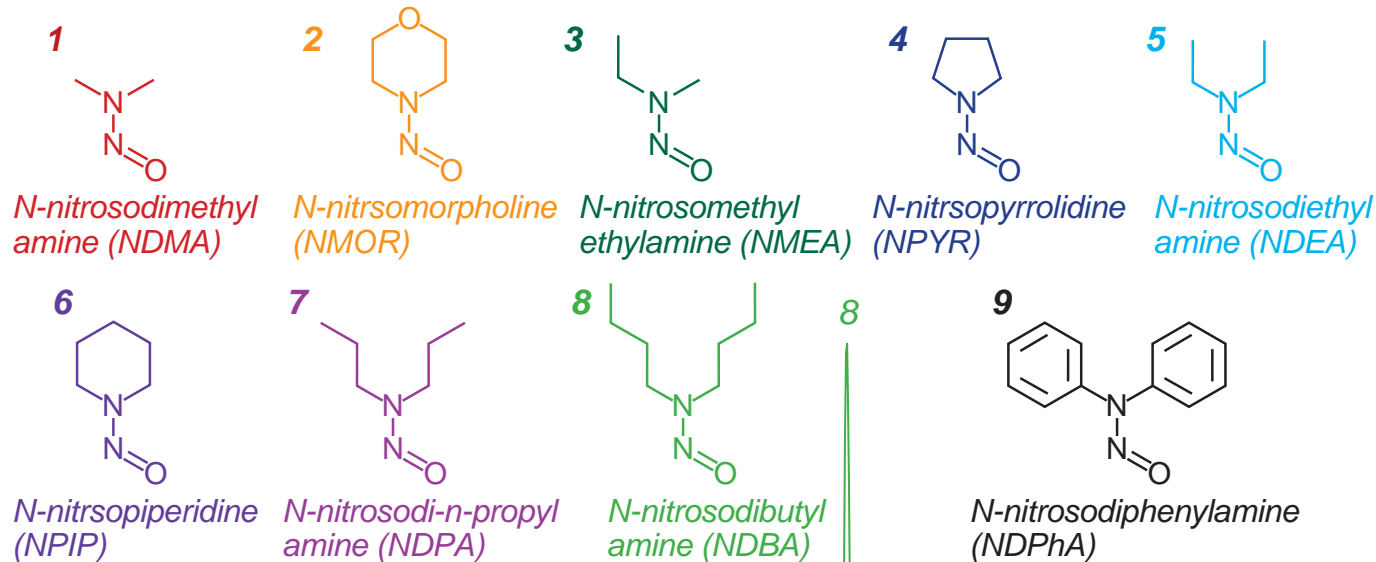

Unison UK-Phenyl MF, 100 x 2 mm
 A: water / formic acid = 100 / 0.1
 B: methanol / formic acid = 100 / 0.1
 0-90%B (0-6min), 90%B (6-10min)
 0.2mL/min (7MPa), 37°C, MS (ESI, SIM, positive)
 5µL (4µg/mL, MeOH/water/HCOOH = 50/50/0.1)