

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## Imatinib ( Antineoplastic )

イマチニブ (抗悪性腫瘍薬 )

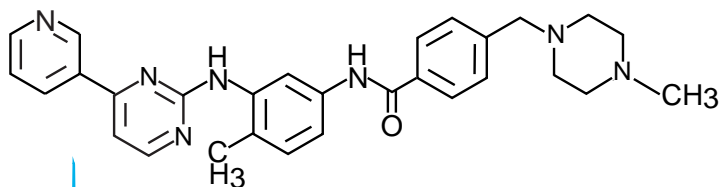

*Imatinib*

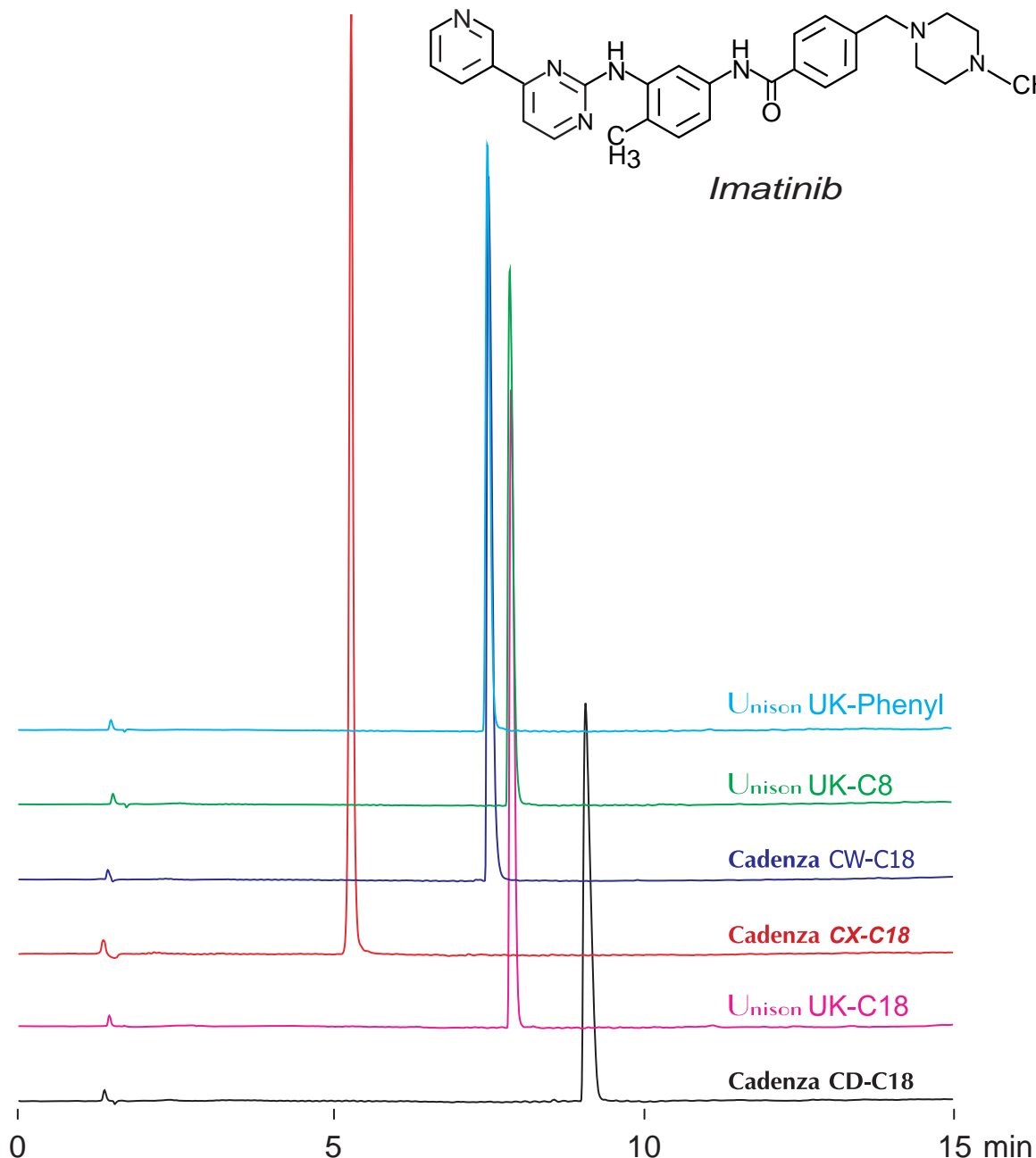

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-85%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (0.05mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)