

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Imatinib (Antineoplastic)

イマチニブ (抗悪性腫瘍薬)

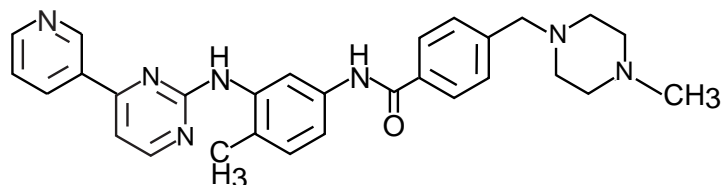

Imatinib

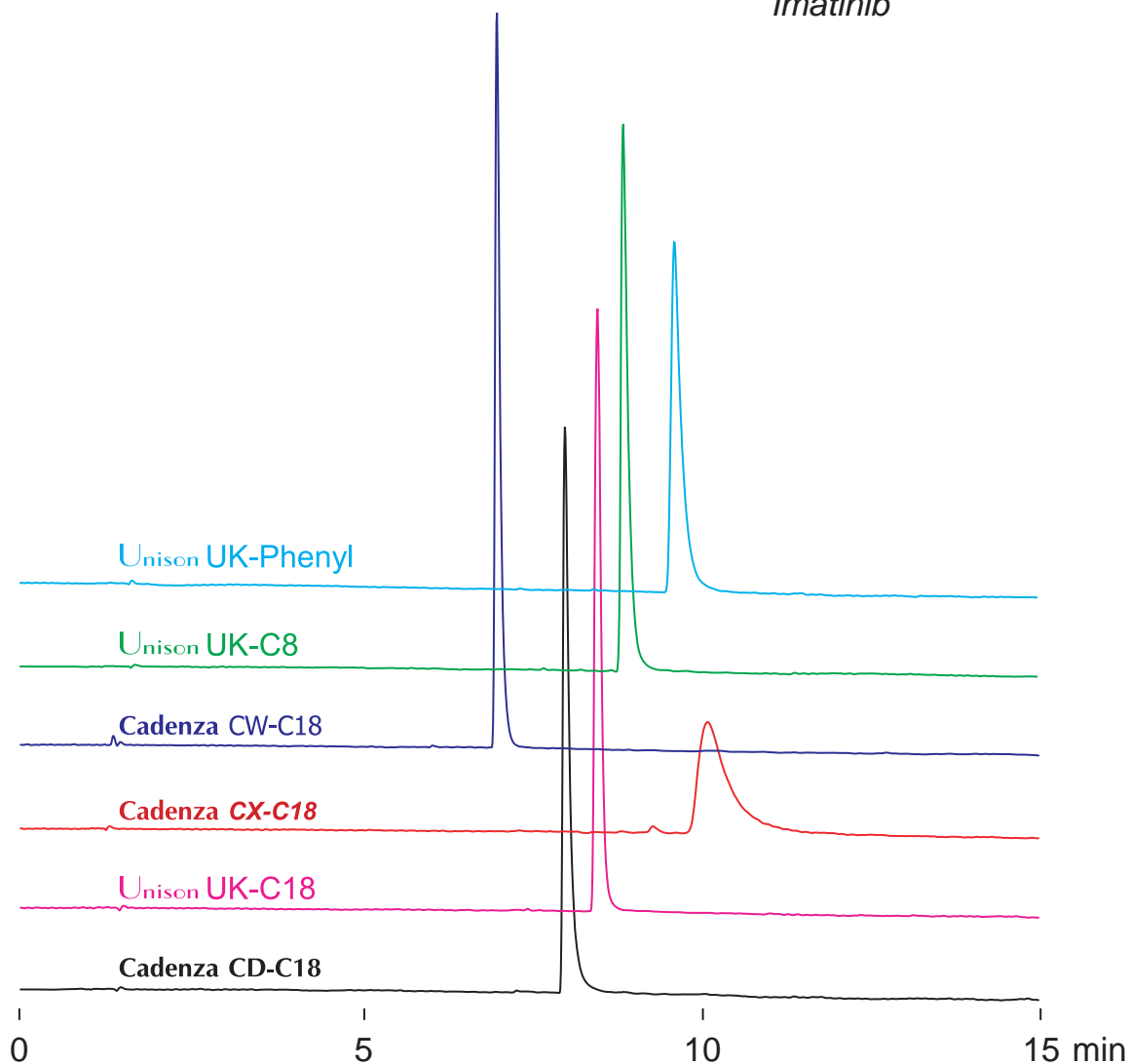

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

15-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (0.05mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)