


Unison UK-C18
Unison UK-Phenyl

100 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Acephate

アセフェートの LC MS 分析

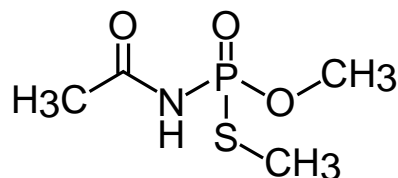

acephate
[M+H]⁺ m/z 184.0

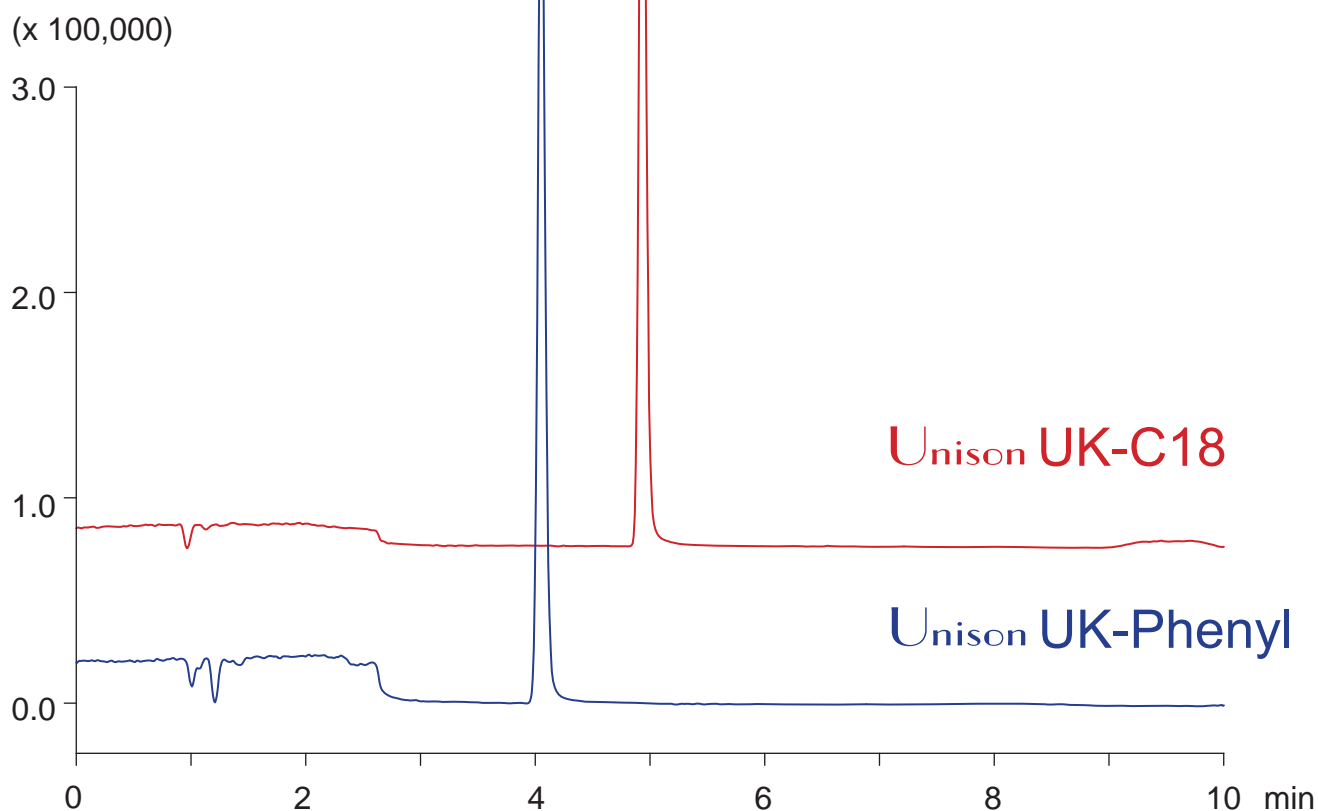

100 x 3 mm
A: water / formic acid = 100 / 0.1
B: acetonitrile / water / formic acid = 20 / 80 / 0.1
0-100%B (0-10min)
0.4mL/min (6MPa), 37°C, MS (ESI, SIM, positive)
5µL (1µg/mL, water / HCOOH = 100 / 0.1)