

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Edoxaban (Antithrombotic)

エドキサバン (抗血栓薬)

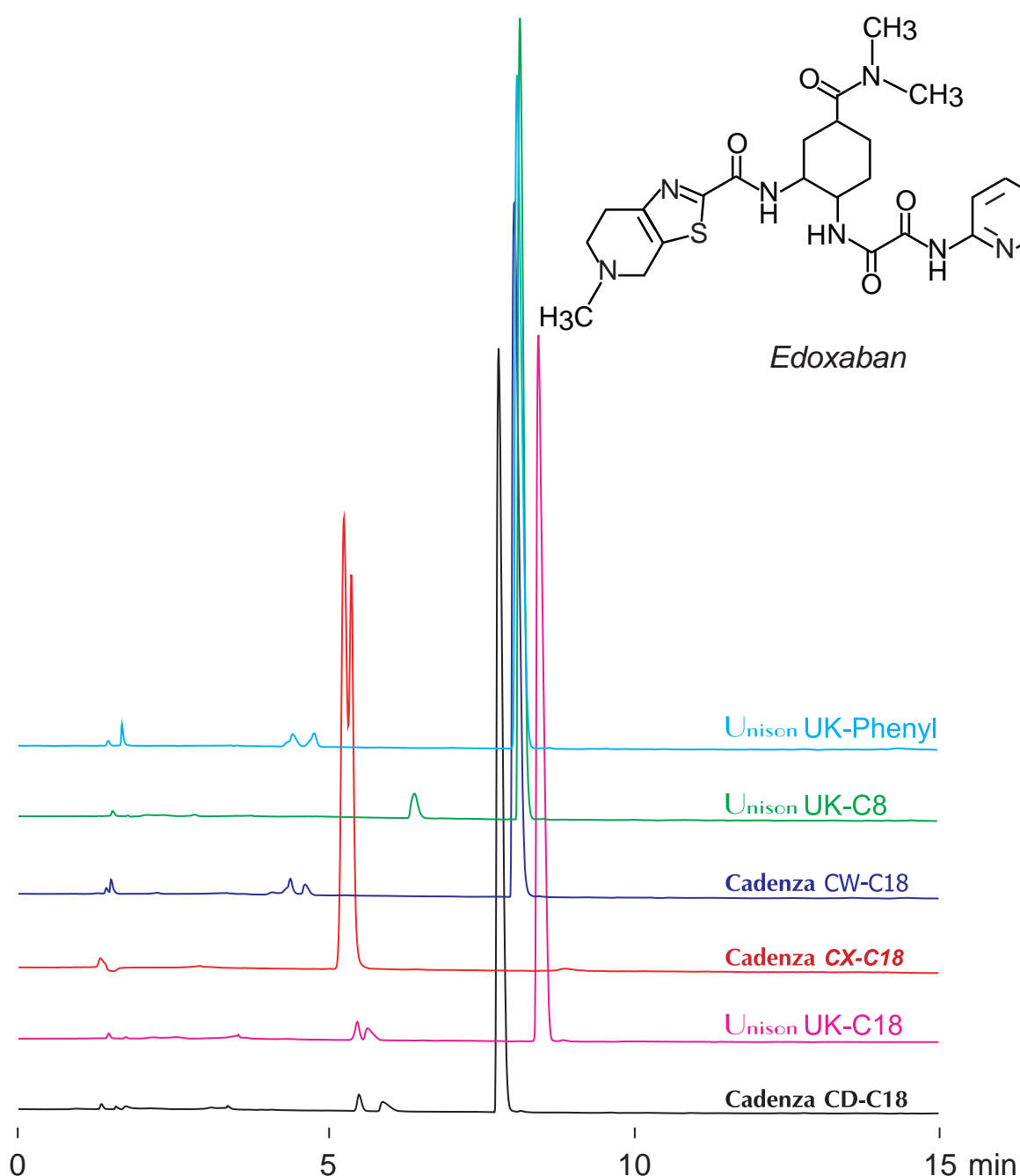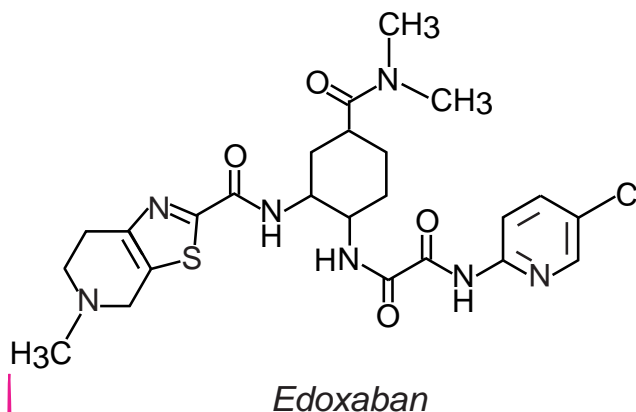

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-60%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)