


Cadenza CD-C18

250 x 4.6 mm



HPLC COLUMN APPLICATION

Polysorbate 20 (Nonionic surfactant)

ポリソルベート20 (非イオン性界面活性剤)

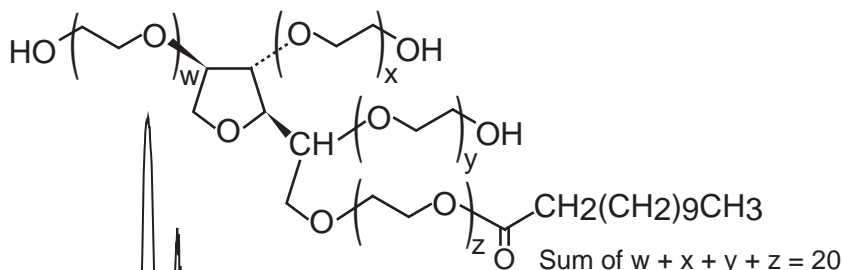

Polysorbate 20
(Tween20)

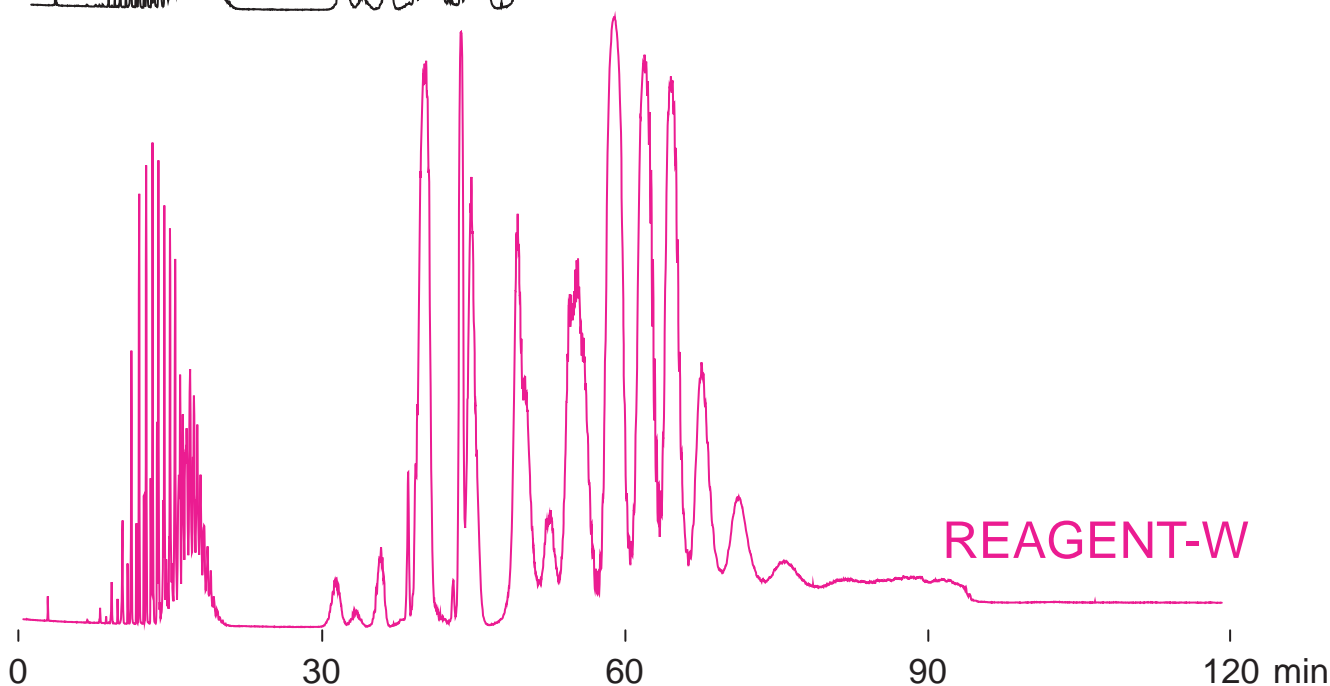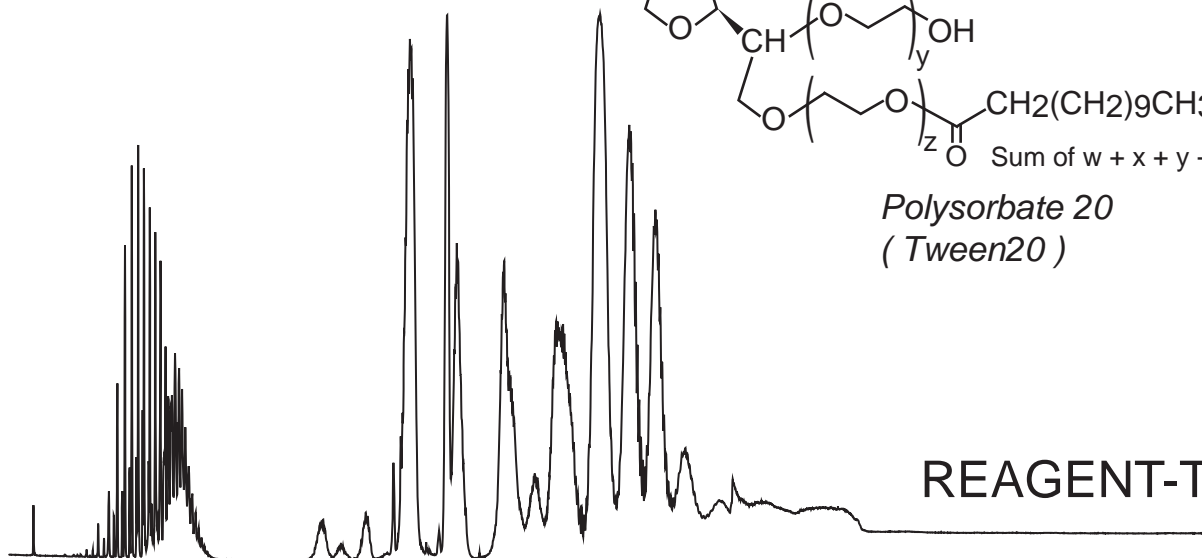

Cadenza CD-C18, 250 x 4.6 mm
 A: water / acetic acid = 100 / 0.1
 B: acetonitrile / acetic acid = 100 / 0.1
 13-100%B (0-60min), 100%B (60-120min)
 1.0mL/min (14MPa), 60°C
 ELSD (spray chamber 70°C, drift chamber 120°C)
 20μL (50mg/mL, 100mM ammonium formate)