

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8
Cadenza CX-C18


Unison UK-C18
Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Fluvastatin (Hypolipidemic agent)

フルバスタチン (高脂血症治療薬)

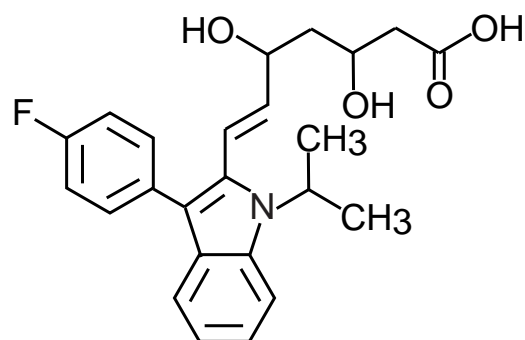

fluvastatin

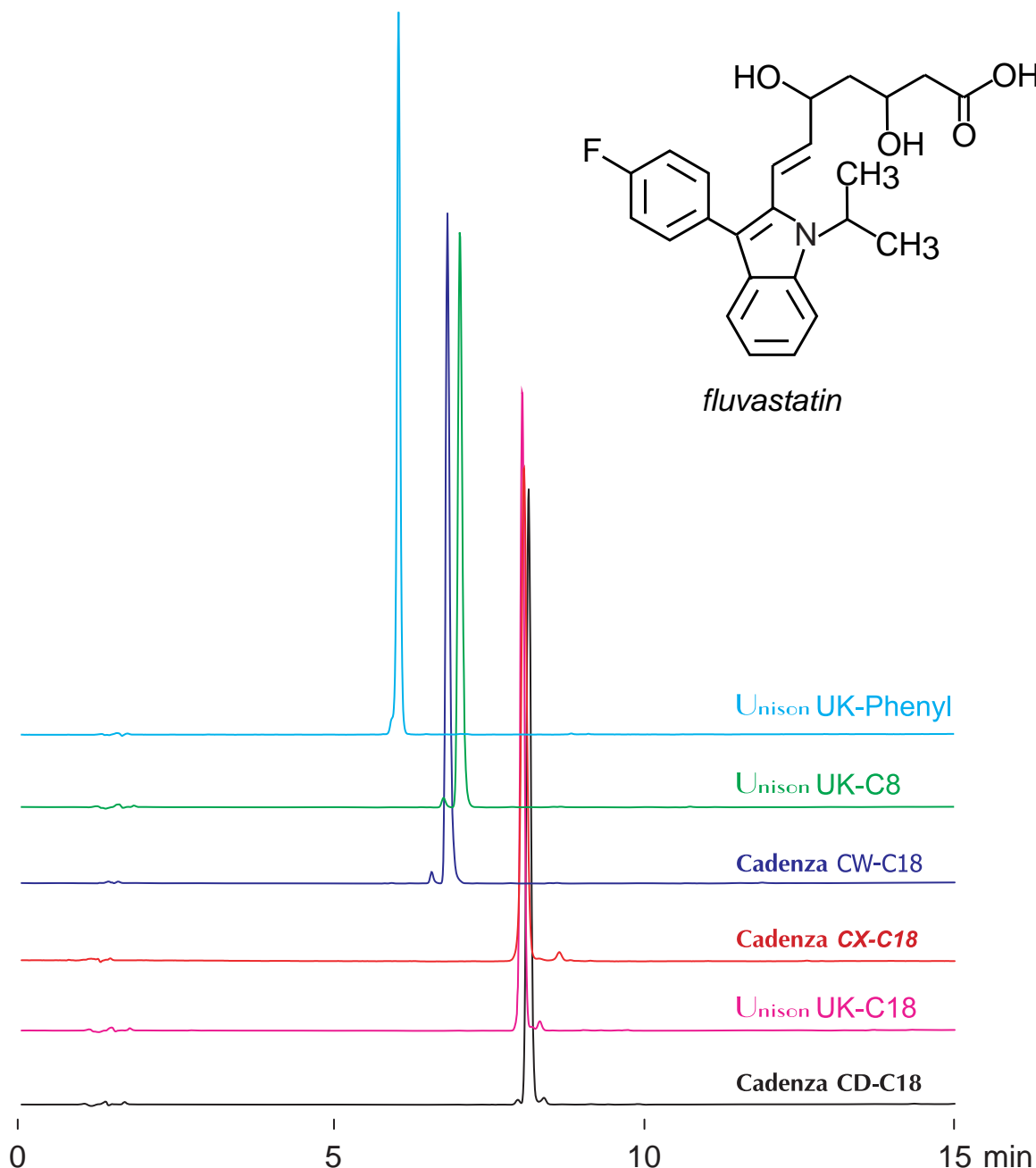

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

30-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10μL (0.25mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)