



Unison UK-C8
Unison UK-C18

250 x 4.6 mm


HPLC COLUMN APPLICATION

Polysorbate 20 (Nonionic surfactant)

ポリソルベート20 (非イオン性界面活性剤)

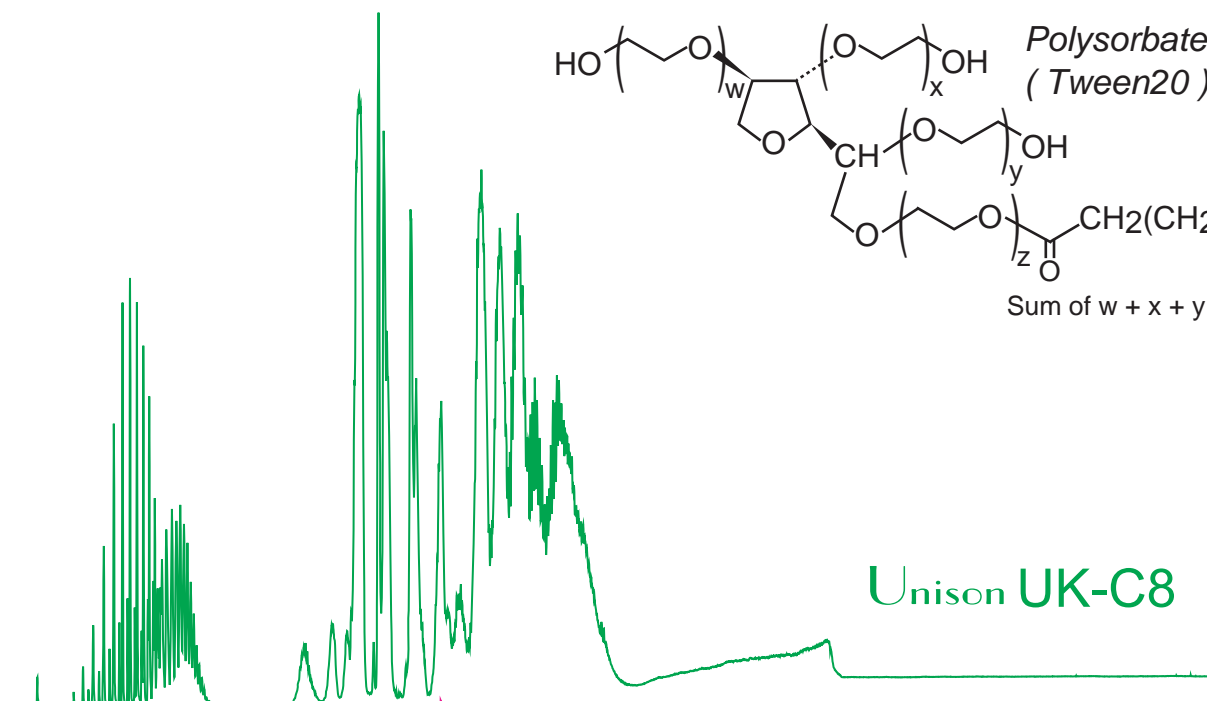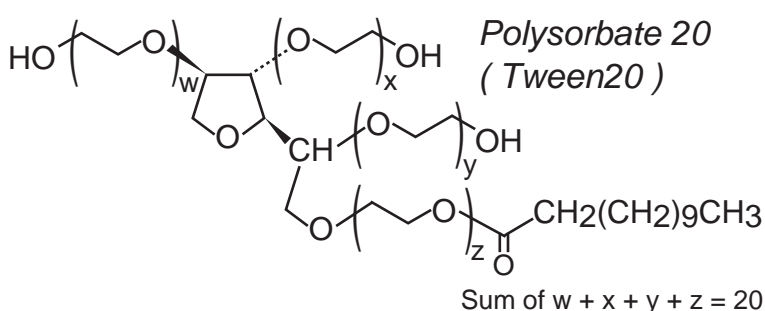

Unison UK-C8

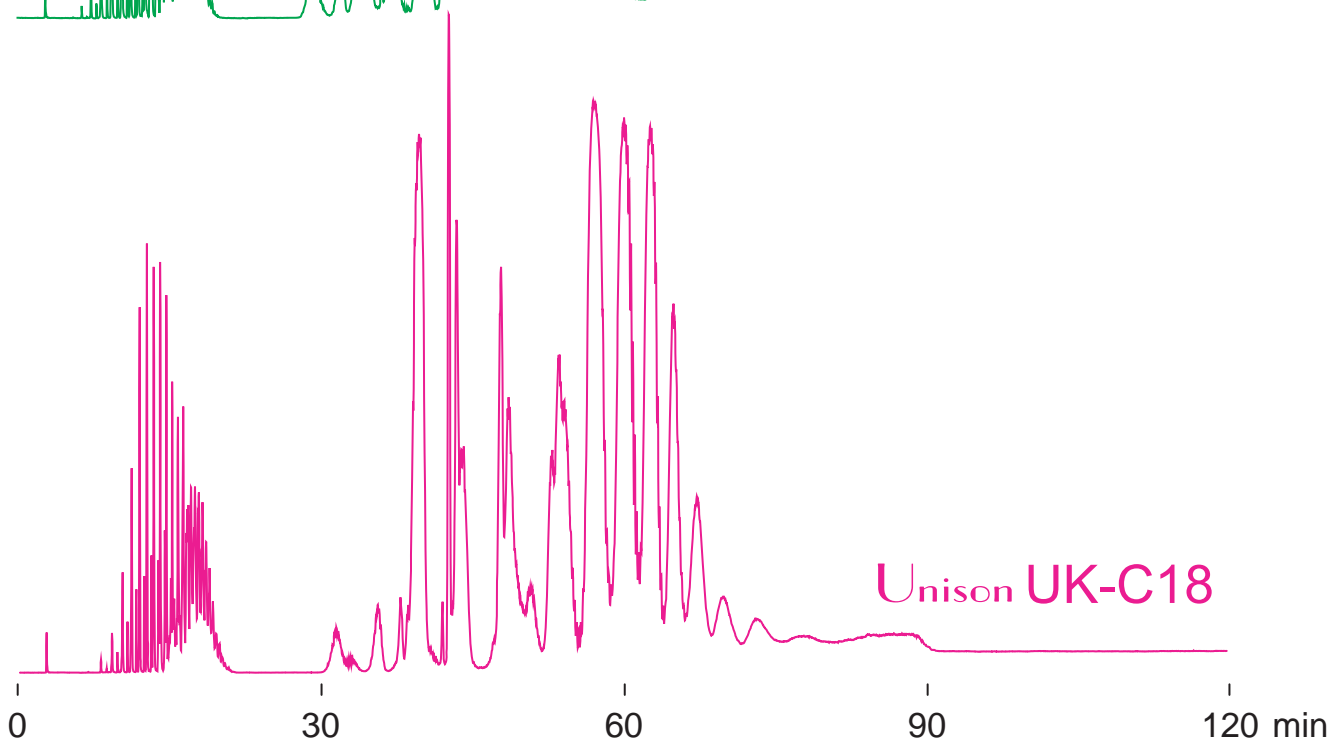

Unison UK-C18

250 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

13-100%B(0-60min), 100%B(60-120min),

1.0mL/min(14MPa), 60°C, ELSD(spray chamber 70°C, drift chamber 120°C)

20µL (50mg/mL, 100mM ammonium formate)