


**Intrada Organic Acid**

150 x 2 mm


**HPLC COLUMN APPLICATION**

**LC-MS: Bialaphos, Glyphosate**


ビアラホス, グリホサートの LC-MS 分析

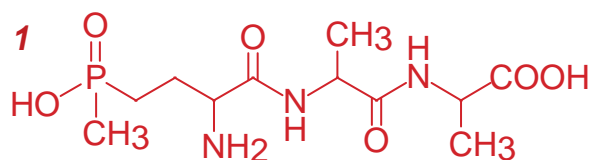

*bialaphos*

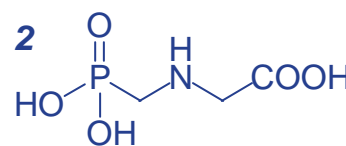

*glyphosate*

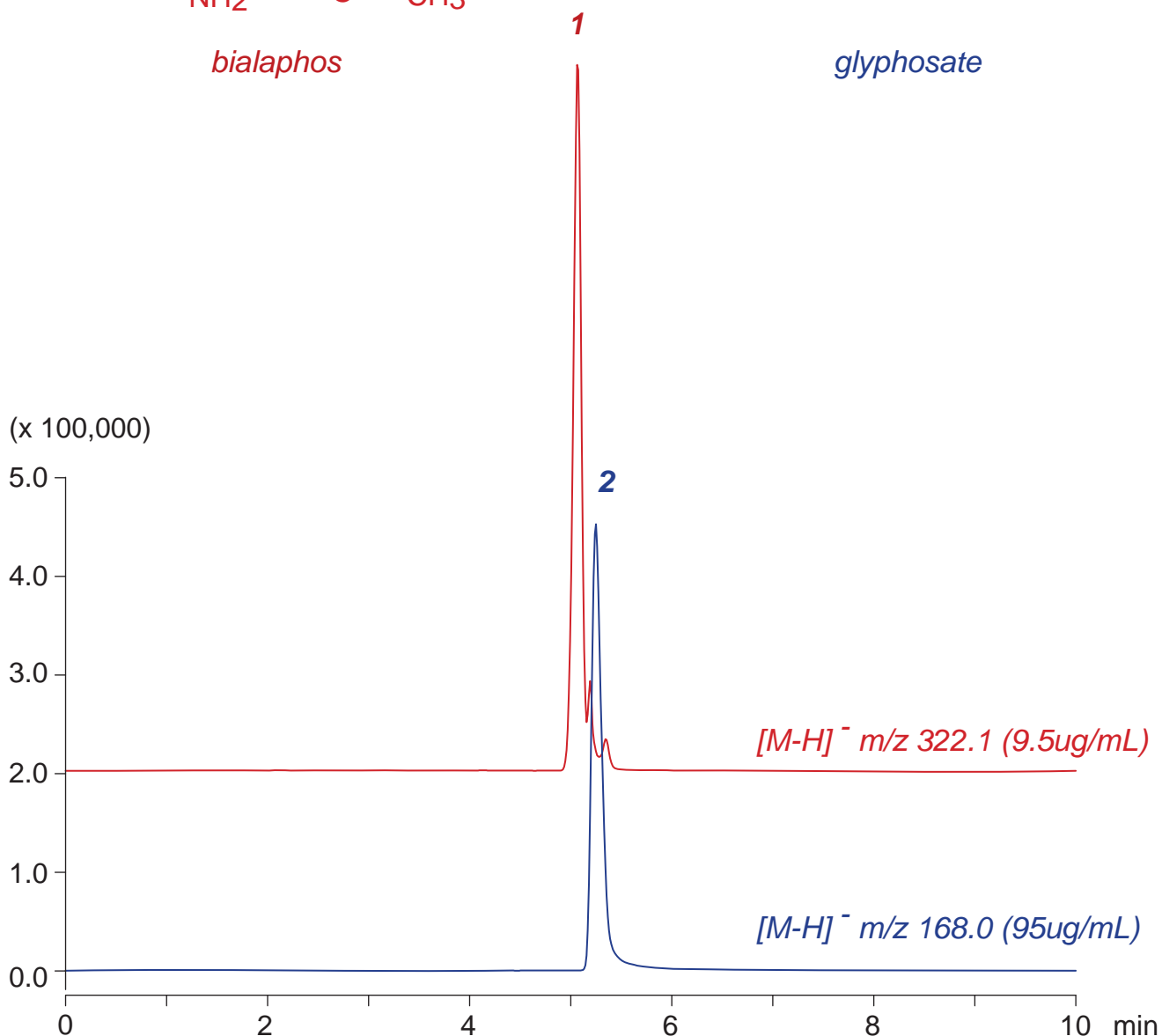

Intrada Organic Acid, 150 x 2 mm  
 A: acetonitrile / 100mM ammonium formate = 95 / 5  
 B: acetonitrile / 100mM ammonium formate = 10 / 90  
 10-100%B (0-7min), 100%B (7-10min)  
 0.2mL/min (6MPa), 37°C, MS (ESI, SIM, negative)  
 5µL (9.5-95µg/mL, water / HCOOH = 100 / 2)