



Cadenza CW-C18 Unison UK-C1  
 Cadenza CX-C18 Unison UK-Phenyl  
 Cadenza CD-C18 Unison UK-C8  
 Unison UK-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Paclitaxel ( Antineoplastic )

パクリタキセル ( 抗悪性腫瘍薬 )

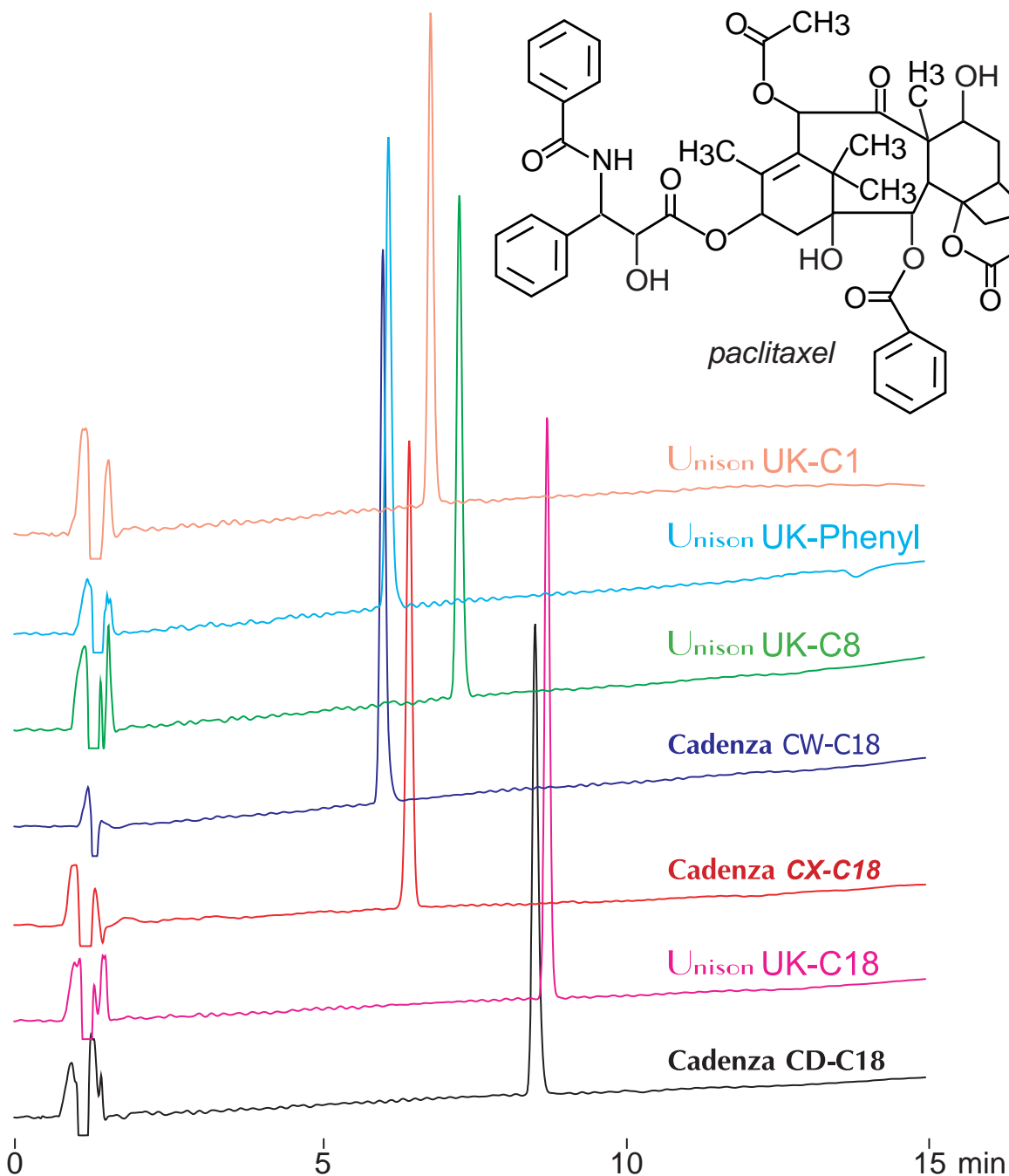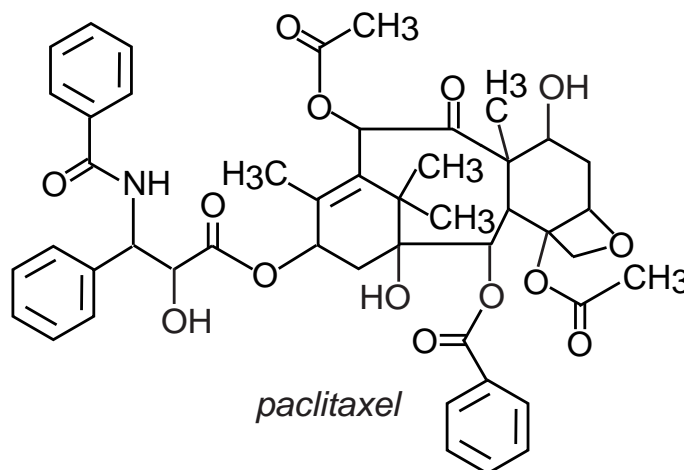

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

40-70%B (0-15min)

1.0mL/min (9MPa), 37°C, 230nm

10µL (0.1mg/mL, methanol)