









Unison UK-Phenyl

100 x 3 mm

HPLC COLUMN APPLICATION

## LC-MS: Neonicotinoid pesticides

ネオニコチノイド系殺虫剤の LC-MS 分析

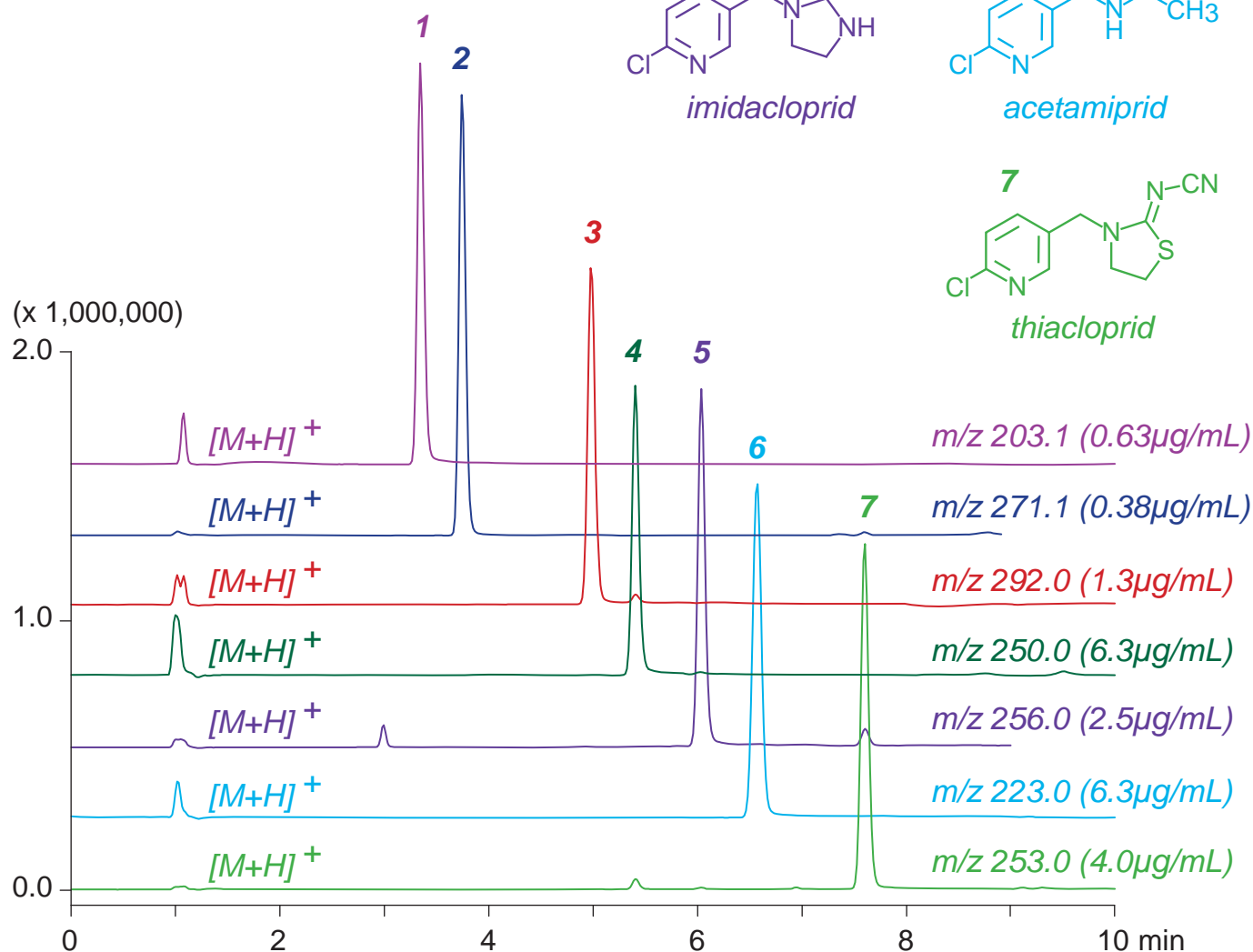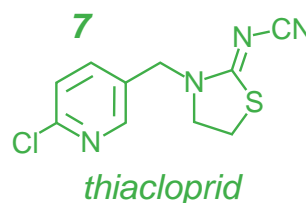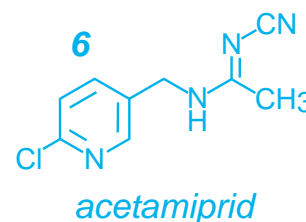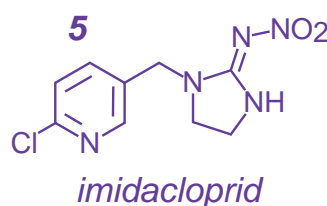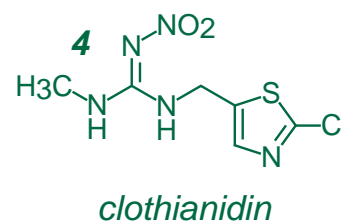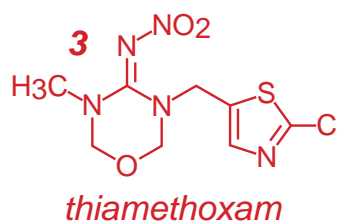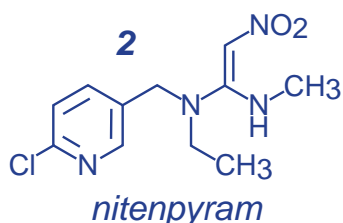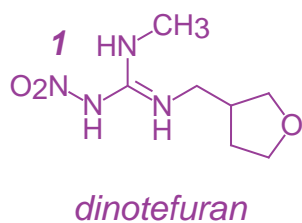

Unison UK-Phenyl, 100 x 3 mm  
 A: water / formic acid = 100 / 0.1  
 B: acetonitrile / formic acid = 100 / 0.1  
 5-50%B (0-10min)  
 0.4mL/min (6MPa), 37°C, Single Quad. MS (ESI)  
 5µL (0.38-6.3µg/mL, 0.1% formic acid)