



Cadenza CD-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Flavonoids in grapefruit juice

グレープフルーツ果汁中のフラボノイド

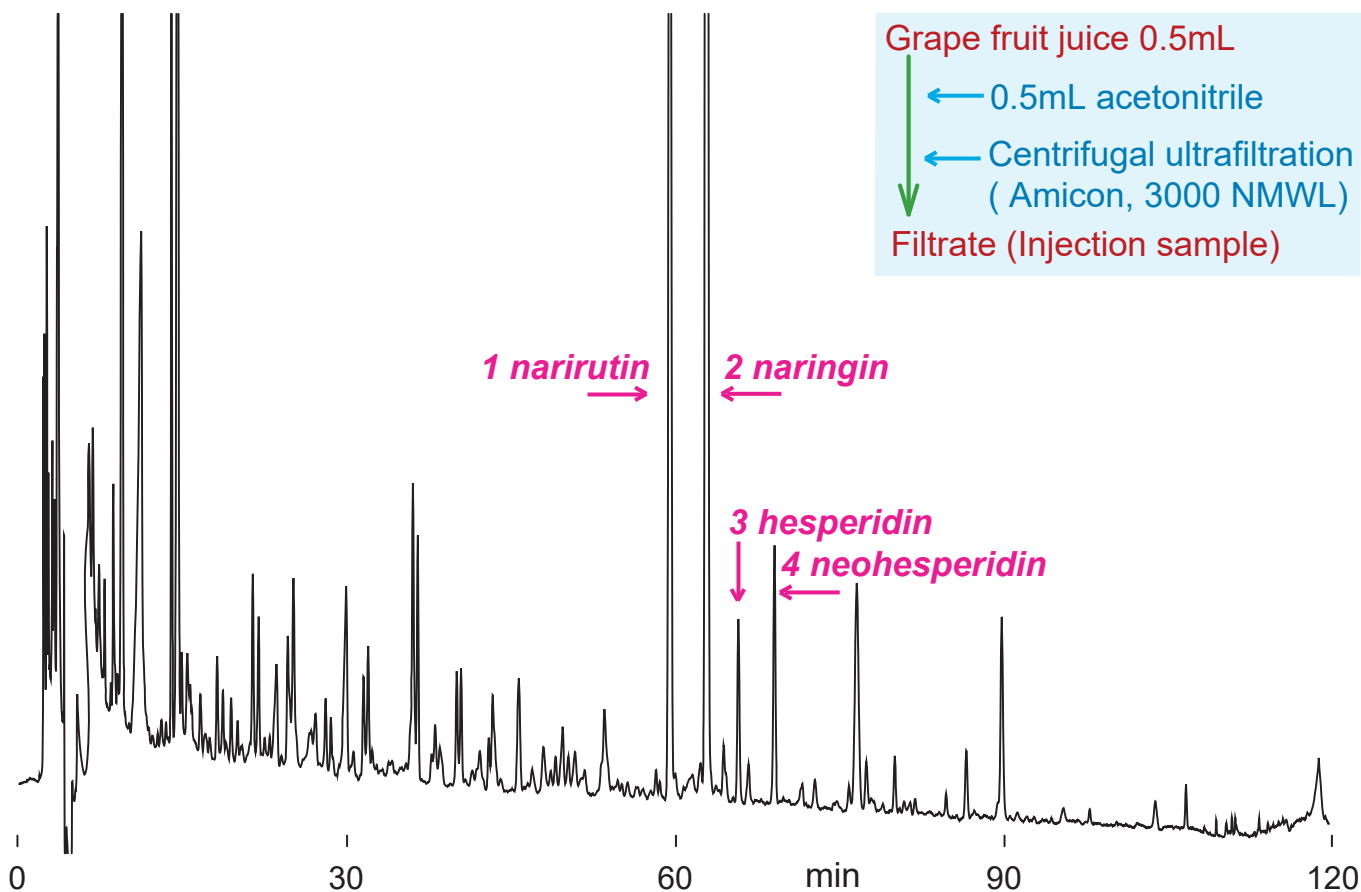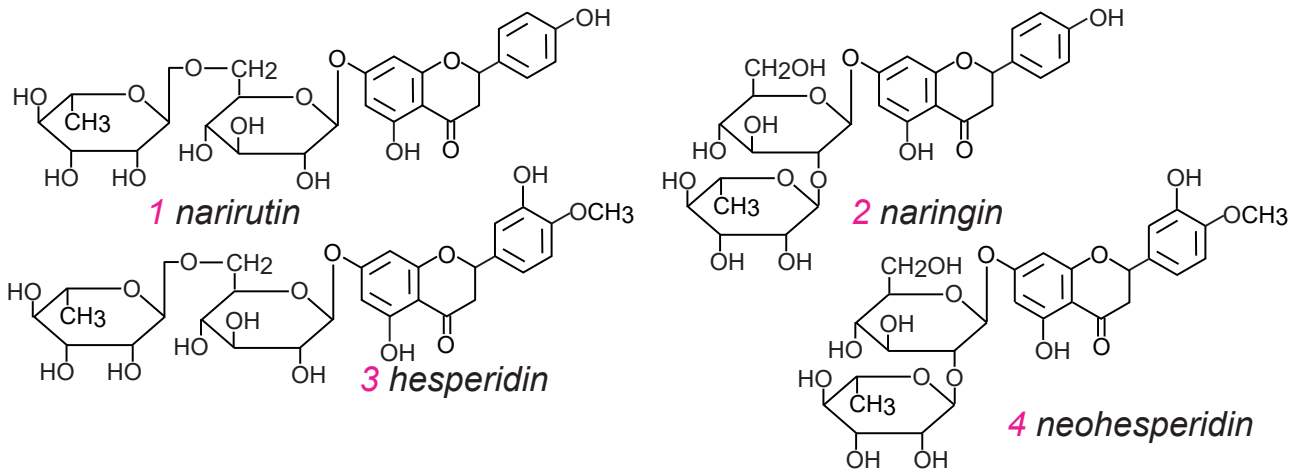

Cadenza CD-C18, 250 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

5-30%B (0-100min), 30-90%B (100-120min)

1.0mL/min (14MPa), 37°C, 280nm

10µL (water / acetonitrile = 50 / 50)