


Cadenza CD-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Organophosphorus pesticides (1)

有機リン系殺虫剤 (1) の LC-MS 分析

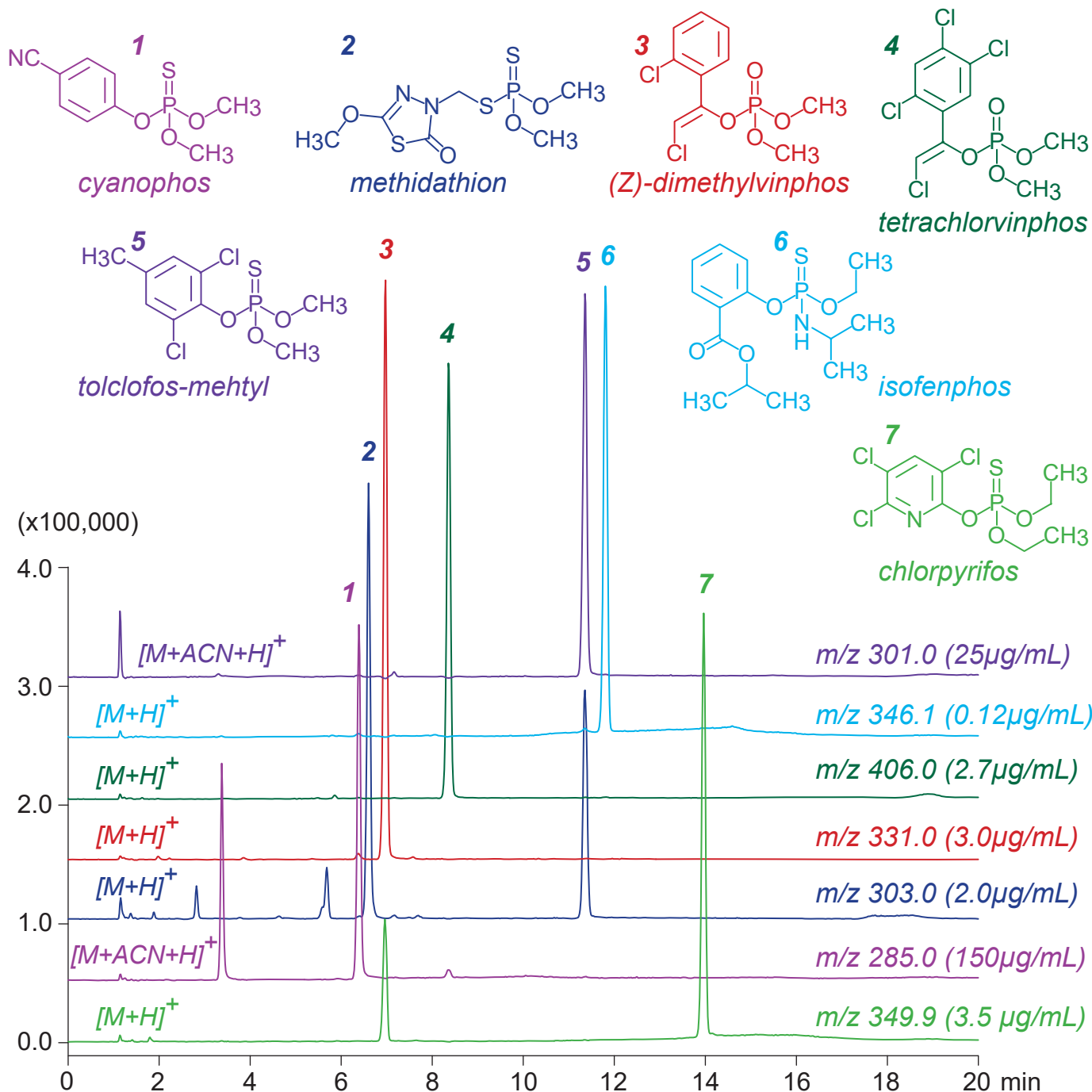

Cadenza CD-C18, 150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

50-95%B (0-15min), 95%B (15-20min)

0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)

5µL (0.12-150µg/mL, ACN / water / HCOOH = 50 / 50 / 0.1)