






Unison UK-Phenyl

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Organophosphorus pesticides (1)

有機リン系殺虫剤 (1) の LC-MS 分析

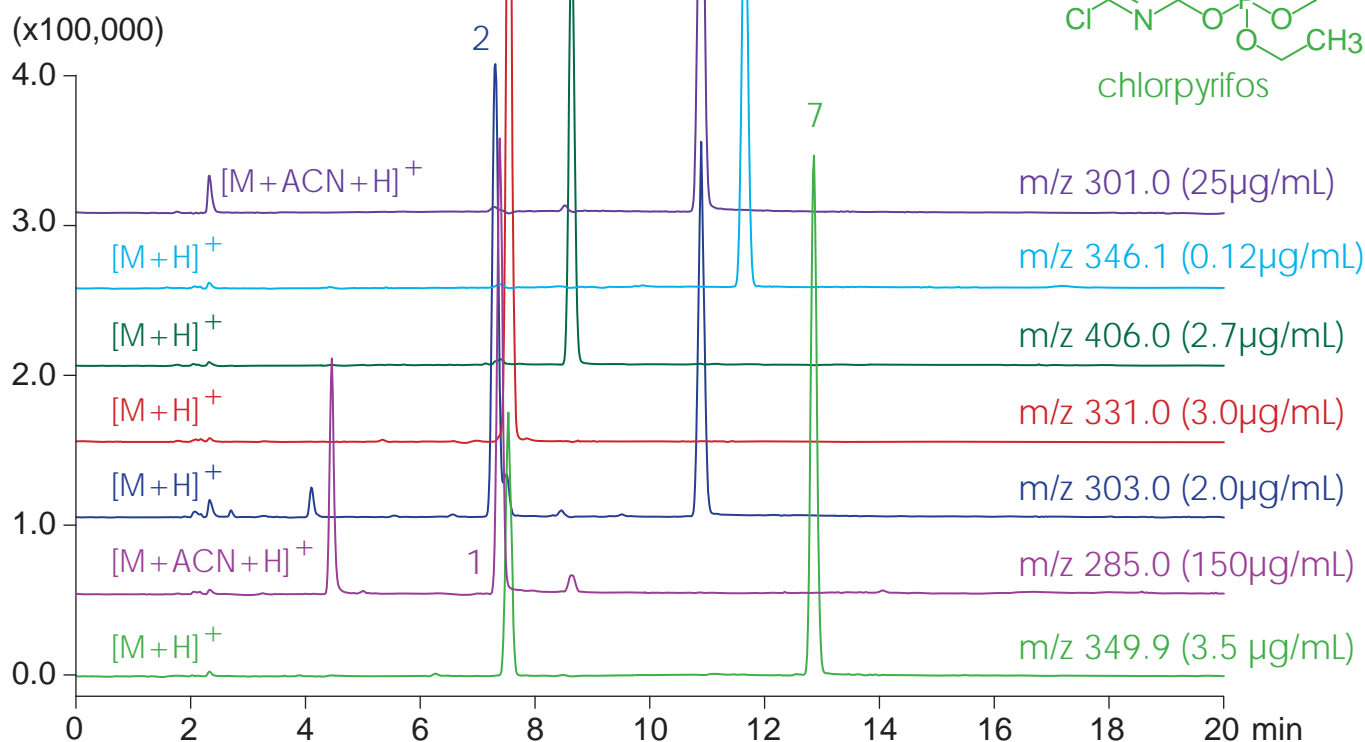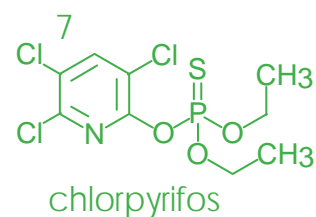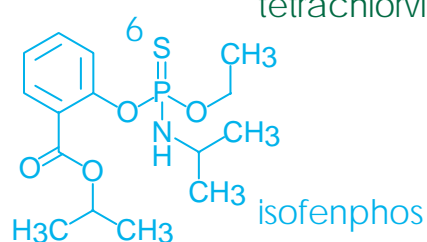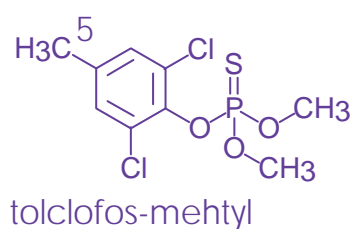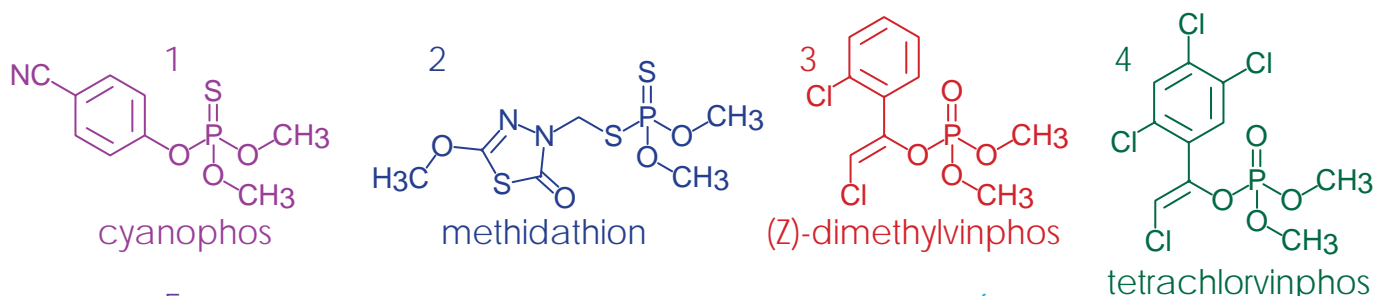

Unison UK-Phenyl, 150 x 3 mm
 A: water / formic acid = 100 / 0.1
 B: acetonitrile / formic acid = 100 / 0.1
 40-70%B (0-15min), 70%B (15-20min)
 0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)
 5µL (0.12-150µg/mL, ACN / water / HCOOH = 50 / 50 / 0.1)