


Unison UK-Phenyl

150 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Organophosphorus pesticides ( 2 )


有機リン系殺虫剤 ( 2 ) の LC-MS 分析

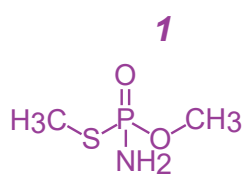


methamidophos

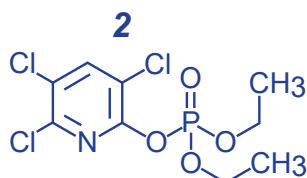


chlorpyrifos-oxon

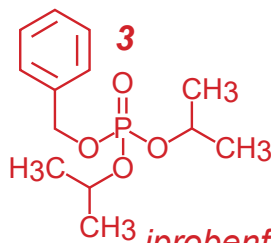


iprobenfos

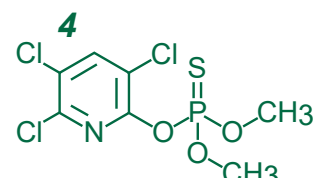


chlorpyrifos-methyl

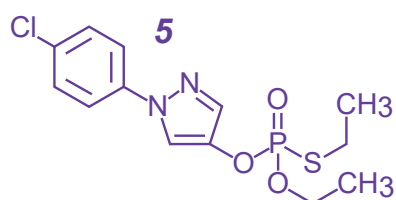

pyraclofos

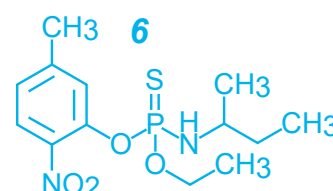

butamiphos

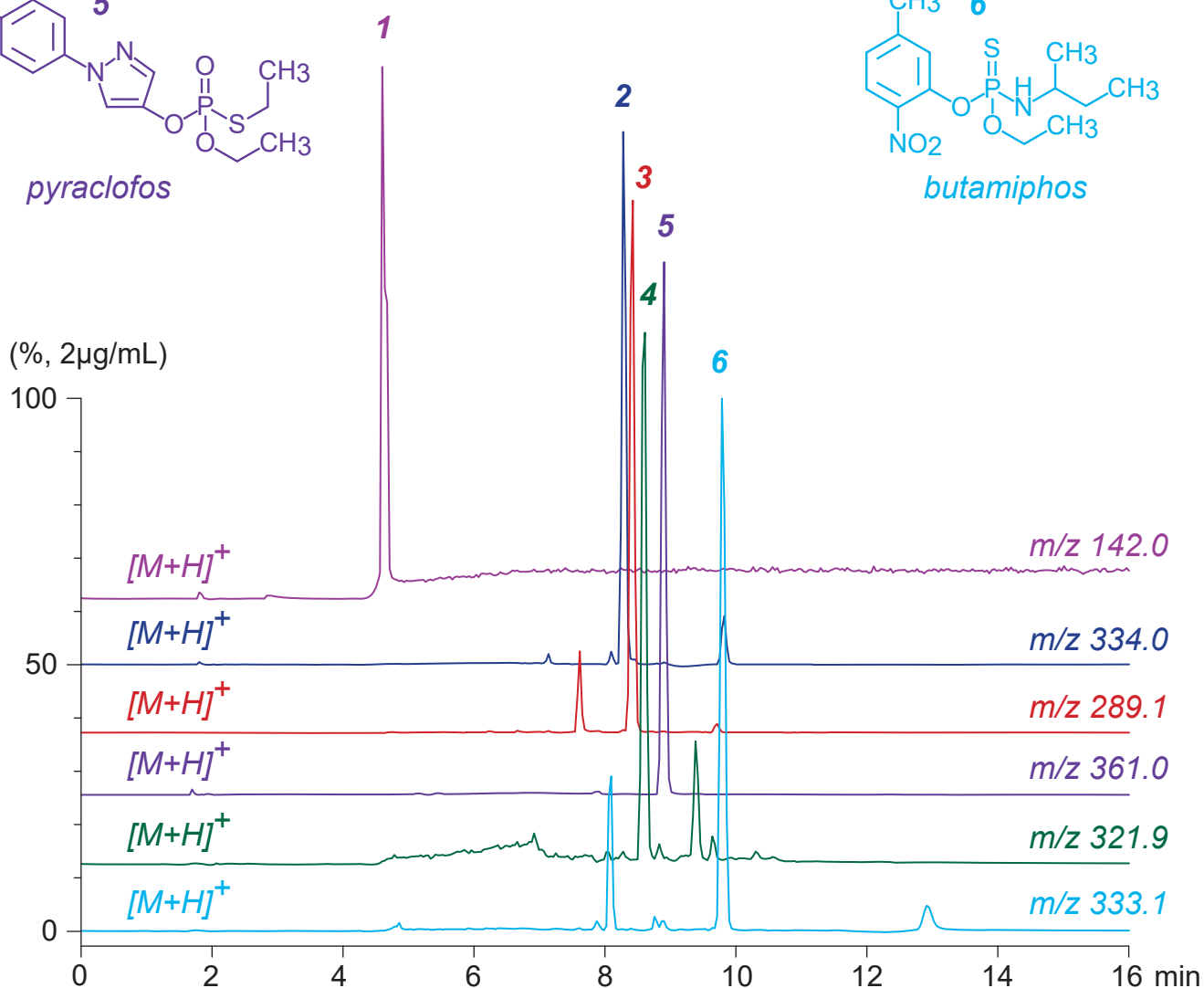

Unison UK-Phenyl, 150 x 3 mm  
 A: water / formic acid = 100 / 0.1  
 B: acetonitrile / formic acid = 100 / 0.1  
 0-65%B (0-5min), 65%B (5-16min)  
 0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)  
 5µL (2µg/mL, ACN / water / HCOOH = 50 / 50 / 0.1)