



Unison UK-Phenyl Cadenza CW-C18  
 Unison UK-C8 Cadenza CX-C18  
 Unison UK-C18 Cadenza CD-C18  
 Unison UK-C1

100 x 4.6 mm

HPLC COLUMN APPLICATION

### Trimetazidine ( Vasodilator )

トリメタジジン ( 血管拡張薬 )

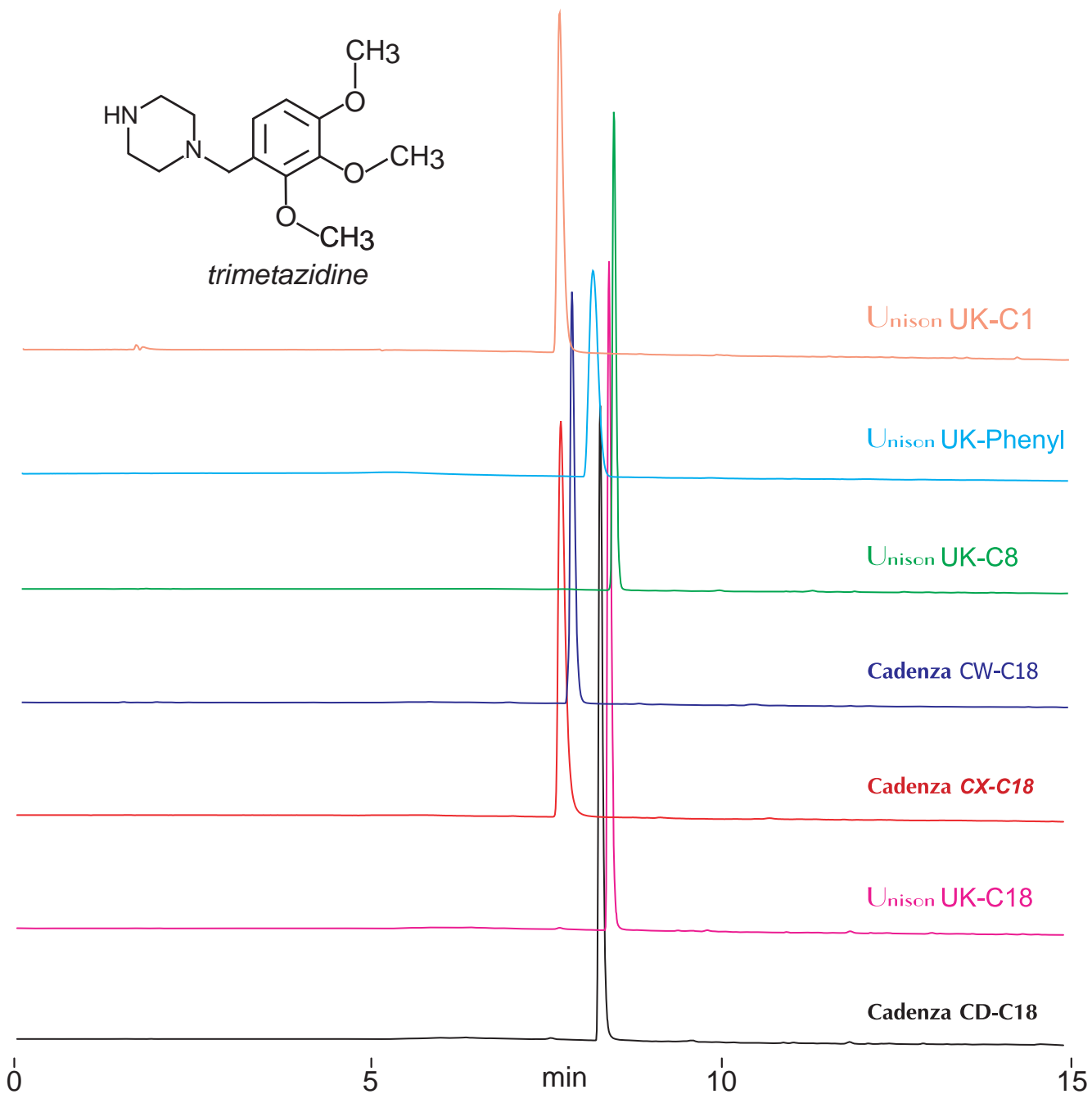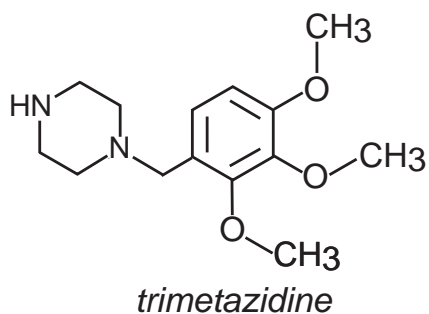

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

0-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 270nm

10µL (1mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)