

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Olmesartan medoxomil (Antihypertensive)

オルメサルタン メドキシミル (血圧降下薬)

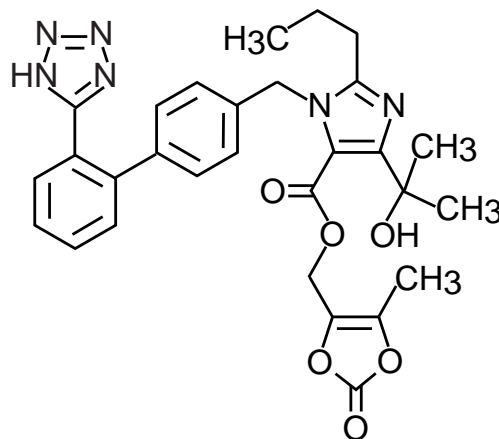

olmesartan medoxomil

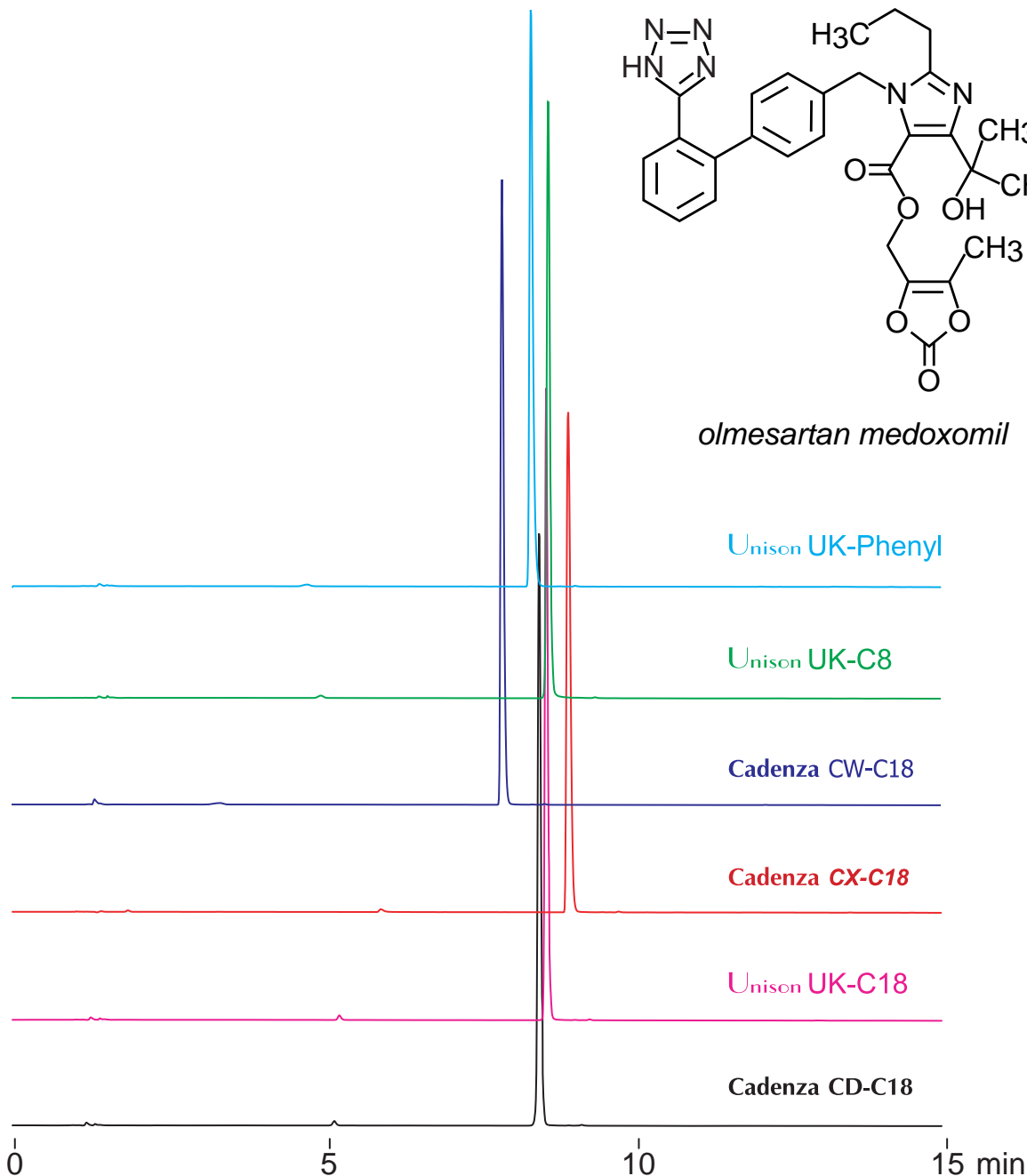

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

20-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1.0mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)