

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Troxeutin (Vascular protectant)

トロキセルチン (血管保護薬)

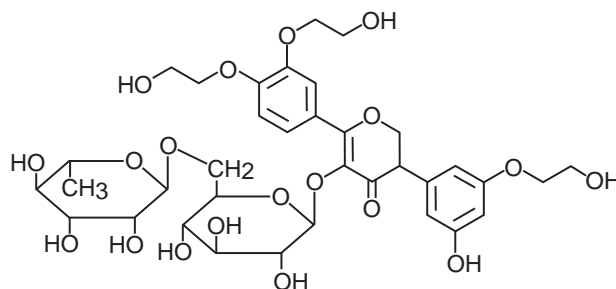

troxeutin

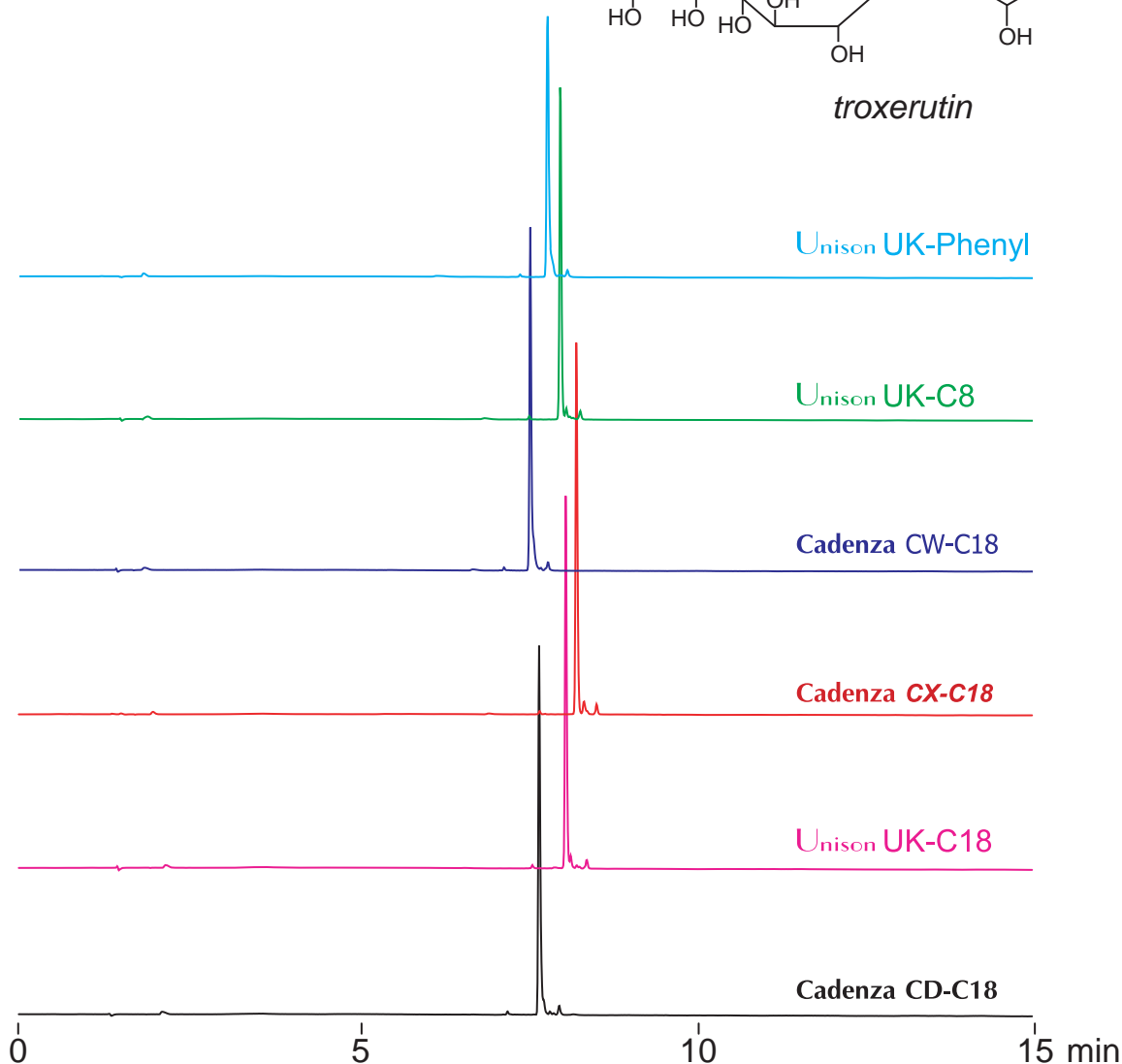

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-85%B (0-15min)

1.0mL/min (12MPa), 37°C, 360nm

10µL (0.1mg/mL, water / formic acid = 100 / 1)