

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

## Troloxerutin ( Vascular protectant )

トロキサセルチン ( 血管保護薬 )

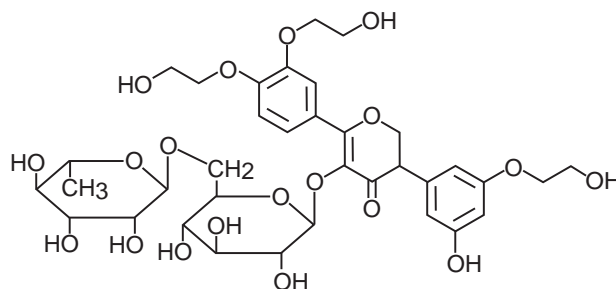

troloxerutin

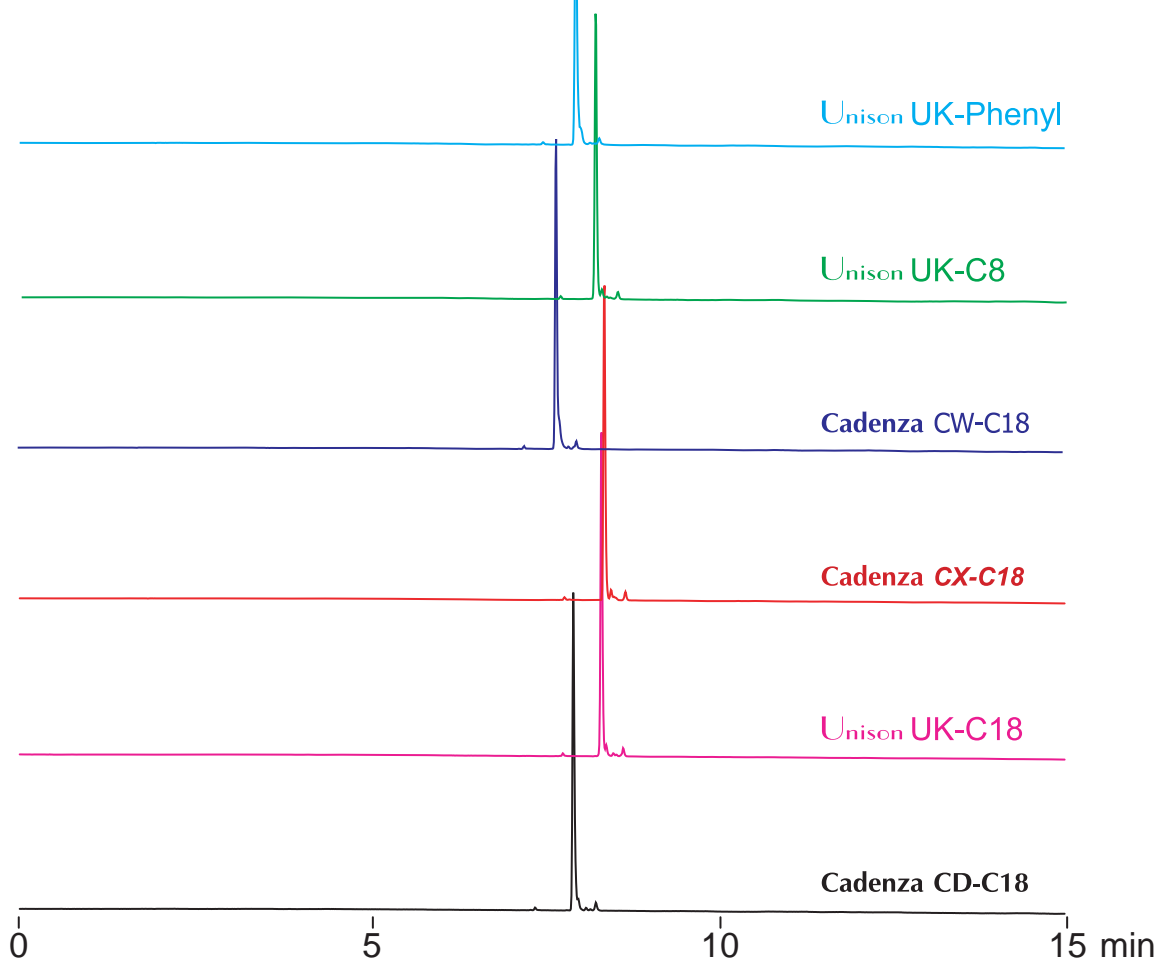

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-85%B (0-15min)

1.0mL/min (12MPa), 37°C, 360nm

10µL (0.1mg/ml, 100mM ammonium formate )